i Neck  
a ir **2**

breo

**breo**

## CONTENTS

Precautions	<b>03</b>
Product configuration	<b>08</b>
Product features	<b>09</b>
Product functions	<b>10</b>
Battery removal steps	<b>19</b>
Troubleshooting	<b>20</b>
Protection	<b>21</b>
Product information	<b>23</b>

## EN Precautions

To ensure safe and correct use of the product, please read the Operation Manual in detail before use. Our company will assume no liability for any consequences arising from use of the product improperly or without observing the instructions.

	Do not use the product during driving.
	Prior to use, please remove any jewelry or object that may affect massage from the massage position.
	Do not use the product in bathroom or other damp environments.
	Keep the product away from water and fire sources and corrosive environments.
	Please observe safety regulations when using the product in specific environments.
	Do not use the product if the massage position has undergone operation.

	Do not use the product if the massage position has skin injury, redness and inflammation, congestion or other discomfort symptoms.
	Persons with heart disease or history of heart disease, or with acute disease, skin disease, malignant tumor or physical abnormality should not use the product. Persons with bleeding tendency or other blood diseases, or with osteoporosis, soft-tissue disorders or apoplexia should not use the product. Persons with slow action, psychosomatic disorders or disturbance of intelligence should not use the product; please follow the doctor's advice if any of these persons intends to use the product. Stop using the product in case of any discomfort during use.
	Women who are pregnant or in menstrual period should not use the massager.

	This device is not intended for use by children under the age of 8. Children aged between 8 and 12 must be supervised while using the device.
	Use alcohol cotton to clean the external appearance; avoid using water for cleaning.
	Do not use any sharp objects to scratch or stab the product.
	<p>The product is an electronic wellness product for home use rather than a medical device. It is not intended for diagnosis, therapy or other commercial purposes.</p> <p>It must only be supplied at safety extra low voltage corresponding to the marking on the appliance.</p> <p>Cut off power supply immediately in case of any abnormality during use. Do not attempt to disassemble or repair the product.</p> <p>Models with moxibustion function will have a heated surface. Persons who insensitive to heat must caution for potential injury when using the appliance.</p>

5

**EN**

## NOTES

**Necessary knowledge****※About use frequency**

Use as needed, 10-30min each time. A proper interval is suggested before reuse.

Overuse may result in physical discomfort, or neck or skin injury.

Therefore, when the product is used for the first time, it is necessary to observe the physical condition and the neck, shoulder and skin changes, and reduce the use time as much as possible.

**※About suitable size**

The product is suitable for neck circumference between 30-40cm.

**※About usage**

When using the product, please close your eyes, lean back in the chair and totally relax your body.

Please remove your necklace and other jewelry, behind-neck

6

headphone, around-the-neck name badge, etc.


If you feel pain or discomfort around your neck or shoulders, please stop using the product and seek medical advice immediately.

#### ※About the status after use

Your neck and the areas around your neck may become red in a short time, but will be recovered after a period of time.

#### ※Miscellaneous

Sounds like "buzz", "puff" and "boom" will be heard during use. These sounds are produced when the massage head rotates and during air flow regulation, which should not be considered as abnormality of the product.

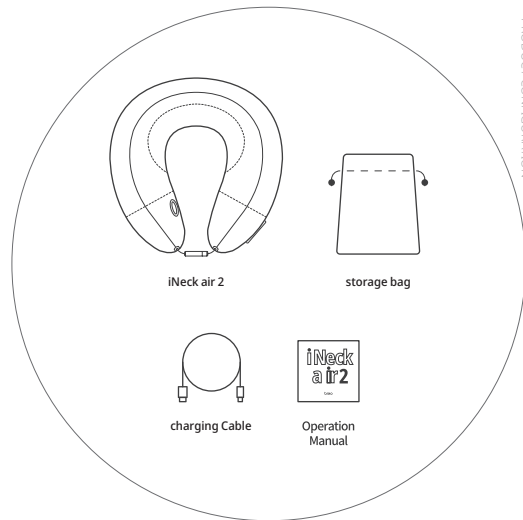Kneading  
massage

Warm  
Compress


Automatic  
inflation

Easy  
folding

9

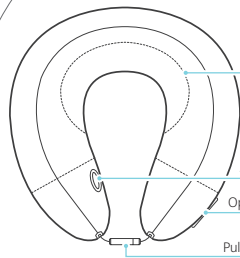
Main unit

Air vent valve


Operation area

Pull strap buckle

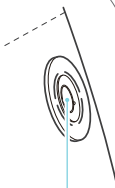Pressure relief valve

10

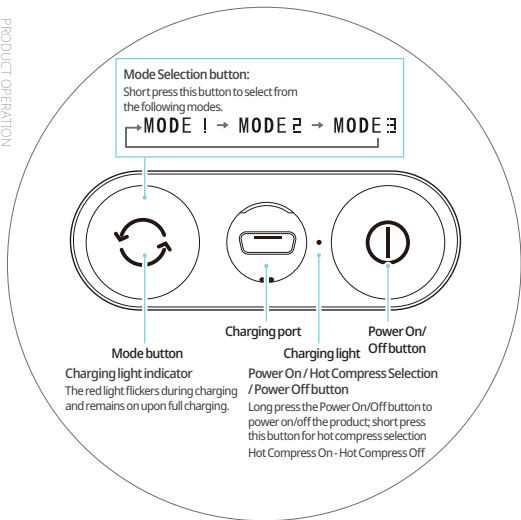11

EN

### Charging method

- 1 Open the external silicone plug in the port.
- 2 Insert the Micro USB connector on the USB connecting cable to the Micro USB port on the main unit.
- 3 Insert the power adapter to the power outlet.

The main unit produces the "beep" sound; then the charging light on the control panel flickers.

After charging, the charging light remains on.

※During charging, it is normal that the main unit and the power adapter get heated slightly.

### Charging time

- 1 It takes about 3.5 hours to fully charge the battery.
- 2 Upon full charging, the product can be used 10 times if it is used for 10 min (MODE 1, MODE 2) or 30min (MODE 3) every time.

(There might be slight difference due to the difference in ambient temperature)

- 3 Please charge the battery immediately when the charging


12

light flickers.

**After charging**

- 1 Unplug the power adapter, and disconnect the Micro USB connector from the main unit.
- 2 Close the external silicone plug in the port.
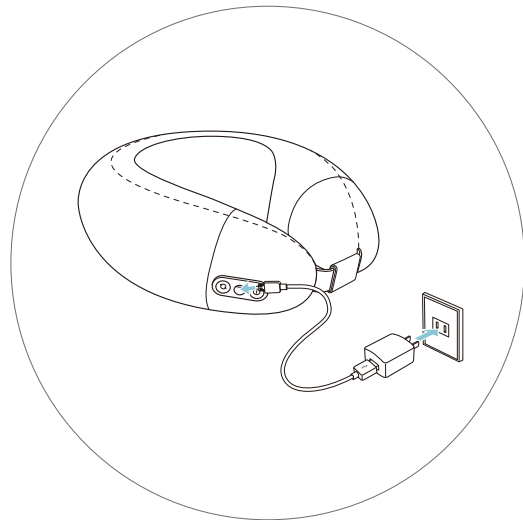
## EN **Wearing method**

### **1 Press the button in the middle of the pull strap buckles to separate the buckles**

The product is suitable for neck circumference up to 40cm.


※ To avoid damage, please do not press the middle button with force when the buckle is not fastened.

### **2 Hold the two ends of the main unit by both hands; pass from the top to the rear of your head to wear the product at your neck**

After wearing, the operation panel may be on your left or right side

You can distinguish the button positions and types by touch.

**Both the Power and Mode buttons use concave design; dot bulge exists beside the Power button.**
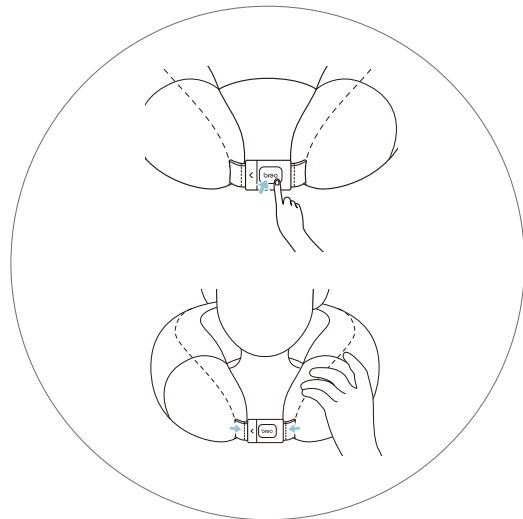
## EN Correct usage

The product should not be used while charging. After charging, unplug the power adapter before using the product.

**1** Long press the Power On/Off button to power on the product.

The main unit produces a "beep" sound; then the light indicator on the control panel comes on; MODE 1 is activated by default, and the Hot Compress function is enabled.

**2** As per your own need, you can select your favorite massage mode by pressing the Power On/Off and Mode buttons.

Upon startup, short press the Power On/Off button to separately disable or enable the Hot Compress function;

Upon startup, short press the Mode button to select among the following three modes;

**MODE 1 The total time is 10min**

massage is enabled; the massage rotation direction is changed at a fixed time;

**MODE 2 The total time is 10min**

massage is enabled; the massage head is retrieved at a fixed

time and is extended again for massage after a short interval;

**The total time is 30min**

massage is disabled; the massage head is retrieved.

**3** Wear the product at your neck. When using the product, please relax your neck and shoulders, lean back in the chair and totally relax your body. Or pull down the pull strap buckles to control the pressure exerted by the massage head on the neck.

**4** A long "beep" sound will be heard after use for 15min or 30min; then the light indicator on the control panel goes out, and power supply is cut off automatically.

**5** How to manually cut off power supply

Long press the Power On/Off button for over 2s; when a long "beep" sound is heard, power supply is cut off automatically.

**6** After use

Press the button in the middle of the pull strap buckles to separate the buckles; hold the two ends of the main unit by both hands to remove the main unit.

※2. Please observe the neck and skin changes when the product is used for the first time or by those with sensitive skin

## EN Operation steps for battery removal

- 1 Cut off the external power supply, and remove the product housing;
- 2 Disconnect the power supply to the battery;
- 3 Remove the battery;
- 4 Place the waste battery at a fixed recycling station;
- 5 Please observe local environmental protection laws and waste disposal guidelines and protect the environment.

### Notes:

- Remove the battery out of the appliance before scrapping the appliance;
- The appliance must be disconnected from power supply before removing the battery;
- The battery should be disposed of safely.

**Except professionals, no one shall attempt to remove the battery!**

## EN Troubleshooting

Please perform the following checks in case of any failure occurring during operation. If the failure cannot be eliminated or the product still cannot work ideally, please cut off power supply and ask a professional technician to repair it. Never attempt to disassemble the product.

Symptom	Possible cause	Basic solution
Power-on failure	Check whether the battery level is too low	If yes, please charge the battery timely
Low massage intensity	Check whether the battery level is low	If yes, please charge the battery timely

If the problem cannot be solved after troubleshooting, please send an email to [support@breocare.com](mailto:support@breocare.com).

EN

PROTECTION

### **Storage**

When the storage period exceeds 3 months, unpack the product, and check and recharge it before use. Do not place the product tightly close to floor, wall and ceiling. The indoor environment should be well-ventilated; keep away from strong sunlight and corrosive gas.

### **Transport**

The product package is provided with simple shockproof measures, therefore is suitable for transport by air, rail, road and sea. Keep away from rain and snow, and avoid inversion and collision.

### **Check of power supply**

Check the power supply; do not power on the product if the supply voltage is beyond the specified range.

### **Maintenance**

Take measures to resist moisture, rats and pollution. Cleaning should be carried out when the product is powered off. After use, put the device in a safe place; do not place any heavy objects on the device. In case of any failure, please

21

ask a professional technician to repair it. Do not attempt to disassemble the product.

### **Disposal and Recycle**

The product should not be abandoned along with domestic wastes. Dispose of the product according to local environmental laws and guidelines.

Please remember to respect the local regulations: hand in the non-working electrical equipments to an appropriate waste disposal centre. The packaging material is recyclable. Dispose of the packaging in an environmentally friendly manner and make it available for the recyclable material collection-service.

This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.

22

EN

Product name	breo iNeck air 2
Product model	neck 2b1
Rated voltage	5V 
Rated power	5W
Battery capacity	2300mAh
Product weight	700g
Product size	L310×W300×H100(mm)

**Shenzhen Breo Technology Co., LTD**

19/F, Hisense Southern Building, No.1777 Chuangye Road, Yuehai Street,  
Nanshan District, Shenzhen 518054, China  
www.breo.com support@breocare.com

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