
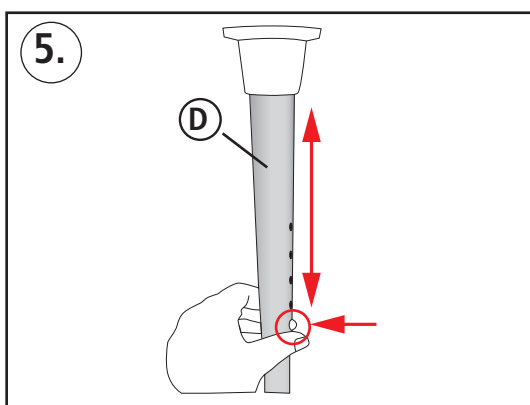
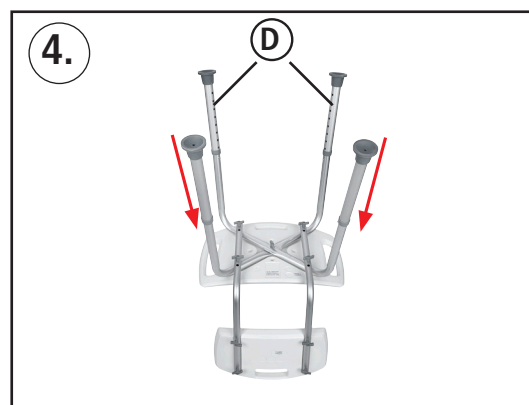
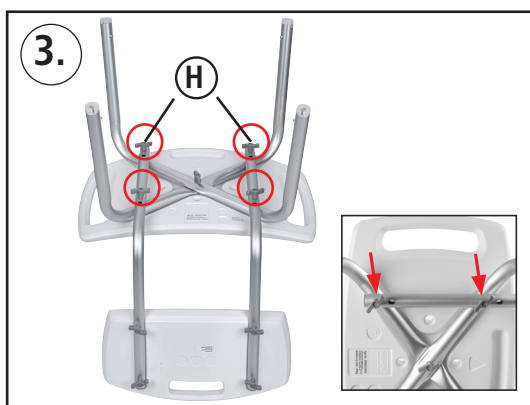
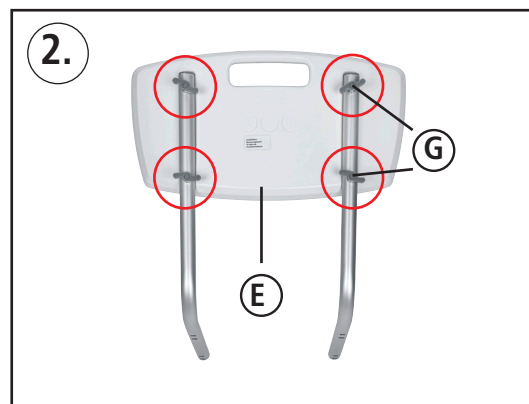
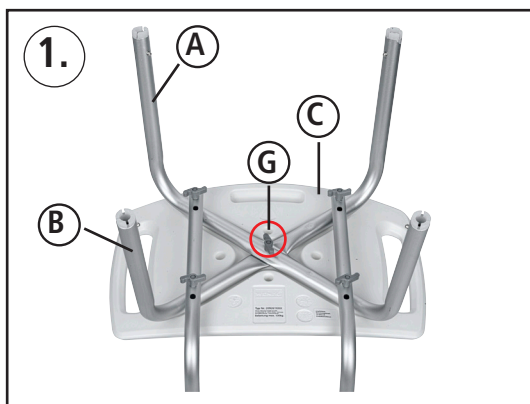
- Ⓛ Montageanleitung
- Ⓜ Assembly Instructions
- Ⓝ Instructions de montage
- Ⓔ Instrucciones de montaje
- Ⓜ Istruzioni per il montaggio
- Ⓝ Montage-instructie

- Ⓟ Instruções de montagem
- Ⓝ Montagevejledning
- Ⓜ Monteringsanvisning
- Ⓟ Instrukcja montażu
- Ⓜ Οδηγία μοντάζ
- Ⓝ Návod k montáži

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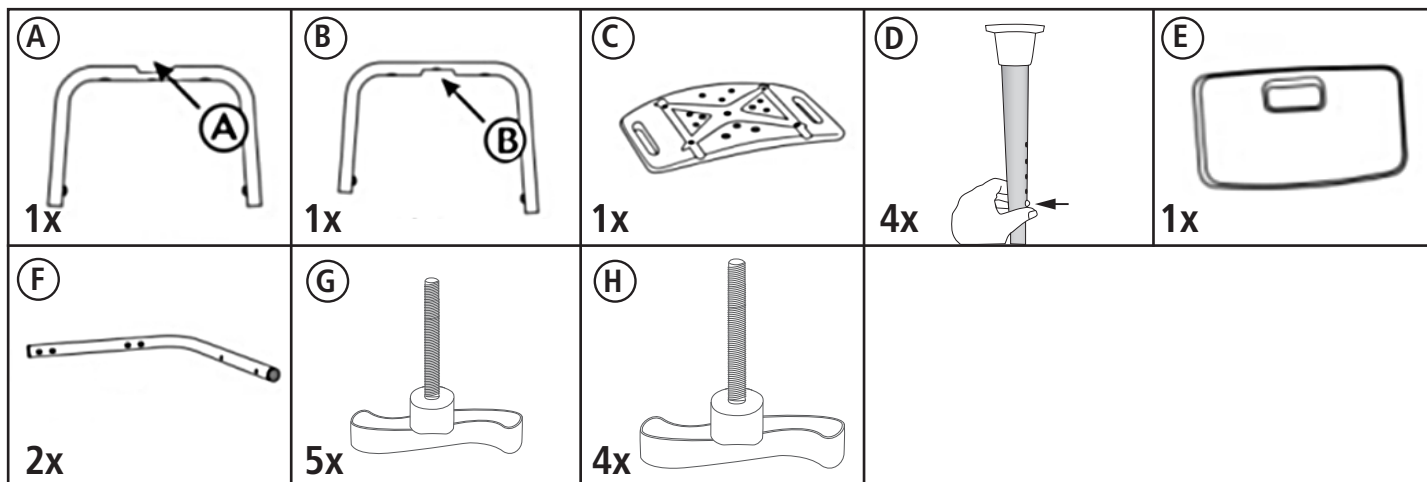Ⓢ Bathroom stool with backrest Model 209261*

Assembly Instructions: 1. Place frame B with notch B into the well on the rear side of the seat section C. Place frame A with notch A onto the frame with notch B and push a short wing screw G through the middle hole, then fasten the two frames to the seat section. 2. Screw the two L-shaped supports F to the backrest E using the short wing screws G. 3. Fasten the backrest to the seat section by aligning the corresponding holes of the L-shaped supports to the holes on the two frames. Now screw in the 4 long wing screws H. You can determine the depth of the seat by selecting one of two holes provided on the L-shaped supports. 4. & 5. Push the legs D into the frame, adjusting the height as you do so, by pressing the buttons and allowing them to engage in the appropriate holes!

CAUTION! Including mounting material. The enclosed mounting elements must be used as described in the

installation instructions. The seat is height adjustable by means of the buttons on the legs. Please make sure that all legs have the same height setting before using the seat. Before subjecting the seat to loads, please ensure that all screws have been tightened and that all locking buttons are engaged. max. load capacity 130 kg. Recommended service life: 15 years from production date.

Areas of use: This product may only be used in private households! Not suitable for use in rehabilitation facilities and for the medical care of persons. This product is intended for use in the wet area of the shower.

Care instructions: Clean with a damp cloth. Do not use any sandy abrasive agents or scrub sponges for cleaning.

Safety instructions: Pay attention to a firm hold! Smooth or wet surfaces can worsen the non-slip properties. Check connecting elements regularly and re-tighten screws if necessary. The seat should be checked regularly for stability and damage, and reuse is only permitted after the appropriate repairs have been carried out. Please keep these instructions, as they contain important information. Do not step on the seat.

Do not use to lift persons! Accidents resulting in serious injuries may occur if the handles are used for incorrect purposes e.g. as gymnastic equipment or steps.

Wenko-Wenselaar GmbH & Co. KG shall not assume any liability whatsoever for incorrect use!

Exclusions: Wenko does not assume any liability for damage to persons and property caused by improper mounting, incorrect use of the handles, in particular by non-observance of the warning and safety instructions listed above.

For after-sales service please contact the manufacturer directly:

Wenko-Wenselaar GmbH & Co. KG
Im Hülsenfeld 10 · 40721 Hilden, Germany
Tel: +49 2103 573-0 · Fax: +49 2103 573-190
service@wenko.de · www.wenko.de

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