PREMIUM NM-089 BEIGE

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ENGLISH

Contents

Introduction	14
Traditional shiatsu massage.....	14
Intended use	14
Product description	14
Parts of the product	15
Product functions.....	16
Product positions	17
Safety precautions	18
Specific group of people	19
Contents included	19
Technical specifications.....	20
Maintenance and storage	20
Electric and electronic devices	21
Customer service	22
Certificates and recommendations	22
Warranty policy.....	23

Introduction

Dear customer,

Thank you for choosing Donnerberg Original Premium neck and shoulder massager.

Please read the manual instructions carefully and save it for later use, be sure to make them accessible to the other users. Please follow the instruction and use the product as described.

Donnerberg Team

Traditional shiatsu massage

Donnerberg Original Premium massager is based on the principles of the Shiatsu massage that originated in Japan. This massage supports the body's natural ability to self-heal. The word Shiatsu means 'finger pressure' in Japanese. Shiatsu uses pressure, touch and manipulative techniques to bring the body and its inner energy paths into balance. Shiatsu massage can be used to reduce stress and promote overall health and well-being.

Intended use



The unit is only intended for domestic /private use, not for commercial use. It is designed for massaging the following body parts: neck, back, shoulders, belly, legs and arms. If you are unsure whether the neck and shoulder massager is suitable for you, please consult your doctor.

Product description


Sit in a comfortable chair, carefully place the massager over your shoulders and switch it on. Put the arms through the loops to facilitate the positioning of the massager and to control the intensity of the massage. Place the massager on the desired areas or parts of the body and adjust the intensity of the massage by pulling harder or softer on the loops.

Each massage cycle lasts 15 minutes but can be shortened on demand. When the 15 minute cycle has elapsed, the massager will automatically switch off. This device must be turned off when not in use.

Parts of the product

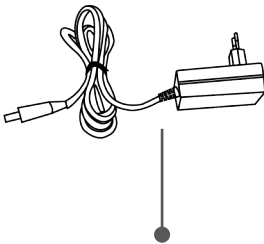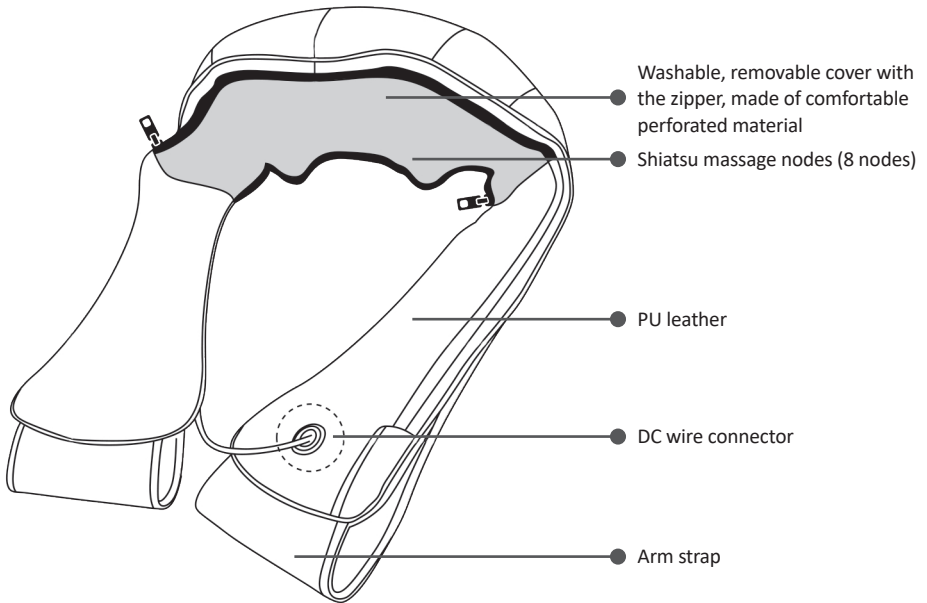
Power adapter

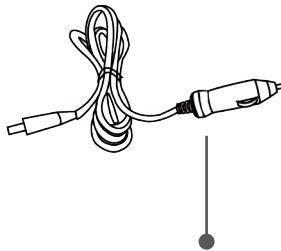Car power adapter

Product functions

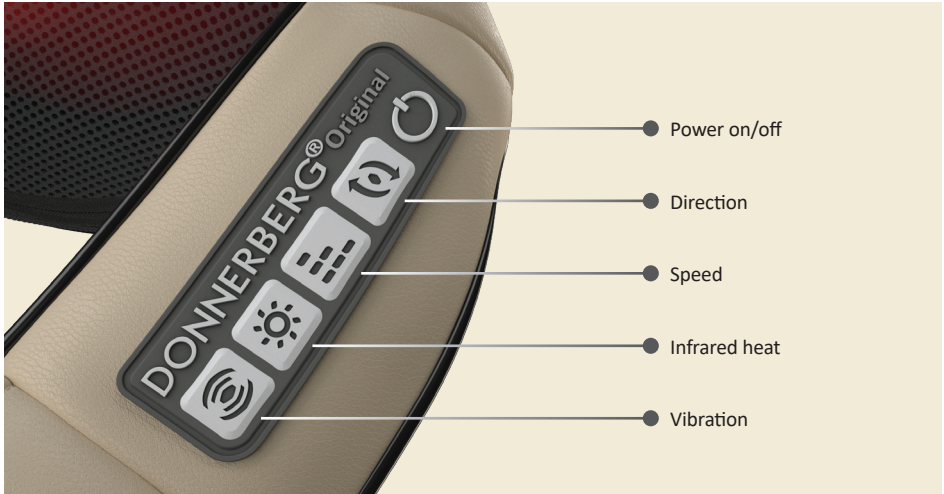1. Power on/off

- Press this button to turn on the device.
- This device is equipped with an infrared heating function which automatically turns on by pressing the “Power” button along with the massage function and rotation.
- Press the “Infrared heat” button to turn off the heating function.

2. Direction

- Press the “Direction” button to change the rotating direction.

3. Speed

- Choose between 3 massage speed modes. The lowest speed is preset by default when the massage device is switched on.
- Pressing the “Speed” button increases the speed of the massage.
 - Press 1x for the medium speed.
 - Press 2x for the highest speed.

4. Infrared heat

- You can choose between two levels of infrared heat:
 - When the device is switched on, the 2nd level of the Infrared heat is on.
 - Press 1x to switch to the 1st level of infrared heat.
 - Press 2x to turn off the infrared heat.
- You can also use infrared heat independently of the massage.
- To activate only infrared heat, switch off the device and press the “Infrared heat” button.

5. Vibration

- Turn on the additional vibration function by pressing the “Vibration” button in case you would like more intensive massage.
 - Press 1x to turn on vibration. Default vibration mode - continuous vibration.
 - Press 2x to change the speed and mode of vibration. Vibration lasts for 5 seconds, then pauses for 2 seconds repetitively.
-
- The vibration function is only recommended on the shoulders, back and leg area.
 - The massage should feel pleasant and relaxing at all times. If the massage is painful or at any time uncomfortable change the pressure using the straps.

Product positions

Please apply on body parts shown in the pictures below.


Neck & shoulders


Back

Belly

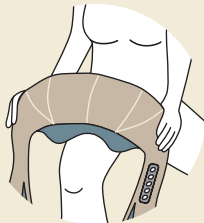Legs

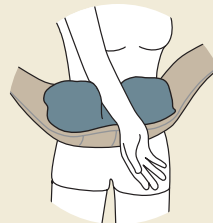Arms

Safety precautions

To prevent the risk of fire, electric shock, or other injury, please follow the instructions:

1. Before connecting the device to your power supply, please ensure that the supply voltage stated on the rating plates is compatible with your mains supply. Only insert the power supply unit into a socket when the appliance is switched off.
2. Do not cover or use the unit under a blanket when operating. Do not use while bathing or showering. The unit must be kept dry and stored in a dry and safe place. The power socket, cable, and appliance must be at appropriate distance from the heat, hot surfaces, moistures and liquids. Connect the unit only in that way, so that the power plug is freely accessible. The power plug should be freely accessible when inserted into the socket. The cable must not be kinked, trapped or twisted.
3. Change the covers only when the device is turned off and the power supply unit unplugged.
4. To use this product safely and correctly, read this manual thoroughly and operate the device as instructed. Never operate this appliance if it has a damaged cord or plug or if it is not working properly. In case the device is not working or has a fault, send it to the service centre for examination. Do not plug in the device in the damaged sockets. If the device has been dropped into water, do not pick it up or touch it. Pull out the power plug from the socket immediately.

Do not operate where aerosol (spray) products are being used or where oxygen is being administered. After you have finished using the unit turn off the power. Before unplugging, the device must be turned off. Do not obstruct with objects or thrust.

Please ensure that the polythene packing is kept away from the reach of children!
Risk of suffocation and choking hazard!

Do not use the massager while driving!

Children are not allowed to play with the device!

Specific group of people

- This appliance can be used by children aged 8 and over and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning the use of the appliance in a safe way and understand the hazards involved.
- Do not use this device to supplement or replace medical treatment.
- Chronic complaints and symptoms may be made worse.

You should refrain from using the Shiatsu massager, or should consult your doctor before using if:

- You are pregnant,
- You have a cardiac pacemaker, artificial joints or electronic implants,
- You suffer from one or more of the following illnesses and complaints: Circulation disorders, varicose veins, open wounds, bruises, broken skin, inflammation of the veins.

- Do not use the appliance near your eyes or near any other sensitive areas.
- Take care when using the heating function. Do not use it if you have circulatory problems or on areas sensitive to heat.
- The surface of the appliance is warm. People who are sensitive to heat must be careful when using the appliance.
- If you experience any pain or discomfort during the massage, stop immediately and consult your doctor.
- Before using the massager, consult your doctor regarding any unexplained pain, if you are undergoing medical treatment or using medical equipment.

Contents included


1	Donnerberg Original Premium Neck & Shoulder massager NM-089
2	Removable covers in different colours
1	Power adapter (length approx 2.8m)
1	Car charger
1	Soft cloth cover for sensitive skin
1	Protective custody - environmentally friendly
1	User manual

Technical specifications

Product name	PREMIUM NM-089
Product model	NM-089 Black
	NM-089 Beige
Voltage supply	Output DC 12V 2A 24W
Input voltage	100-240 V, 50/60 Hz
Automatic switch off	After approx. 15 minutes
Net Weight	Approx. 1.85 kg
Gross Weight	Approx. 2.46 kg
Product dimensions	Approx. 50 x 43 x 17.6 cm

Maintenance and storage

- Before cleaning, make sure the switch is turned off and power plug is disconnected from the outlet. Do not use abrasive cleaners such as benzene or thinner to clean the product, otherwise the colour of the product may fade. Do not immerse the appliance in water or the other liquids during the cleaning! Small marks can be removed by using a damp cloth. Please remember to clean only the top side of the device.
- You may only clean the massager yourself. If the device malfunctions please do not repair it yourself as this leads to automatic expiry of all guarantee rights.
- Thanks to the integrated zipper, the cover can be easily removed and washed in the washing machine at 30° C. Simple and easy, just unzip the zipper and replace the cover with another one.

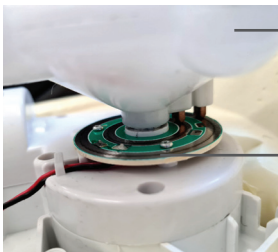● Massage node

● Contact plate

- We would kindly advise you to clean your massager regularly by taking off the cover and removing dirt (hair, dust) that gathered between the massage node and the contact plate with cotton buds. Regular removal of contamination can prevent malfunctions of certain components of infrared heat and massage.

Electric and electronic devices

Information for private households

The Waste Electrical and Electronic Equipment Regulations 2018 (WEEE) in the United Kingdom contains several requirements concerning the handling of electric and electronic devices as waste. In the following paragraphs, we list the most important requirements for you.

1. Separated recycling of „old devices“

Electric and electronic devices that have become waste are called old devices. Owners of such old devices have to bring these to recycling facilities, where these are handled and recorded separately from urban solid waste. Old devices must not end up in the household refuse. In fact they have to be collected in special collecting and recycling systems.

2. Possibility for disposal of „old devices“

Owners of old devices from private households can dispose these at the public waste disposal authorities or bring them to producer or seller related disposal locations. You can find these disposal locations online: <http://www.recycle-more.co.uk>

3. Privacy information

Old devices often contain private data, especially appliances from data processing and communication equipment, as smartphones or hard discs. In your own interest, please arrange to eliminate all personal and private data from your disposed equipment.

4. Meaning of the symbol „crossed out waste container“

The labelled symbol showing the „crossed out waste container“ on electric and electronic equipment indicates that the concerning equipment has to be disposed separately from urban solid waste.

5. The Waste Electrical and Electronic Equipment Regulations 2018 (WEEE)

Following the implementation of the Directive 2012/19/EU and its British transportation WEEE Regulations 2018, we as a producer of EEE are registered in the Public register of Electrical and Electronic Equipment Producers. The register is held in England and Wales by the Environment Agency.

Customer service

We have a reliable customer service team that is committed to providing solutions to any issue you might have. If you have any questions or concerns about the products, warranty or exchange do not hesitate to contact us at:

E-mail	de.info@donnerberg.net
Telephone	+ 49 89 215 400 980 / + 49 95 619 769 199
<p>Company Donnerberg Von-Werthern-Str. 4a 96487 Dörfles-Esbach Germany</p>	

Certificates and recommendations

 	<p>Recommended by the German Olympic Sports Confederation</p>
	<p>TÜV SÜD (Technical Inspection Association)-the product was tested for electrical safety, hazardous material components, and unwanted noise production. Certificate No.: Z2 091698 0029 Rev.00</p>
	<p>RoHS (Restriction of Hazardous Substances) - restricts the use of specific hazardous materials found in electrical and electronic products 2011/65/EU</p>
	<p>CE (Conformite Europeenne)- indicates that the product conforms with the relevant EU directives regarding health and safety or environmental protection (765/2008)</p>

Warranty policy

The company Donnerberg, Mr. Svetozar Grbic, Von-Werthern-Str. 4a, 96487 Dörfles-Esbach, Germany grants a valid manufacturer's warranty for the territory of the European Union and the civil-law defect liability. The warranty starts from the time that the device was delivered to the customer. It is not extended by a previous warranty case. The warranty is granted only to the buyer and is not transferable to third parties.

In addition to the statutory liability for defects, we also repair or replace equipment whose defects are not initial defects unless they are caused by improper use or are insignificant. The choice between replacement and repair is ours. The warranty is only granted if the defect is notified immediately upon its occurrence and is reported.

You may report your warranty claim to one of the following:

**Company Donnerberg
Von-Werthern-Str. 4a
96487 Dörfles-Esbach
Germany**

**e-mail:
de.info@donnerberg.net**

**or telephone:
+ 49 89 215 400 980
+ 49 95 619 769 199**

Upon your notification, our customer service team will send you a claim form.

You must include the claim form and the copy of your purchase receipt as the proof of warranty. Please note that we will not process warranty claims without this claim form.

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