POWER GROW COMB

by Velform®

Instruction Manual

Mode d'emploi

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Manual de instrucciones

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**IMPORTANT INSTRUCTIONS
RETAIN FOR FUTURE USE.**

The power Grow comb by Velform® includes –
 A - The power Grow Comb: Hair growth unit.
 B - Massage comb: Relaxing rubber massage comb to use in the shower.
 C - Manicure set: Luxury pampering travel kit.
 D - Hair brush: Professional hair salon brush.
 E - Instructions

Introduction

Congratulations you have purchased the most advanced hair growth comb on the market today. Please read and follow the instructions carefully.

Information about the Power Grow Comb by Velform®

The Power Grow Comb by Velform® uses Lasers and infrared LED lights to help stimulate scalp cells. This ground breaking therapy also helps improve blood circulation. The combination of Lasers and infrared-red lights helps to nourish hair follicles.

The Power Grow Comb by Velform® also stimulates hair follicles by vibration technology. The patented vibration therapy will help increase blood flow to your hair follicles and in turn help improve hair growth.


The attached comb also enables the Laser and light therapy to reach your hair follicles. The teeth of the combs are positioned before or after the Laser and LED's to create perfect hair separation.

Through scientific investigation, engineering and design the Power Grow Comb by Velform® will help grow your hair!

INSTRUCTIONS


Inserting the Batteries

1. Open the battery cover.
2. Insert 4 AAA size batteries as shown:

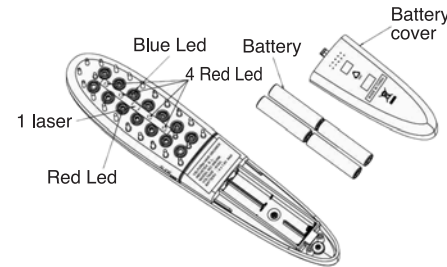Turning the Power Grow Comb by Velform® on.

To turn on the unit press the power switch by the logo. You will hear a short beep and the LED by the power switch will glow red. This means the Power Grow Comb by Velform® is on.

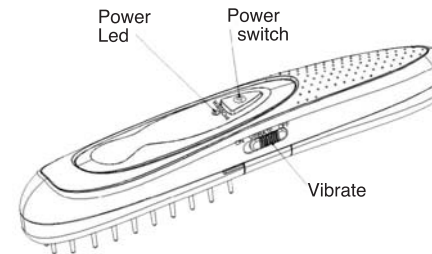
Turning on the vibration mode.

The Power Grow Comb by Velform® also features a vibrate mode. On the side of unit, move the VIBRATE switch to the 'ON' position. The vibration mode will then be activated. Note – the vibration mode will only work when the unit is on and the power button has been pressed.

Detaching the Hair Comb

If you desire, it is possible to detach the comb of the Power Grow

Comb by Velform® and apply the laser and infrared lights directly over your scalp. To detach the comb please locate the 'slide' labels at either side of the comb. Once located please slide the comb away from the unit. To re-attach the comb please click the unit back into its original position. You will hear a 'clicking' sound.

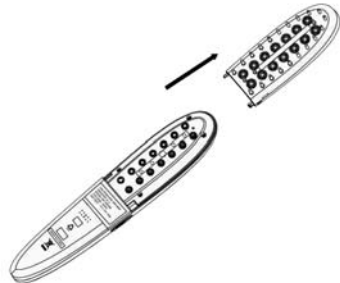Directions of Use

You should use the Power Grow Comb by Velform® without applying any lotions or gels.

Use Power Grow Comb according to your own schedule and convenience: during 5 or 10 minutes a day, or during 10 to 15 minutes every other day. You will begin to notice results in approximately 8-12 weeks.


Please move the unit slowly around your head in the direction of hair growth, from front to back, and bottom to top. Please make sure you cover all of your head and move in a slow controlled motion.

Massage Comb

The Massage Comb is especially designed to help circulate the blood flow in your scalp and improve your hair growth.

It is possible to use the Massage Comb with either wet or dry hair. For an exhilarating experience you may also use the Massage Comb in the shower to help apply shampoo.

Please place your fingers in the holes of the Massage Comb. (See Figure)
Move the massage comb around your scalp in

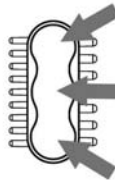a circular and controlled motion.

SPECIFICATION

Power Grow Comb by Velform® (model n°: WB/096)

Input voltage

6V DV (4 x AAA size battery)

Laser wavelength

650nm +/- 5nm

Laser class

4.5 mW Class 3R certified

Red LED wavelength

660 nm +/- 5nm

Blue LED wavelength

470 nm +/- 5nm

Output

420mW without vibration

600m W with vibration

The Power Grow Comb by Velform® complies to the following regulations -

ROHS Certificate.

CE Certificate.

EMC: EN 55014-1 y EN 55014-2

WARRANTY

This product is covered by a warranty against manufacturing defects subject to the time periods stipulated in the legislation in force in each country.

ATTENTION

DISPOSING OF MATERIALS

The symbol of a crossed out rubbish bin with wheels indicates that you must be aware of and follow local regulations with regard to the disposal of this kind of product.

Do not dispose of this product the way you would dispose of normal household waste.

You must follow the relevant local regulations.

Electric and electronic products contain hazardous substances that have harmful effects on the environment and human health and should be recycled properly.

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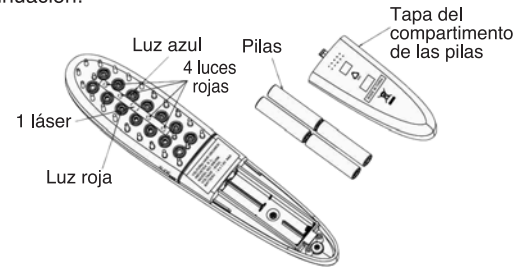
Encender Power Grow Comb by Velform®

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