## Digitale Personenwaage Gina

Digital bathroom scale Gina | Báscula de baño digital Gina |  
Pèse-personne numérique Gina | Bilancia digitale pesapersona  
Gina | Cyfrowa waga osobowa Gina

Deutsch

English

Español

Français

Italiano

Polski

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BE 2011

Operating instructions | Instrucciones de operación | Mode d'emploi  
Istruzioni per l'uso | Instrukcja obsługi

# Liebe Kundin, lieber Kunde !

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# Dear Customer!

You have chosen to purchase a high-quality ADE brand product, which combines intelligent functions with an exceptional design. This bathroom scale makes it possible for you to determine and monitor your body weight. The wealth of experience of the ADE brand ensures a high technical standard and proven quality.

We hope you enjoy weighing yourself and wish you success in reaching your goal weight!

Your ADE Team

## General information

---

### About this operating manual

This operating manual describes how to operate and maintain your product safely. Keep this operating manual in a safe place for future reference. If you pass this product on to someone else, pass on this operating manual as well.

Failure to observe this operating manual may result in injuries or damage to your product.

### Explanation of symbols

This symbol combined with the word DANGER warns of serious injuries.

This symbol combined with the word WARNING warns of moderate and mild injuries.

This symbol combined with the word NOTE warns of material damage.

This symbol indicates additional information and general notes.

# Table of Contents

---

General information .....	12
Safety .....	14
Intended use .....	15
Overview .....	16
Scope of delivery .....	16
Start-up .....	17
Inserting/changing batteries .....	17
Switching on voice playback.....	17
Choosing the unit of measurement .....	17
Setting up the scale .....	18
Weighing .....	18
Resetting the scale.....	18
Cleaning .....	18
Fault/Remedy.....	19
Technical data.....	19
Declaration of Conformity.....	20
Warranty .....	20
Disposal.....	20

# Safety

---

- This product may be used by children aged 8 and above, as well as by persons with reduced physical, sensory, or mental capacities or lacking experience and knowledge if they are supervised or have been instructed in the safe use of the product and understand the risks that may arise.
- Children must not play with the product.
- Cleaning and user maintenance must not be carried out by children, unless they are 8 years of age or older and are being supervised.
- Batteries, when swallowed, can be life-threatening. You should therefore keep the batteries and scale out of reach of small children. If a battery has been swallowed, medical assistance must be sought immediately.

## **DANGER for children**

- Keep packaging material away from children. There is a risk of choking if swallowed. Also pay attention that children do not pull the packaging material over their heads.

## **RISK of fire/burning and/or explosion**

- There is a risk of explosion if the batteries are replaced improperly. Replace the batteries only with the same or equivalent battery type. Pay attention to the "Technical data".
- Batteries must not be recharged, reactivated by other means, dismantled, thrown into fire, or short-circuited. Otherwise, there is an increased risk of explosion and leakage; furthermore, gases may escape!
- Avoid contact of the battery fluid with skin, eyes, and mucosa. In case of contact, immediately rinse the affected areas with plenty of clean water and consult a doctor.

## **WARNING of risk of injury**

- Do not step onto the outer edge of any one side of the scale. Risk of tipping!
- Do not step onto the scale with wet feet or when its surface is damp. Risk of slipping!
- If the scale falls down and shatters, there is a risk of injury from glass splinters.

## BEWARE of material damage

- Place the scale on a stable, even surface.
- Do not overload the scale. Pay attention to the “Technical data”. Place yourself on the scale exclusively for weighing purposes.
- Protect the scale from strong shocks and vibrations.
- Remove the batteries from the scale when they are depleted or when you are no longer using the scale. You will thus avoid damage that may result from leakage of battery acid.
- When inserting and changing batteries, pay attention to the correct polarity (+/-).
- Do not subject the batteries to any extreme conditions, by – for example – storing them on radiators or under direct sunlight. Increased risk of leakage!
- Clean the battery and device contacts, as needed, before inserting the battery.
- Do not make any modifications to the product. Only allow repairs to be carried out by a specialist workshop. Repairs carried out improperly can considerably endanger the user.
- Do not place the scale in the vicinity of devices that may emit electromagnetic radiation, such as mobile phones, wireless devices, or the like. Doing so may lead to incorrect display indications or malfunctions.
- There are feet on the underside of the scale. Since floor surfaces are made from various materials and are treated with diverse care products, it cannot be completely excluded that some of these substances contain ingredients that may adversely affect and soften the feet of the device. If necessary, place a non-slip pad under the scale.

## Intended use


---

The bathroom scale is designed for the weighing of people in the home. It is not suitable for use in clinics, medical practices, etc.

The maximum load capacity of 250 kg must not be exceeded.

# Overview

---

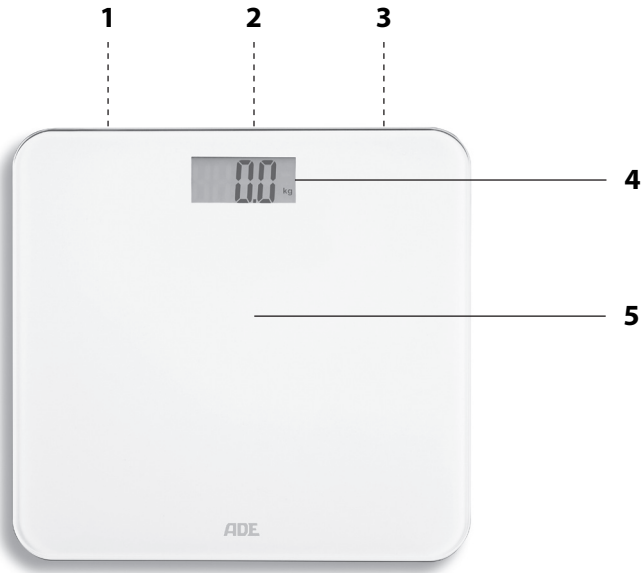- 1 Battery compartment (underside)
- 2 Adjusting voice playback and volume
- 3 Speaker
- 4 Display
- 5 Weighing surface

## Scope of delivery

---

- Bathroom scale
- Battery type LR03 (AAA), 1.5 V  (3x)
- Operating Manual

# Start-up

## Inserting/changing batteries

For use you need 3 batteries of type LR03 (AAA), 1.5 V (included). The battery compartment is located on the underside of the device.

1. Push the locking tab of the battery compartment cover towards the battery compartment and remove the cover.
2. Insert the batteries as indicated at the bottom of the battery compartment. Pay attention to the correct polarity (+/-).
3. Replace the battery compartment cover.  
The locking tab must snap audibly into place.

Suggestions for changing the battery:

- A change of batteries becomes necessary when “Lo” appears on the display, or when nothing at all is displayed.
- Only use the type of battery that is specified in the “Technical data”.
- Always change all batteries at the same time.
- Dispose of old batteries in an environmentally appropriate way (see “Disposal”).
- Calibrate the scale, see “Setting up the scale”.

## Switching on voice playback

The measured weight can also be voiced at the end of the weighing process. You can choose among the following languages:

German (DE), English (EN), French (FR), Spanish (ES), and Italian (IT)

1. To switch on the scale, step on the middle of the weighing surface.
2. To switch on voice playback, turn the small dial beyond the resistance point (“CLICK”).
3. Select the desired volume level using the scale .
4. Tap the small button below the language abbreviations to select the desired language. Note the display.

## Choosing the unit of measurement

The unit of measurement is permanently linked to the set language.

Language	Indicator	Unit of measurement
DE, FR, ES, IT	kg	Kilogram
EN	lb	Pound

## Setting up the scale

1. Place the scale on an even, firm and dry surface. Do not place the scale on a carpet.
2. First, briefly step onto the middle of the weighing surface to calibrate the scale. As soon as the display shows "0.0", calibration is complete.  
This calibration must be performed after each battery change.

## Weighing

---

### **WARNING** against injuries

- Do not step onto the scale with wet feet or when its surface is damp.  
Risk of slipping!
  - Position yourself on the scale. Stand still and do not hold on to anything.  
The measured weight will be displayed shortly thereafter.  
A few seconds after you step off the weighing surface, the display will return to "0.0" and you can weigh yourself again. Otherwise, the scale will switch off automatically shortly afterwards.
- Weigh yourself with as little clothing on as possible.
  - Try to always weigh yourself at the same time and under consistent conditions in order to obtain comparable values.

## Resetting the scale

---

Reset the scale to its default settings if unusual or obviously incorrect values are shown on the display.

1. Remove the batteries from the battery compartment and wait a moment.
2. Reinsert the batteries.  
Pay attention to the correct polarity of the batteries (+/-).

After switching the scale on, you can weigh yourself as usual.

## Cleaning

---

### **NOTE** on material damage

- Do not immerse the scale in water or other liquids.
  - Never use harsh, scouring or abrasive cleaning agents or the like to clean the scale. Otherwise, the surface could be scratched.
1. Wipe the scale as needed with a slightly damp soft cloth.
  2. Disinfect the weighing surface from time to time using a suitable disinfecting agent (available at specialist retailers and pharmacies).

## Fault/Remedy

Fault	Possible cause
No function	<ul style="list-style-type: none"> <li>▪ Are the batteries empty or inserted the wrong way around?</li> </ul>
The scale appears to be displaying wrong weighing results	<ul style="list-style-type: none"> <li>▪ Is the scale located on a firm, even surface?</li> <li>▪ Is a device that can cause electromagnetic disturbances located in the vicinity of the scale?</li> <li>▪ Remove the batteries from the scale and wait a moment. Reinsert the batteries.</li> </ul>
"Err" is shown on the display	<ul style="list-style-type: none"> <li>▪ The maximum capacity of the scale of 250 kg has been exceeded.</li> </ul>
"Lo" is shown on the display	<ul style="list-style-type: none"> <li>▪ The batteries are empty (see "Inserting/ changing batteries").</li> </ul>

## Technical data

Model:	2011 "Gina" digital bathroom scale
Battery:	3x LR03 (AAA), 1.5 V ===
Rated current:	< 150 mA
Measuring range:	max. 250 kg
Variance:	1% +/- 100 g
Display scaling:	100 g
Scale dimensions:	approx. 310 x 310 x 28 mm
Distributor (no service address):	GREND S GmbH Stahlwiete 23 22761 Hamburg, Germany

Our products are constantly being developed and improved.  
For this reason, design and technical modifications are possible at any time.

# Declaration of Conformity

---

GREND'S GmbH hereby declares that the BE 2011 digital bathroom scale complies with the 2014/30/EU and 2011/65/EU directives. This declaration will lose its validity if any modification is made to the scale without our consent.

Hamburg, March 2022  
GREND'S GmbH

## Warranty

---

GREND'S GmbH warrants for 2 years from the purchase date the cost-free rectification of defects in material and workmanship through repair or exchange.

In the event of a warranty claim, please return the product with the proof of purchase (specifying the reason for claim) to your dealer.

## Disposal

---

### Disposing of the packaging

Dispose of the packaging with similar materials. Dispose of cardboard with waste paper and of films with recyclable materials.

### Disposing of the product

Dispose of the product in accordance with the regulations that apply in your country.

#### **Devices must not be disposed of in normal household waste.**

At the end of its life, the product must be disposed of in an appropriate way. In this way, valuable materials contained in the device will be recycled and the burden on the environment avoided. Bring the old device to a collection point for electronic waste or to a recycling depot. For more information, contact your local waste disposal firm or your local administrative body.

#### **Batteries and rechargeable batteries should not be disposed of in the household waste.**

All batteries and rechargeable batteries must be brought to a local or district collection point or to the retailer. Such batteries and rechargeable batteries can thus be disposed of in an environmentally friendly way.

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GREENDS GmbH