




⌚ 2 min

⌚ 2 min
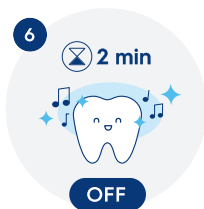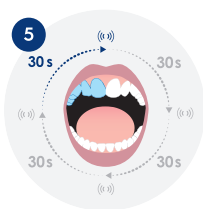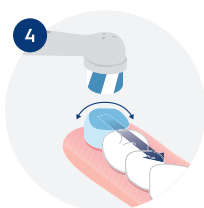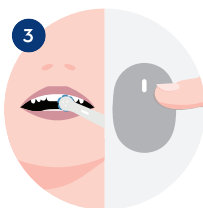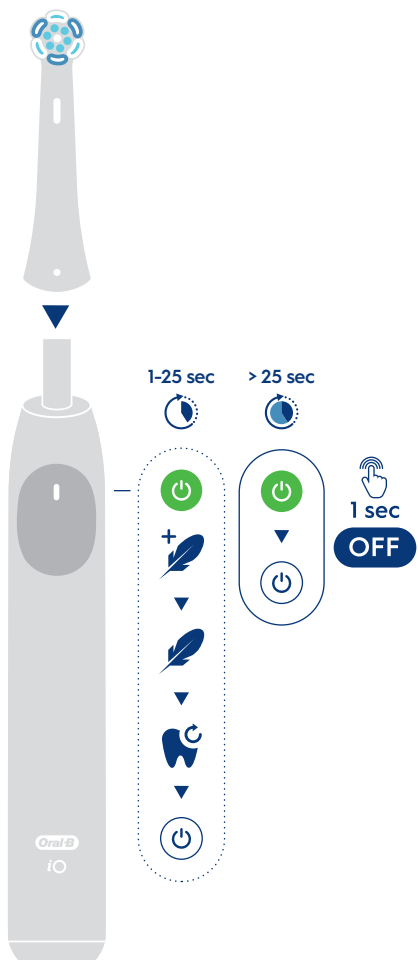
## MODES

MODE 1


MODE 2

MODE 3

## ★ MODE

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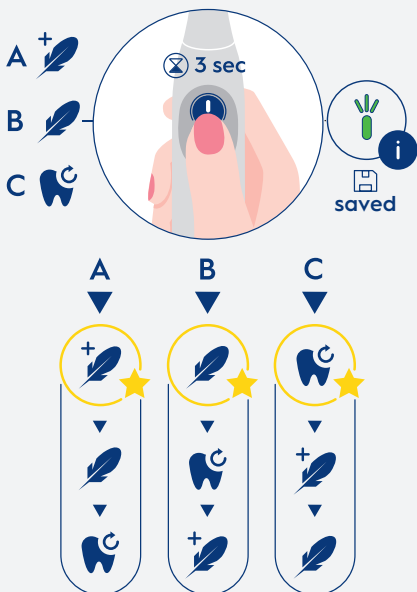## **Welcome to Oral-B!**

Before operating your Oral-B iO toothbrush, please read these instructions and save this manual for future reference.

Rechargeable electric toothbrush.

## **IMPORTANT**

- Toothbrushes can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge, if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children shall not play with the appliance.
- Cleaning and user maintenance shall not be made by children without supervision.
- If the product/cord is damaged, take it to an Oral-B Service Centre.
- Only use charger provided with your appliance.

## **Warning**

- Do not place the charger, smart plug or charging travel case in water or liquid or store where it can fall or be pulled into a tub or sink. Do not reach for it when fallen into water. Unplug immediately.
- This appliance contains batteries that are non-replaceable. Do not open and disassemble device. For battery recycling please dispose whole device according to local environmental regulations. Opening the handle will destroy the appliance and invalidate the warranty.
- Usage by children under age 3 is not recommended.

- Do not chew or bite on brush. Do not use as a toy or without adult supervision.
- This toothbrush is a personal care device and is not intended for use on multiple patients in a dental practice or institution.
- For hygienic reasons we recommend the use of individual brush handles and brush heads for each family member.
- Small parts may come off, keep out of reach of children.
- Clean brush head properly after each use (see «Cleaning Guidelines» paragraph). Proper cleaning ensures the safe use and functional life of the toothbrush.
- When unplugging, always pull the power plug instead of the cord. Do not touch the power plug with wet hands. This can cause an electric shock.
- For tooth brushing only.

To avoid brush head breakage which may generate a small part choking hazard or damage teeth:

- Prior to each use, ensure the brush head fits well. Discontinue use of the toothbrush, if brush head no longer fits properly. Never use without a brush head.
- If the toothbrush handle is dropped, the brush head should be replaced before the next use even if no damage is visible.
- Replace brush (head) every 3 months or sooner if brush head becomes worn.
- Periodically check the entire product/ cord/ accessories for damage. A damaged or non-functioning unit should no longer be used.
- Do not modify or repair the product. This may cause fire, electric shock or injury.
- Never insert any object into any opening of the appliance.
- Avoid sharp bending, jamming or squeezing of the cord.
- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.
- Use this product only for its intended use as described in this manual. Do not use attachments which are not recommended by the manufacturer.


- If stickers are provided, follow instructions on the sticker sheet. Ensure valve on back of the handle is not covered.
- This toothbrush has a waterproof handle, is electrically safe and designed for use in the bathroom.
- Periodically check the operation of the Smart Pressure Sensor by pressing moderately on the brush head during use.

## Important Information

- This toothbrush and pouch (accessory) contain NdFeB magnets and generate magnetic fields. If you use a medical device or have an implanted medical device, consult with the medical device manufacturer or a physician, if you have questions about safely using your Oral-B iO toothbrush.

## Specifications

For voltage specifications refer to the bottom of the charging unit.

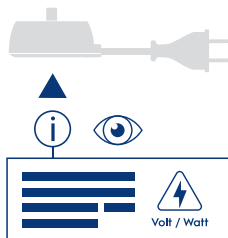
Noise level:  $\leq 66$  dB (A)

## Charging and Operating

To charge the toothbrush battery place the toothbrush on the plugged-in charger.

During charging the on/off button is blinking green. Once it is fully charged the indicator light turns off.


A full charge typically takes 24 hours.

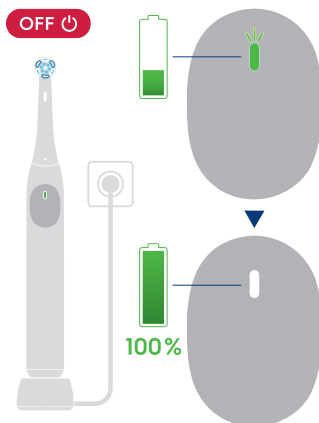English

If the battery is running low the low charge indicator blinks in red after switching off.  
Once the battery is empty, the motor stops.  
To charge the battery please put the toothbrush on the charger.

**NOTE:** In case the battery is empty, charge for at least 30 minutes.

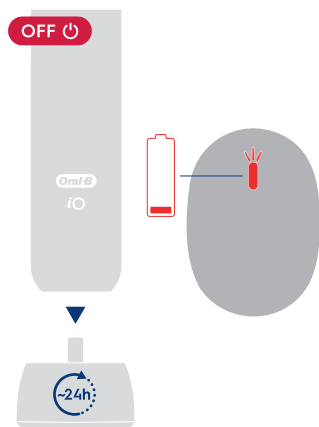You can always store the toothbrush on the plugged-in charger to maintain it at full charge; a battery overcharge is prevented by the device. However, for environmental purposes we recommend to keep charger in the power socket during charging only.

The toothbrush switches off when exceeding the outside operating limits and/or standard battery current. To enable operation the toothbrush needs to be put on the charger.

**Note:** Store the toothbrush at room temperature for optimal battery maintenance.

**CAUTION:** This toothbrush contains a Li-Ion battery. Do not expose the toothbrush to temperatures higher than 40 °C / 104 °F. Recommended ambient temperature for charging is 5°C to 35°C (41°F to 95°F).

## Using your iO toothbrush

### Start brushing


1. Place the brush head on the toothbrush handle. A small gap should remain between handle and brush head.

2. Wet brush head and apply toothpaste.  
To avoid splashing, guide the brush head to your teeth before switching the toothbrush on by pressing the on/off button.
3. Guide the brush head from tooth to tooth.
4. Brush the outsides, the insides and the chewing surfaces. Brush all four quadrants of your mouth equally. You may also consult your dentist or dental hygienist about the right technique for you.


In the first days of using any electric toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or dental hygienist.

If you have sensitive teeth and/or gums, Oral-B recommends using the «Sensitive» or «Super Sensitive» modes (optional in combination with an Oral-B «Gentle Care» brush head). Availability of modes may vary depending on model purchased.
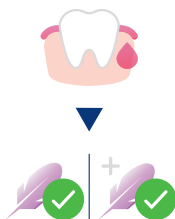
### **Toothbrush Timer (Stutter timer)**


The iO toothbrush is equipped with timer features, which help you to achieve the professionally recommended 2-minute brushing time. The timer signals with a short stuttering sound at 30-seconds intervals to move to the next quadrant of your mouth. A short melody indicates the end of the professionally recommended 2-minute brushing time. Once reached the professionally recommended 2-minute brushing time, the iO toothbrush switches off automatically.

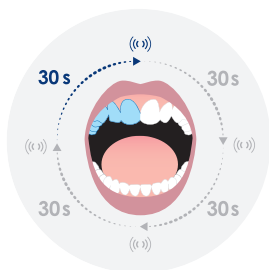The timer memorises the reached brushing time, even when the handle is briefly switched off during brushing. The timer resets if pausing longer than 30 seconds or placing the handle on the plugged-in charger.

### **Pressure Sensor**

The iO toothbrush is equipped with a pressure Sensor. If too much pressure is applied, the toothbrush slows down the oscillating rotation of the brush head.
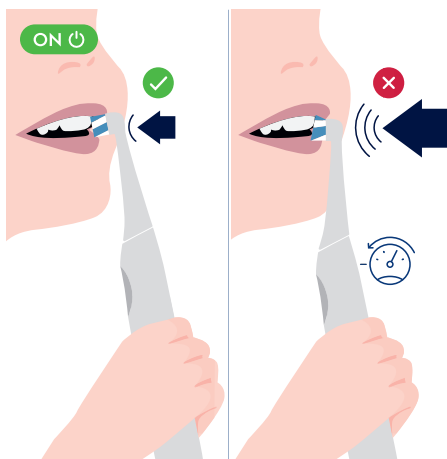
## Brushing modes (depending on model)

Your toothbrush offers different brushing modes:

«**Super Sensitive**» – Super gentle cleaning for extra sensitive areas (LED blinks in green once)

«**Sensitive**» – Gentle, yet thorough cleaning for sensitive areas (LED blinks in green twice)



«**Daily Clean**» – Standard mode for daily cleaning (LED blinks in green three times)

English

ON 


|                                                                                     |                                                                                         |
|-------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------|
|    | 1 x    |
|   | 2 x   |
|  | 3 x  |

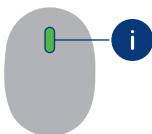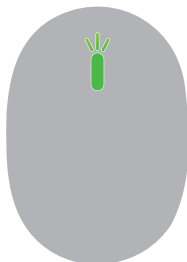**Note:** The sequence and availability of brushing modes may vary depending on model purchased.


### Mode setting:

1. The toothbrush automatically starts in the «Super Sensitive» mode.
2. To switch to other modes, successively press the on/off button in the first 25 seconds during use. After 25 seconds the toothbrush switches off by pressing the on/off button.
3. To switch off the toothbrush at anytime, press and hold the on/off button (approximately 1 sec) until the motor stops.
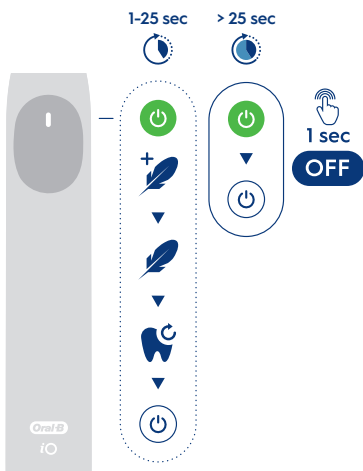
### Favorite mode setting:

1. You might choose your favorite mode as first mode.
2. To setup your favorite mode, switch on the toothbrush by pressing the on/off button and select your favorite mode.
3. Press and hold the on/off button for 3 seconds. The motor will stop (approximately after 1 sec -> continue pressing) and the LED blinks green confirming the successful setup of your favorite mode.
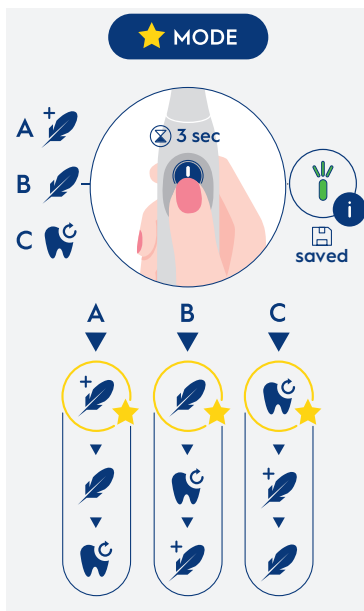
## Brush heads

Oral-B offers you a variety of different Oral-B iO brush heads that only fit your Oral-B iO toothbrush handle. These unique brush heads are dentist-inspired to leave your mouth feeling clean and refreshed.

Brush head «**iO Ultimate Clean**» - for 360° Advanced Deep Cleaning

- Thousands of bristles angled at 16° for superior plaque removal
- Spiraled bristle pattern for better adaptability
- Dual-length bristles for cleaning even in hard-to-reach areas

Brush head «**iO Gentle Care**» - for Healthier Gums\* as of Day 1

- 4000+ soft bristles for effective cleaning
- Curved bristle pattern for gentle gliding effect
- Ultra-thin adaptable bristles – even for sensitive teeth & gums

\*vs. a regular manual toothbrush

Brush head «**iO Radiant White**» - for Whiter Teeth\* as of Day 1

- Thousands of densely arranged bristles for thorough cleaning
- Curved & angled outer bristles for better reach in-between teeth
- 5 polishing petals to remove everyday stains (like the ones from coffee, wine, or smoking)

\*by removing surface stains


Brush head «**iO Interdental Clean**» - for Extra Clean in Hard-To-Reach-Areas

- Dual-length high-density bristles clean in between teeth
- Outer bristles angle in to support central bristles
- Use as Step 2 after your usual brushing routine – even implants & braces

**Note:** iO Brush heads are available on [www.shop.oralb.de](http://www.shop.oralb.de) and cooperating stores.


### When to exchange brush head

Most iO Oral-B brush heads feature INDICATOR® bristles to help you monitor brush head replacement need. With thorough brushing, twice a day for two minutes, the bristle colour will fade half-way approximately within 3 months, indicating the need to exchange your brush head. If the bristles splay before the colour recedes, you may exert too much pressure on teeth and gums. In general, we recommend replacing the brush head every 3 months.
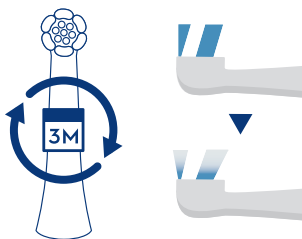
## Cleaning Guidelines

1. After each brushing, rinse your brush head under running water with the handle switched on.
2. Remove the brush head and clean the handle and brush head separately. Make sure water runs inside the top of the handle while the handle is switched on, as well as inside the brush head.
3. Wipe all parts dry and let them air-dry before reassembling your toothbrush again.

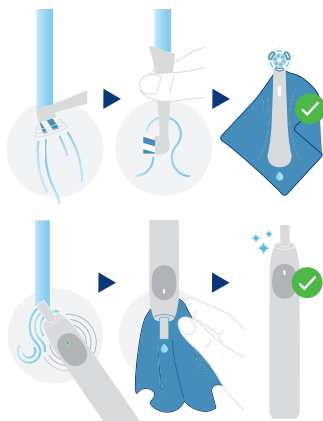4. Unplug the charger before cleaning.
5. Clean the charger regularly with a damp cloth.

6. Please check individual cleaning instructions for each accessory in its «Accessories» paragraph.

## Accessories (depending on model)

### Travel Case

Store the clean and dry iO toothbrush/ brush heads in the travel case. The outer and inner case should be cleaned with a damp cloth only.


**Note:** Clean stains on travel case surfaces immediately. Store travel case in a clean and dry place.

### Brush head holder stand / Brush head holder with cap

You may store your cleaned and dry brush heads on the pins of the brush head holders. Make sure your brush heads are dried before placing them in the holder and/or closing the cap. Brush head holders and cap are dishwasher-safe.

### Zipper Pouch / Magnetic Pouch

You may store your cleaned brush heads and toothbrush in the pouches. Make sure your brush heads and toothbrush are dried before placing them in the pouches. The outer and inner surfaces can be cleaned under running water. Store pouches in a clean and dry place.

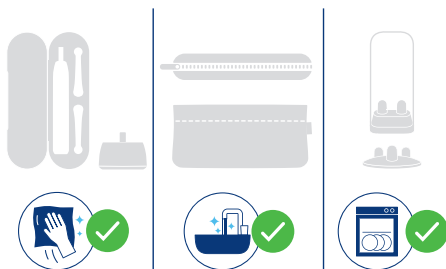English

**Note:** Clean stains on pouch surfaces immediately.

Subject to change without notice.


## Environmental Notice


Product contains batteries and/or recyclable electric waste. For environment protection do not dispose in household waste, but for recycling take to electric waste collection points provided in your country.

## Battery removal at the end of the product's useful life

Do not disassemble the product except when disposing of the battery. Run down battery before removing. When taking out the battery for disposal, use protective clothing (gloves, safety glasses, etc.) to omit skin contact whilst leaking batteries. Use caution not to short the positive (+) and negative (-) terminals. Insulate terminals (e.g. with tape) and put the battery in a plastic bag for disposal.

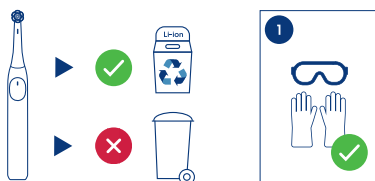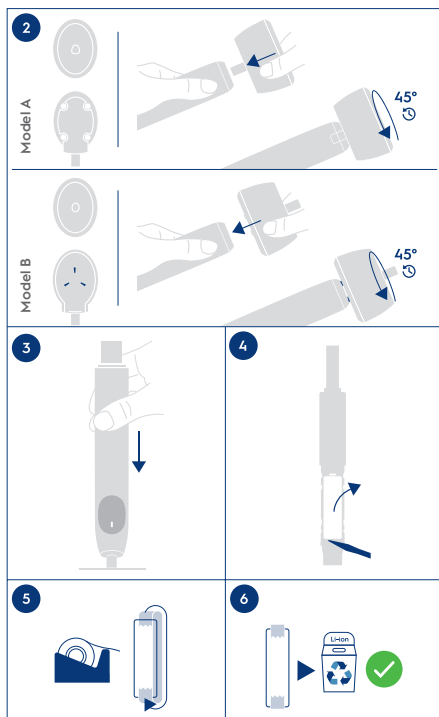Open the handle as shown, remove the battery and dispose of it according to local environmental regulations.

**CAUTION! Opening the handle will destroy the appliance and invalidate the warranty.**

## Guarantee

We grant 2 years guarantee on the product commencing on the date of purchase.

Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship,

either by repairing or replacing the complete appliance as we may choose. Depending on availability replacement of a unit might result in a different colour or an equivalent model.

This guarantee extends to every country where this appliance is supplied by Braun / Oral-B or its appointed distributor. This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brush heads, as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorized persons and if original Braun / Oral-B parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorized Oral-B Braun Customer Service Centre.

### **Warranty replacement brush heads**

Oral-B warranty will be voided if electric rechargeable toothbrush handle damage is found to be attributed to the use of non-Oral-B replacement brush heads. Oral-B does not recommend the use of non-Oral-B replacement brush heads.

- Oral-B has no control over the quality of non-Oral-B replacement brush heads. Therefore, we cannot ensure the cleaning performance of non-Oral-B replacement brush heads, as communicated with the electric rechargeable handle at time of initial purchase.
- Oral-B cannot ensure a good fit of non-Oral-B replacement brush heads.
- Oral-B cannot predict the long-term effect of non-Oral-B replacement brush heads on handle wear.

All Oral-B replacement brush heads carry the Oral-B logo and meet Oral-B's high quality standards. Oral-B does not sell replacement brush heads or handle parts under any other brand name.

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Braun GmbH