# VELFORM® DERMAPLANING SYSTEM

## English

**Velform® Dermaplaning System**  
Thank you for purchasing the Velform® Dermaplaning System, the lighted facial exfoliator that exfoliates, removes hair and smooths the skin. Please read this user manual carefully and pay close attention to all the safety instructions before using the product for the first time. Keep this manual on hand for future use.

1. Safety Warnings
2. Product Components
3. Basic Directions
4. Maintenance

### 1. Safety Warnings

**DO NOT USE Velform® Dermaplaning System** around the lips or nose, or directly on the eyelids. Velform® Dermaplaning System is for use ONLY on the forehead, cheeks, upper lip area and chin, and above the upper brow. **DO NOT USE** on skin showing dryness, acne, sunburn, psoriasis or eczema, or skin that has just received medical treatment. **DO NOT USE** if you have diabetes or circulatory or bleeding disorders. Discontinue use if skin irritation occurs. **DO NOT USE** on moles, sores or lesions. **DO NOT** rinse the handle or product under water, or immerse it in liquid. Keep the Velform® Dermaplaning System dry. Do not store in wet or damp places or near containers of water that it could fall into.

**DO NOT** dispose of batteries in fire, as they may explode. Dispose of old batteries in accordance with all applicable local regulations. Keep your Velform® Dermaplaning System out of reach of children. Not recommended for children under age 12.

- ### 2. Product Components
1. Head
  2. Top of the handle
  3. Head ejection button
  4. ON/OFF button
  5. Handle
  6. Tray with six replacement heads
  7. Battery cap

### Getting Started BEFORE FIRST USE

- Test on a small, discreet area.
- Discontinue use immediately if sensitivity or an allergic reaction occurs.

The Velform® Dermaplaning System comes ready to use right out of the package.

Recommended for exfoliation/hair removal on the forehead, cheeks, upper lip area and chin, and above the upper brow.

### BATTERY INSTALLATION / REPLACEMENT

- Gently remove the battery cap to reveal the battery compartment. (Picture 8)
- Insert one AAA battery into the compartment. Make sure the positive end (+) is facing inside the compartment. (Picture 9)
- Both the product and the cap have a mark. Make sure both marks are lined up (Picture 10).
- To secure the cover on the handle, first line up both marks and then slightly turn the cover clockwise to ensure it is firmly closed. (Picture 11).

### Instructions for Use DRY USE ONLY.

We recommend using the Velform® Dermaplaning System at night, before bed, to give the skin time to recover.

Be sure to space treatments one to two weeks apart.

Select a new, clean Velform® Dermaplaning System head for each use.

- With one hand holding the handle, gently slide the safety guard off with the other hand (Picture 2).
- Press the ON button to illuminate the area (Picture 3).
- For best results, hold the skin taut with your free hand and the unit at a 45° angle with the other.
- Use short, feathery strokes and gently skim across skin, beginning at the hairline and moving inward following the contours of the face (Picture 4).

- Do NOT press the head down into the skin. No pressure should be used. Pass over the same area only ONCE during each use.
- The Velform® Dermaplaning System should be used by following the recommended pattern (Picture 5).

### PROFESSIONAL TIPS

After using the Velform® Dermaplaning System, let your skin rest.

Do not use creams or make-up for several hours after use.

- ### AFTER EACH USE
- Let the skin breathe for several hours before applying lotions, serums or make-up.
  - Applying make-up immediately after use is NOT recommended.
  - Do not sunbathe after using the product.

Note: After use, the skin may be slightly pink or red. This is normal after treatment and should fade within a few hours.

If redness or burning increases over time, discontinue use of the product.

### Replacement Heads

Use a new, clean Velform® Dermaplaning System head for each use.

- Make sure the unit light is off.
- Remove the head by sliding the ejection button down (Picture 6).
- Discard the used head.
- Hold the tray securely with one hand.
- Grasp the handle with the other.
- Hold the top of the handle over the new head (Picture 7A).
- Angle the top of the handle down towards the top of the new head (Picture 7B).

- Press the top of the handle down into the top of the new head (Picture 7C).
- Rock the handle backward while pressing down until you hear the new head click into place. Now you can lift the handle (Picture 7D).
- The unit is now ready for use (Picture 7E).

NOTE: DO NOT USE THE PRODUCT if the replacement head appears to be damaged.

### Care and Maintenance

Keep the Velform® Dermaplaning System safety guard on when handling the unit. Remove the safety guard when you are ready to use the unit. For best results, discard the head after each use. Make sure the unit is turned off after use. Never immerse it in water or any other liquid. Keep the product in a safe, dry place.

### Quality Warranty







This product is covered by a warranty against manufacturing defects subject to the time periods stipulated by the legislation in force in each country. This warranty does not cover damage resulting from inadequate use, negligent commercial use, abnormal wear and tear, accidents or improper handling.

Australia & New Zealand only  
This product has the benefit of certain Consumer Guarantees. These are prescribed by the Australian Consumer Law & the New Zealand Consumer Guarantees Act 1993 both of which provide protection for consumers. There is no express warranty for this product in Australia or New Zealand. The above paragraph refers to other countries.

### CAUTION DISPOSAL OF MATERIALS

A symbol of a crossed-out wheeled bin means you should find out about and follow local regulations about disposing of this kind of product.

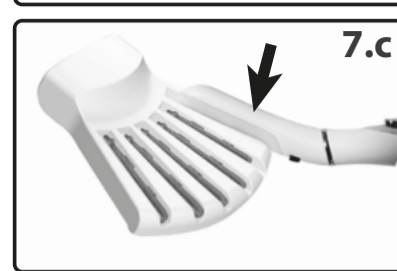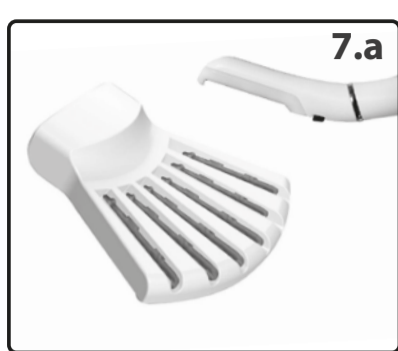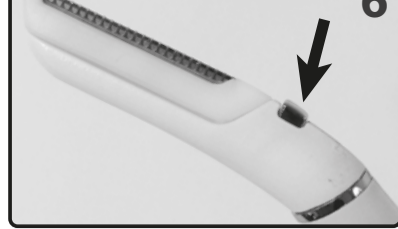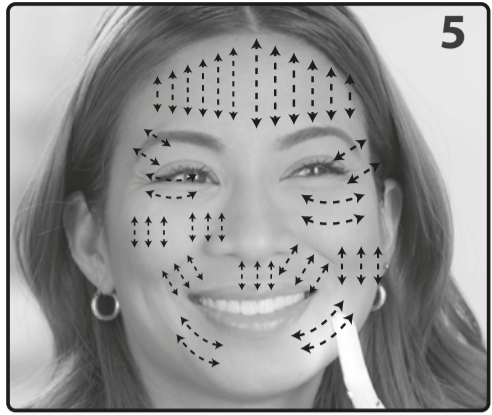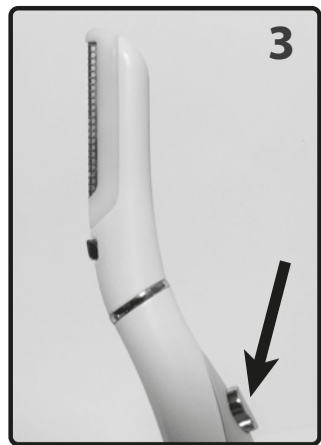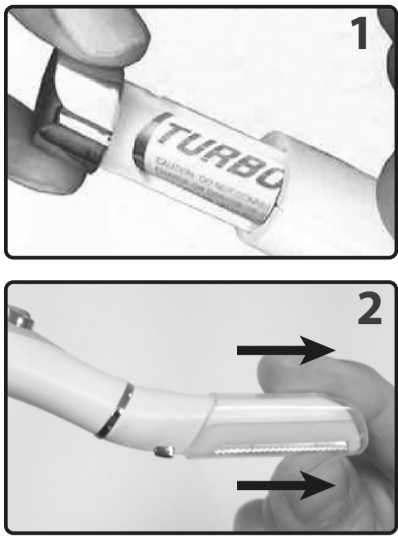Do not dispose of this product as you would other household waste. Dispose of this product in accordance with the corresponding local regulations. Electrical and electronic appliances contain hazardous substances that can have harmful effects on the environment and/or human health and should be recycled properly.

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