Sport
MG**2**

breo

CONTENTS

Caution	06
Product Configuration	10
Product Characteristics	11
Product Operation	12
Troubleshooting	17
Disposal Instructions	18
Protection	20
Product Information	22

EN

EN

C a u t i o n

Precautions

1. Please read the instructions carefully.
2. To avoid pain caused by excessive massage, do not use this device for a long time.
3. Do not use this device while driving.
4. Do not use this device on the head or other privates. Do not use the product on surgical area. Do not use the product on areas of bruises, damaged skin, swelling, or inflammation.
5. People with cardiovascular (or once had), acute, or dermatologic diseases, malignant tumors, or other abnormal symptoms, including hematologic disorders with hemorrhagic symptoms, osteoporosis, soft tissue diseases, and strokes, should not use this device.
6. In case of restricted physical capabilities, physical or mental problems, or mental disabilities, do not use the device without medical consultation first.
7. Do not use the device during menstruation or pregnancy. Consult the medical professional before using

EN the device.

8. Do not use it for entertainment purposes.
9. This device is not intended for use by children under the age of 8. Children aged between 8 and 12 must be supervised while using the device.
10. Do not scratch or stab the device with sharp objects.

Cautions For Use

1. This device is a household electronic health care device, not as a medical device, and is not intended for diagnosis, treatment and other commercial purposes.
2. Please remove any accessories you might be wearing or other objects that might affect the massaging process before use.
3. Stop at the first sign of pain or discomfort of the massage area.
4. If the device works abnormally, please shut it down immediately, cut off the power supply, and contact the customer service team. Do not disassemble or repair it by

EN yourself.

Product Protection

1. Only use the original charger and the adapter conforming to the standard to charge, and pay attention to the electric safety.
2. Please keep this device away from water sources, fire sources and corrosive environments.
3. This device contains a non-replaceable battery. Do not disassemble the battery.

***Frequency**

1. Do not operate continuously for more than 10 minutes each time.
2. Overuse might cause discomfort or neck and skin injuries.
3. When you use the massager for the first time, be mindful of how you feel and observe your skin and body; do not massage for too long.

EN

*** Before Use**

1. Before using this device, please completely relax your body.
2. The skin may become a little red after massaging. This phenomenon is normal and will quickly disappear.


***After Use**



The massage area may turn red for a short period of time, but will return to its original state after a while.


***Product Cleaning**

Clean the surface with a dry cloth dipped in water or 75% alcohol. Do not immerse the device in water.

Product Configuration

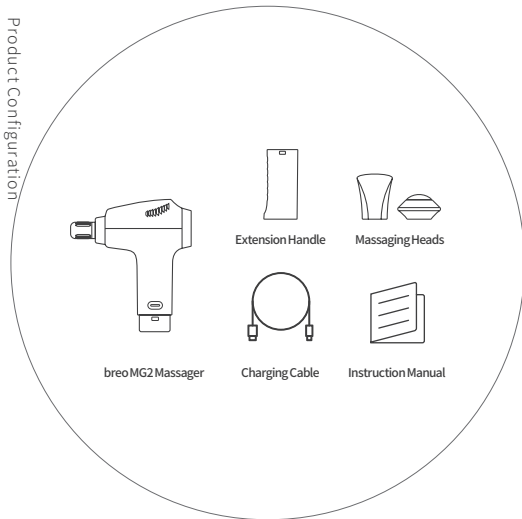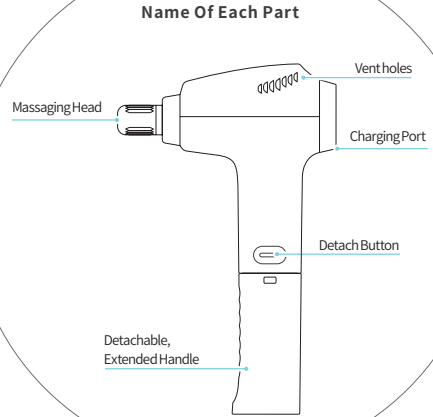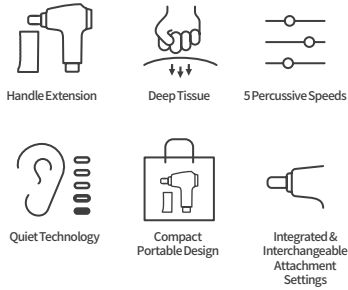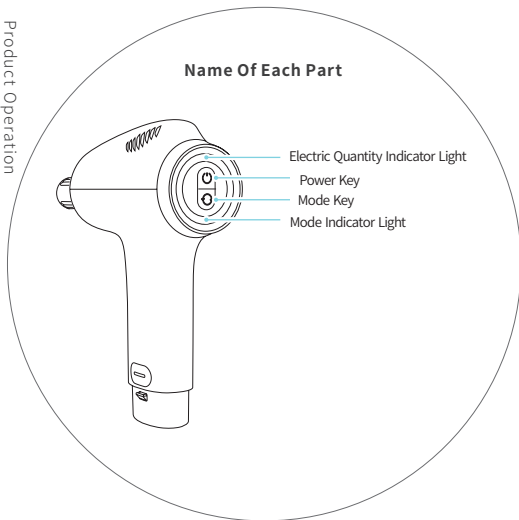**EN****Proper Usage**

Preparation before use

Please check the device and accessories for completeness and damage.

1. Turn on

Firmly press and hold the power button for 2 seconds to turn on the device.

2. Mode selection

The default is the first gear speed. Press the mode key to switch among the speeds.

To let the device enter standby status, press the power button during operation. Press the power button again to resume operation.

3. Turn off

The device will automatically turn off after 10 minutes of use.

To turn off the device manually, firmly press and hold the power button for 2 seconds until the device shuts down.

EN 4. The skin may become a little red after massaging. This phenomenon is normal and will quickly disappear.

Charging Method:

1. This device cannot be used when charging.
2. Connect the charging plug to the device's charging port.
3. The self-protection program will be activated and stop working automatically if the device is connected to the power supply for charging when the device is working.
4. A red light indicates the battery power is less than 20%.
5. A white light indicates the battery is fully charged.

* It is normal for the device and the power adapter to heat up a little while charging.

Charging Time:

The charging time is about 3-4 hours.

End Of Charge:

When the battery is fully charged, unplug the power

EN adapter.

Troubleshooting

In case of any failures during operation, please perform the following checks. If the fault cannot be eliminated or the device still cannot operate satisfactorily, please cut off the power supply, and professional technicians shall be consulted for maintenance. Do not disassemble the device by yourself!

Symptom	Possible Cause	Basic Solution
Power-on failure	Check whether the battery level is too low	If yes, please charge timely
Weak massage intensity	Check whether the battery level is too low	If yes, please charge timely

If the problem cannot be solved after troubleshooting, please send an email to support@breocare.com or call us at +1(866)-546-3587 (USA). We are always happy to help.

Disposal Instructions

- The battery must be removed from this device in accordance with regulations before it is discarded;
- When removing the battery, the device must be powered off;
- The battery should be disposed of safely in accordance with laws and regulations.

Do not disassemble the device if you are non-professional.

Operating steps for battery removal:

1. Cut off the external power supply and remove the device shell;
2. Unplug the power supply on the battery;
3. Remove the battery;
4. Please put the waste battery to the fixed recycling point;

EN 5. Please abide by local environmental protection laws and discard guidelines and take good care of the environment.

EN

Protection

Storage

When the device is stored for more than 3 months, check the device and recharge it before using. Do not store the device closely against the ground, wall or ceiling. Store in a well-ventilated place and avoid direct sunlight and harmful erosion.

Transport

The device is packed with simple anti-vibration measures for transportation by air, railway, highway and sea. Avoid rain, snow and collision shall be avoided and do not leave it upside down.

Power Supply Check

Check the power supply. Do not start the machine if the power voltage exceeds the specified range.

Maintenance

Take measures against moisture, rodents and pollution.

EN Turn the device off before cleaning. Store the device in a safe place after use. Do not place heavy items on the device. If any fault occurs, call professional and technical staff for repair. Do not disassemble the device if you are non-professional.

Disposal And Recycling

This device should not be discarded with household waste and should be disposed of safely in accordance with local environmental laws and guidelines.

EN

Product Information

Product Name	Breo Deep Tissue Massager
Product Model	MG2
Rated Voltage	5V
Rated Power	10W
Battery Capacity	2400mAh
Product Weight	1.72lb
Product Size	L7.28×W2.28×H8.15 (inch)
Extension Dimensions	L7.28×W2.28×H11.69 (inch)
Place Of Production	China

Shenzhen Breo Technology Co., Ltd.

Floor 19, Hisense South Building, No. 1777 Chuangye Road, Nanshan District, Shenzhen, China

Tel: +1(866)-546-3587 (USA)

Customer Service Support: support@breocare.com (Global)

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