**LEBEN
LANG**

USER MANUAL

KÖRPERANALYSE- WAAGE

BODY COMPOSITION SCALE

LBL701

Download on the
App Store

GET IT ON
Google Play

Deutsch	3
English	17
Français	29
Español	41
Italiano	53
Nederlands	65

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Thank you for your trust in LEBENLANG. We strive for the highest standards in quality and design and hope you enjoy your new product.

Supported devices:

- iOS 8.0 or higher
- Android 6.0 or higher

Install the Fitdays App:

1. Search for **“Fitdays”** in the App Store or Google Play Store or scan the QR code next to it:
2. Download and install the app on your device.
3. After installation is complete, the app icon will appear on your phone or tablet.

Scan the QR code to
download the app

GENERAL SAFETY INSTRUCTIONS

When using electrical appliances, basic safety precautions must always be observed, including the following:

- Please read this user manual to ensure effective use. If you pass the product on to someone else, please include this manual.
- LEBENLANG devices comply with relevant safety regulations. Repairs and replacement of the power cord must only be carried out by specialists authorized by the manufacturer. Improper repairs can cause significant risks to the user.
- The scale is intended exclusively for private use and not for commercial purposes.
- Children should not play with the scale. Cleaning and maintenance should not be performed by children.
- The scale can be used by children aged 8 and above, and by individuals with reduced physical, sensory, or mental abilities, or lack of experience/knowledge, provided they are supervised or instructed on the safe use of the scale and understand the risks involved.
- Keep the scale away from heat, direct sunlight, moisture, sharp edges, etc.
- Check the scale regularly for damage. Do not use the device if it is damaged.
- Do not use the scale outdoors.
- Do not attempt to repair the scale yourself; contact an authorized expert.
- If the scale is damaged, it must be replaced by the manufacturer, their service representative, or similarly qualified personnel to avoid hazards.
- Always place the scale on a hard, dry, and flat surface before taking a measurement.

SPECIAL SAFETY INSTRUCTIONS

- DO NOT use with medical implants such as pacemakers.
- DO NOT stand on the edge of the scale or jump on it.
- DO NOT overload the scale (Max. 180kg).
- DO NOT drop the scale or drop objects on it, as this can damage the sensors.
- DO NOT immerse the scale in water or use chemical cleaning agents. It is sufficient to clean the scale with a slightly damp cloth.
- Not recommended for pregnant women.
- Not recommended for infants, toddlers, and children under 10 years old.
- Every measurement made with this scale is for reference purposes only and should not be considered as medical advice.
- Consult your GP or doctor before making any changes to your diet, exercise plan, or physical activities.
- Ensure your feet are dry before stepping on the scale.

CONTENTS & TECHNICAL DETAILS


- 1 x Smart Body Fat Scale (max. 180kg)
- 1 x User Manual
- 1 x Type-C Charging Cable (Universal adapter not included)
- Weight unit: kg/lb/st (units can be changed via APP settings)

BATTERY INFORMATION

- Lithium battery capacity: 500 mAh
- Charging time: 30 minutes
- Charging cycles: 99
- Operating temperature: 10°C - 40°C
- Operating humidity: 20% - 90%



FITDAYS APP FUNCTIONS

Body weight	Body fat	Visceral fat	Protein mass	BMI (Body Mass Index)
Muscle weight	Bone mass	Water percentage	BMR (Basal Metabolic Rate)	Body age
Muscle rate	Ideal weight	Weight without fat	BFR	Protein rate
Subcutaneous fat	Height recording (<16 years)	Body fat measurement by body part	Partial muscle rate	Obesity level

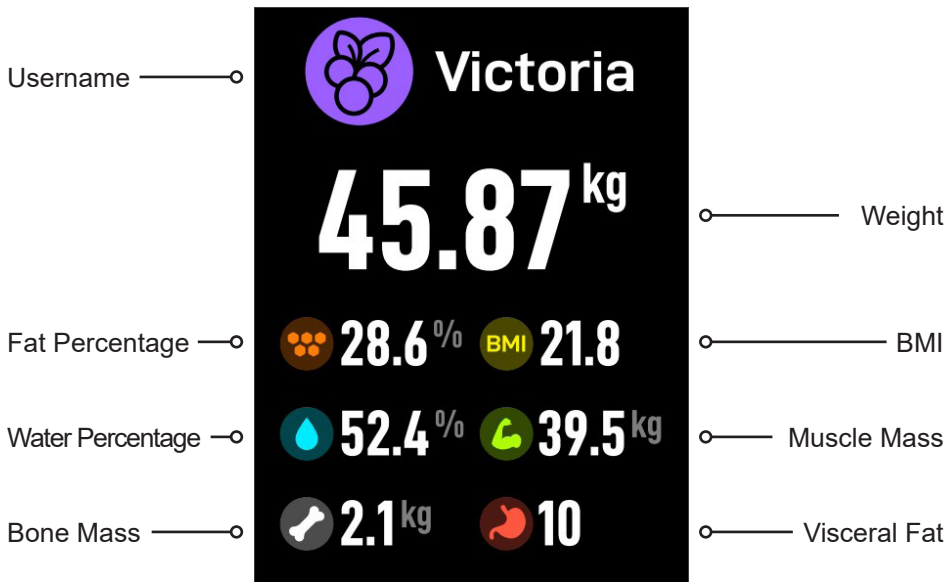APPLICATION / OPERATION

1. Open the Fitdays APP.
2. Press the power button next to the USB charging port.
3. Place the Smart Body Fat Scale on a hard, flat surface & turn on Bluetooth on your smartphone or tablet.
4. Register in the Fitdays app with your own ID, via email or social media account.

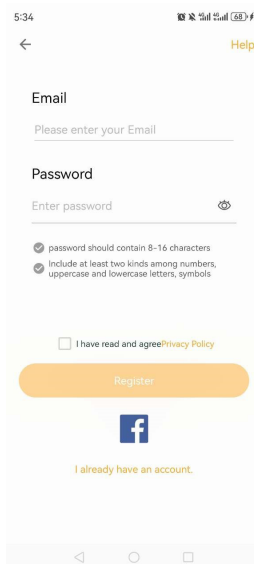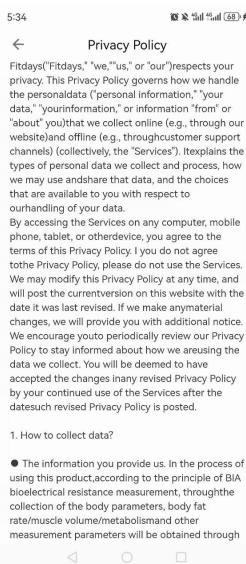5. Enter the required information and confirm with “OK.”
- **Note:** Incorrect information about gender, age, and height can lead to inaccurate measurements.

5:35 personal information

Personal information is only processed as health data. Fitdays will keep your information confidential.

Male Female

Nickname 94109 >

Birthday Jan. 1, 1995 >


Height 165cm >

Target weight 56.1kg >




Mode selection Standard Mode >

For the accuracy of the measurement results, it is recommended that you fill in the real body data


Confirm

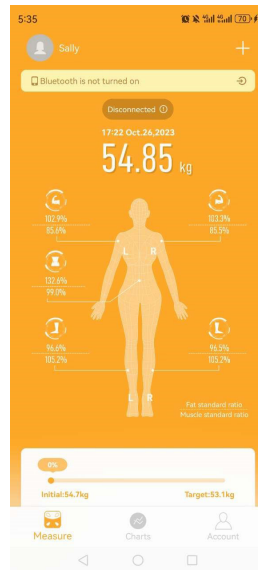
6. Pair the Smart Body Fat Scale via Bluetooth. Click on Account -> Device -> “+ Add Device” -> The pairing search will now start.
- The devices found by the app will now appear.
 - » **To find the scale, it must be switched on!**
 - After successfully connecting to your device, the active connection will be displayed on the main page.
 - If “Disconnected” is displayed, pairing failed, please repeat the pairing steps.

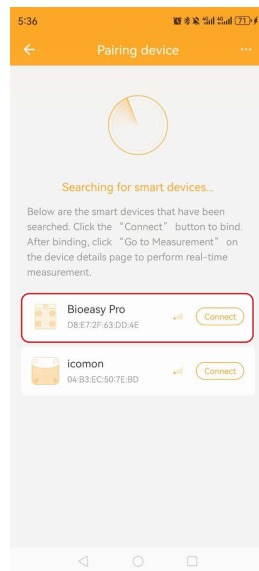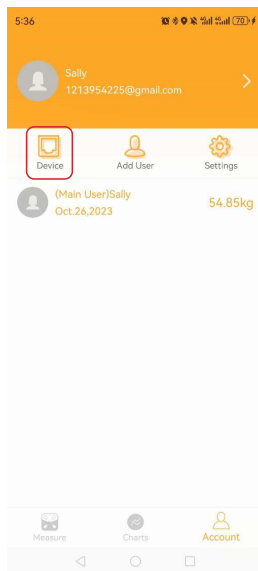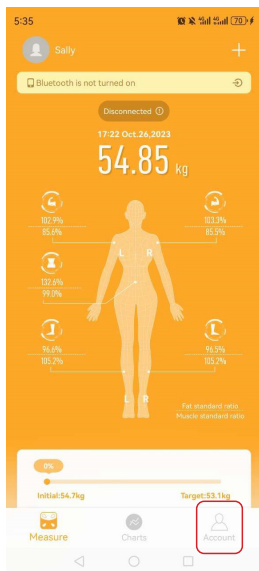7. Tap the Smart Body Fat Scale with your foot to activate the display.
 - Wait until the display shows “0.00” (kg/lb/st).
 - The Bluetooth connection symbol will now appear on the LED display of the scale.
 - The connection status will also be displayed on the main page of the Fitdays app.
8. Stand barefoot on the Smart Body Fat Scale.
 - Make sure your fingers fully touch the four electrodes on the extendable handle (please note to hold the left and right side of the handle in the correct hand) and at the same time stand on the two metal electrode plates with your feet according to the instructions below.

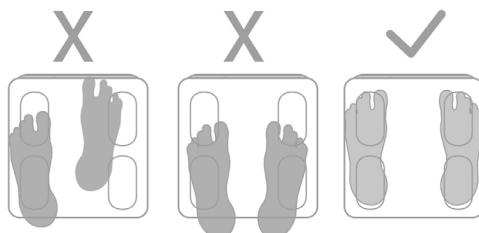
- Pay attention to the handle grip, especially the position of the thumb and four fingers, ensure both make full contact with the metal electrodes and are on the correct side!

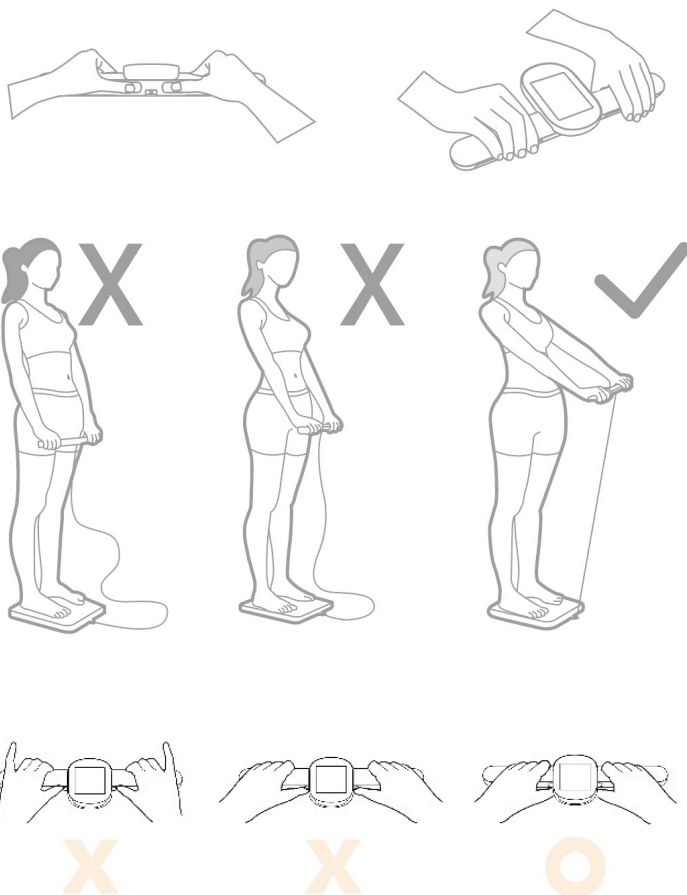
Note:

If you step on the scale without being barefoot or your body posture is incorrect, only the weight and BMI will be measured and displayed on the main page. For full data, stand barefoot on the scale and hold the handle in the correct position as shown above.


9. The number on the display will blink while the Smart Body Fat Scale weighs and analyzes.
 - Please stay still and in the same position for about 30 seconds.
 - The app will show the measurement progress and play music during the body fat measurement.
 - After the measurement is complete, you can check the additional measured data in the app.

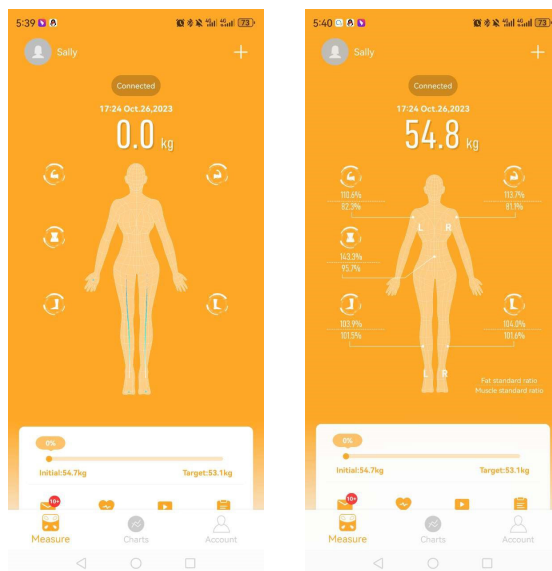
During the body data measurement, the scale displays the following:

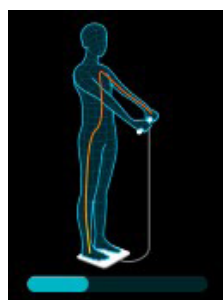
Measurement

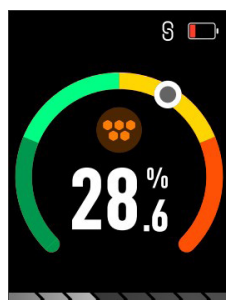
Body Fat Percentage

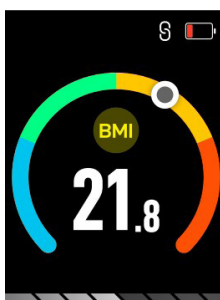
BMI

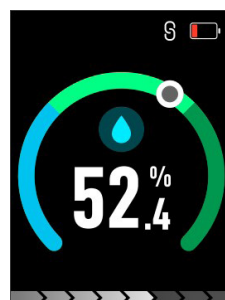
Body Water Percentage

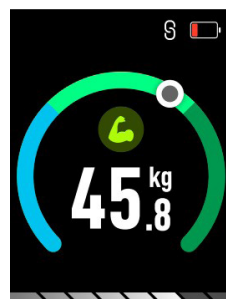Muscle Mass

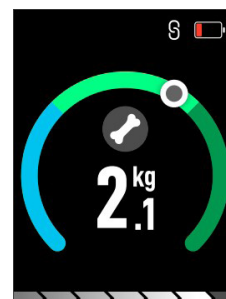Bone Mass

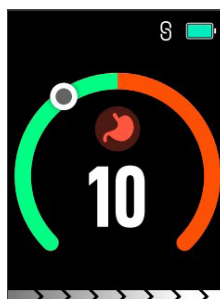Visceral Fat

Summary


After completion, you can see the 7 basic values on the TFT screen of the scale, including body fat percentage, BMI, weight, body water percentage, bone mass, muscle mass, and visceral fat.



The screen also displays the weight difference compared to the last weight measurement.


10. You can click “More Data” to get additional parameter data.
- If you click the “Collapse Data” button at the top right, you will receive a detailed body composition analysis report.

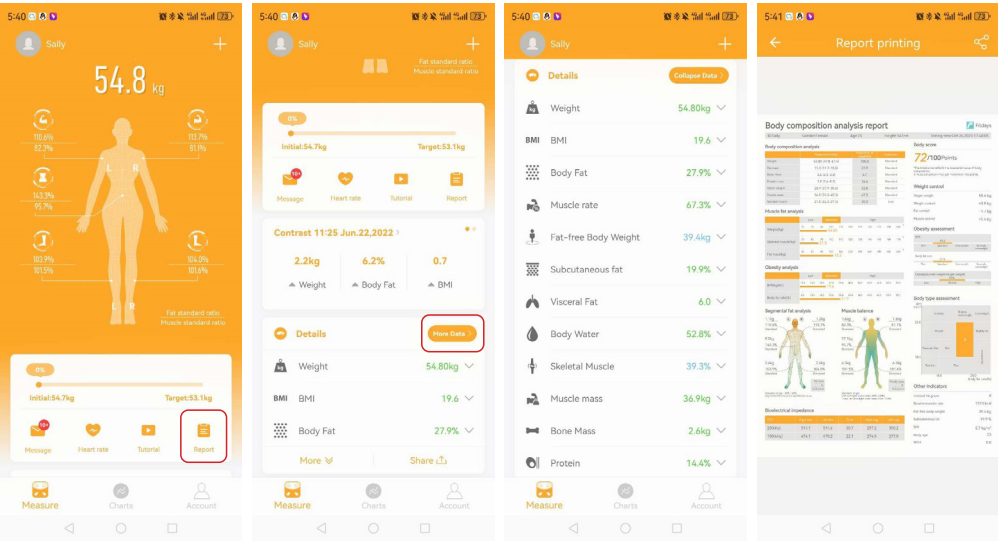11. By clicking “Record,” you will see your long-term progress report.

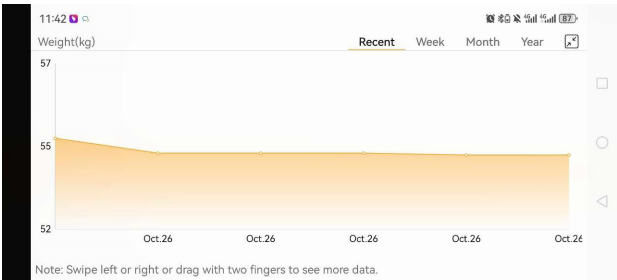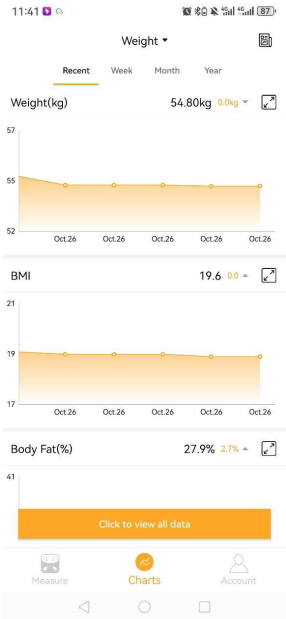12. Add users or switch users by clicking >> ⓘ Add User or go to the user account and select an existing user.
- If multiple users have similar values (e.g., family members with similar weights), the scale may sometimes show options for the user to manually confirm which profile should be used.
 - Up to 24 profiles can be stored on each account!

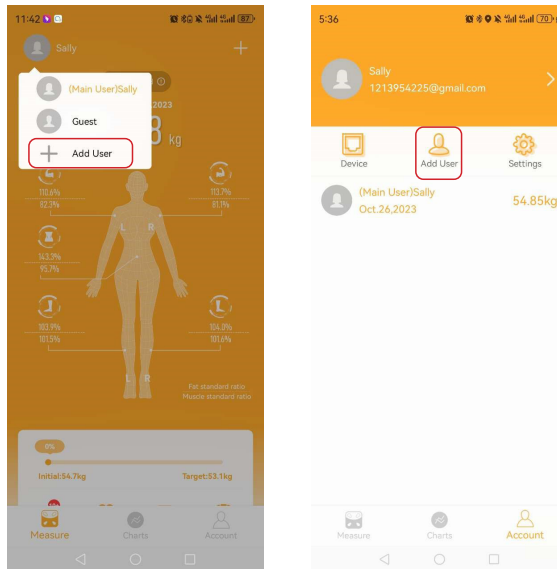Note:

- This Smart Body Fat Scale offers offline functionality. Setting up the scale and connecting it to the app is necessary for the first use to analyze body fat and other parameters. After the initial setup, you can conveniently measure body fat, BMI, and body types without needing to connect the app each time. However, it is important to note that offline measurement will not be activated if there is a weight difference of more than 2kg compared to the user's last measurement.
- In offline mode, all successful measurement results will automatically sync with your app the next time the scale is connected to the app. This ensures your data is seamlessly updated and accessible for further analysis and tracking.

Display symbols of the scale:

LO	Low battery – Please charge via Type-C-USB cable
Err	Overload – The scale shows Err (the app cannot display this symbol)
C	An error occurred during measurement. Step off the scale and then step on again to repeat the measurement process.

TROUBLESHOOTING

1. Bluetooth not connected?

- Ensure both the scale and the smartphone have Bluetooth enabled.
- Check the iOS version of your smartphone, it must be iOS 8.0 or higher/ Android 6.0 or higher.
- Click on Account -> Device -> "+ Add device" -> The search for Bluetooth devices will now begin. The scale must be switched on to be found!

2. No body fat data measured when weighing?

- The user must stand barefoot on the Smart Body Fat Scale.
- Ensure proper grip during measurement.
- Make sure both your feet and the scale are dry.
- The scale will not finish weighing if the number display is flashing. Remain standing on the scale for about 30 seconds until the number stops flashing.
- Bluetooth is not turned on.
- Bluetooth is paired with another scale.

3. The scale provides inconsistent weight readings?

- Ensure the scale is placed on a hard, flat surface.
- Check each sensor foot on the back to ensure nothing is stuck on the bottom.
- Press the power button again.

4. The app shows incorrect fat-muscle ratio?

- The body parameter setting is incorrect. Check if the correct gender, height, and age have been entered.

5. How do I recalibrate the scale after moving or changing batteries?

- Step on the scale to turn it on. It will automatically calibrate to 0.0 kg.

6. The scale will not turn on?

- Ensure the battery is not fully discharged.

7. The scale only shows weight and BMI?

- For first use by each user, the online mode must be used and the app connected for successful measurement.

- For every use, you must hold the correct grip, paying attention to the thumb and finger positions.
- Keep feet apart and stand barefoot on the scale platform, facing the display.
- Remain still for about 20-30 seconds to complete the full measurement.

- This symbol indicates that this product must not be disposed of with household waste within the EU. Disposal can take place via local collection and return points.
- Batteries should be removed before disposal and disposed of separately from the device. To protect the environment, non-rechargeable and rechargeable batteries must not be disposed of with regular household waste but must be taken to suitable collection points. Please also observe the applicable legal regulations regarding battery disposal.

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LEBENLANG GmbH