
# VENDOMNIA

ORIGINAL INSTRUCTION MANUAL / ORIGINAL-BEDIENUNGSANLEITUNG  
MODE D'EMPLOI ORIGINAL

EN DE FR

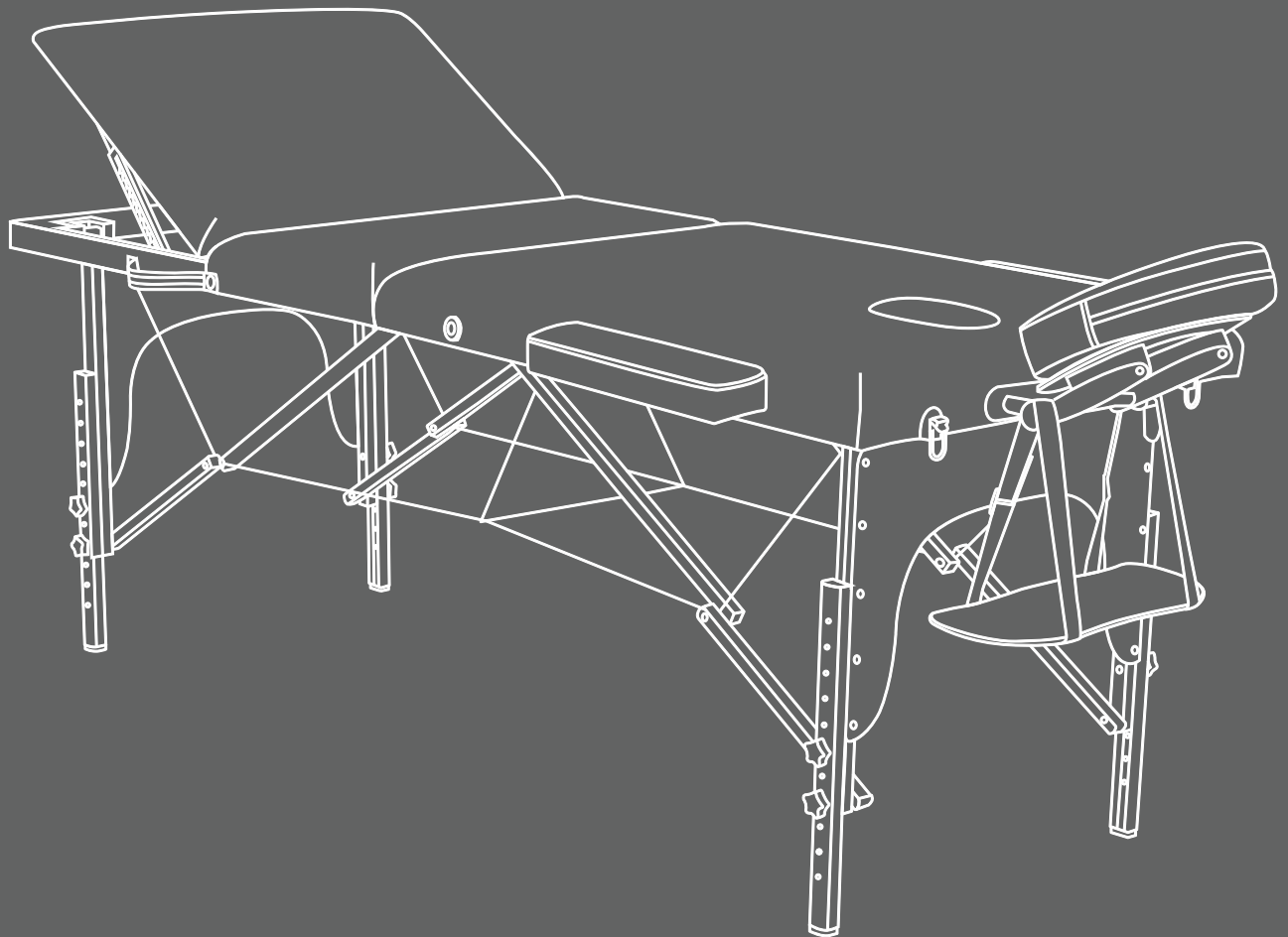ART. NR.: MSSLG

# PORTABLE MASSAGE TABLE

## DEAR CUSTOMER,

THANK YOU FOR CHOOSING OUR PRODUCT.

PLEASE READ THIS MANUAL CAREFULLY BEFORE ASSEMBLING OR USING THE PRODUCT TO AVOID DAMAGE CAUSED BY IMPROPER USE. PLEASE PAY SPECIAL ATTENTION TO THE SAFETY INSTRUCTIONS. IF THE PRODUCT IS PASSED ON TO THIRD PARTIES, THIS MANUAL HAS TO BE PASSED ON ALONG WITH THE PRODUCT.

PLEASE KEEP THIS MANUAL FOR FUTURE REFERENCE!

## SAFETY INSTRUCTIONS

1. EXAMINE THE MASSAGE TABLE FOR DAMAGE BEFORE EACH USE.
2. USE THE MASSAGE TABLE ON A LEVEL SURFACE. OTHERWISE THE TABLE WILL BE RAMSHACKLE.
3. DO NOT USE ALCOHOL BASED CLEANERS WHEN CLEANING THE VINYL. THIS COULD CAUSE DAMAGE TO THE VINYL MATERIAL. USE WARM SOAPY WATER AND WIPE DRY.
4. MAKE SURE WHEN USING THE TABLE THAT ALL FOUR LEGS SHOULD BE ADJUSTED TO THE SAME LEVEL. USING THE LEGS AT DIFFERENT LEVELS CAN CAUSE THE MASSAGE TABLE TO BE UNSTABLE AND CAUSE INJURY.
5. DO NOT ALLOW CHILDREN TO PLAY ON THE MASSAGE TABLE.
6. DO NOT OVERLOAD THE TABLE BEYOND THE RECOMMENDED WEIGHT LIMIT. THIS COULD RESULT IN SERIOUS INJURY TO THE USER.
7. **MAXIMUM WEIGHT CAPACITY - SAFE WORKING LOAD: 250 kg.**
8. KEEP THE TABLE FAR FROM NAKED LIGHT OR SMOKING. USING THE TABLE UNDER THE TEMPERATURE 40°C.
9. IF THE FRAME/LEGS/TABLE STRUCTURE IS DAMAGED IN ANYWAY, THEN YOU SHOULD NOT USE THIS TABLE. USING THIS TABLE WHEN DAMAGED MAY CAUSE HARM OR SERIOUS INJURY TO YOU OR YOUR PATIENT. YOU SHOULD SEEK REPAIR IMMEDIATELY.

## 10. THE INTENDED USER:

- THE TABLE IS AVAILABLE FOR NORMAL ADULT.
- FOR DISABLED INDIVIDUALS, PLEASE USE THE TABLE UNDER THE GUIDANCE OF MASSAGER.
- THE TABLE IS NOT AVAILABLE FOR PREGNANT WOMAN AND CHILDREN.

## TECHNICAL DATA

PRODUCT MODEL	MSSLG01 / MSSLG02	MSSLG03 / MSSLG04
SIZE WITHOUT ARM-REST AND HEADSET	LENGTH X WIDTH: 185 × 60 cm	LENGTH X WIDTH: 180 × 60 cm
SIZE WITH ARMREST AND HEADSET	LENGTH X WIDTH: 212 × 80 cm	LENGTH X WIDTH: 208 × 80 cm
SIZE FOLDED	93 × 60 × 16 cm	90 × 60 × 16 cm
HEIGHT ADJUSTABLE FROM	61 - 84 cm	59 - 79 cm
NET WEIGHT	13.5 kg / 15.6 kg	12 kg / 14.3 kg
WORKING WEIGHT	250 kg	250 kg



## PARTS LIST

1	HANDLE	2	HINGE	3	BUCKLE
4	END PANEL	5	ADJUST LEG	6	FOOT RUBBER
7	SIDE ARMRESTS	8	FACE CUSHION	9	HEADREST CRADLE
10	CABLE	11	FACE PLUG	12	PLASTIC TRAY
13	KNOB	14	ARMREST		

## ASSEMBLY

1. RELEASE THE BUCKLES AND THEN OPEN THE TABLE SO THAT IT IS STANDING ON THE SIDE WITH THE RUBBER FEET. TAKE OUT THE FACE CUSHION, U SHAPE HEADREST, ARMREST AND SIDE ARMRESTS FROM UNDERNEATH OF THE TABLE.

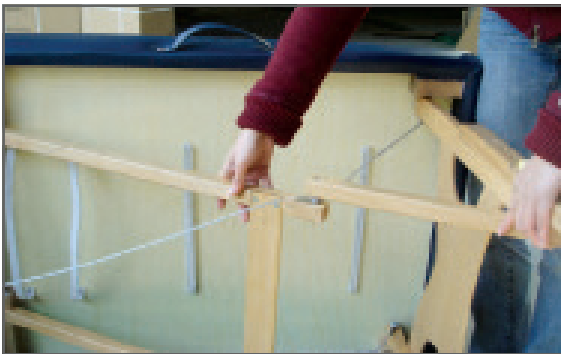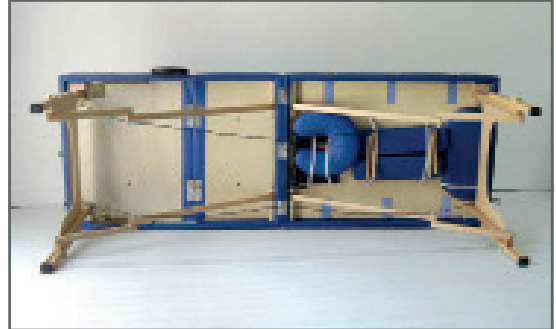
2. PULL OUT THE LEFT AND RIGHT PAIR OF LEGS AND STRAIGHTEN EACH JOINT.

3. LIFT THE TABLE UP UNTIL ALL FOUR LEGS ARE ON THE FLOOR.

4. LIFT ONE END OF THE TABLE AND PULL OUT TO OPEN TABLE FULLY.

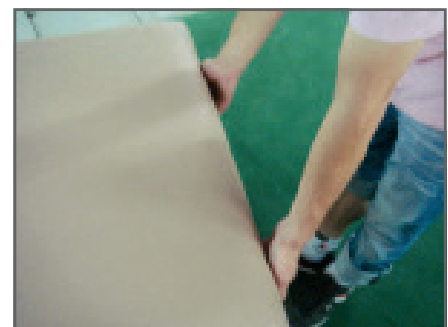5. INSERT BOTH SIDE ARMRESTS INTO THE SIDE ARMREST BUSHINGS OF THE FRONT BED.

**WARNING:**

THE SIDE ARMREST IS FOR LIGHT USE ONLY (MAX. 5 KG)!

PLEASE DO NOT PLACE HEAVY WEIGHT OR STAND UP BY SUPPORTING ON THE SIDE ARMREST.

6. PUT THE ARMREST TO THE HEADREST CRADLE. THEN INSERT BOTH RIGHT/LEFT HEADREST STICKS WITH HEADREST INTO THE HOLES OF THE HEADREST BUSHING PLATE.

7. PLACE THE FACE CUSHION ONTO THE U SHAPED HEADREST BY LINING UP THE VELCRO PADS.

**WARNING:**

HEADREST IS FOR LIGHT USE ONLY (MAX. 25 KG)!


PLEASE DO NOT PLACE HEAVY WEIGHT OR STRONG PRESSURE WHEN RESTING.

8. NOW THE MASSAGE TABLE IS READY FOR USING.


9. IF YOU WANT TO USE THE BACKREST OF THE TABLE, PLEASE RELEASE THE BUCKLE FIRST. AND THEN LIFT THE BACKREST.
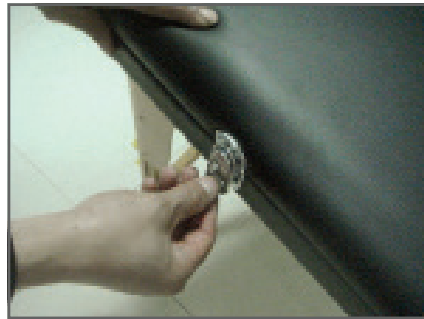
#### ATTENTION:

PLEASE LIFT THE BACKREST TO THE SAME LEVEL. THERE ARE 10 DIFFERENT LEVELS FOR CHOICE.

## ADJUSTMENT

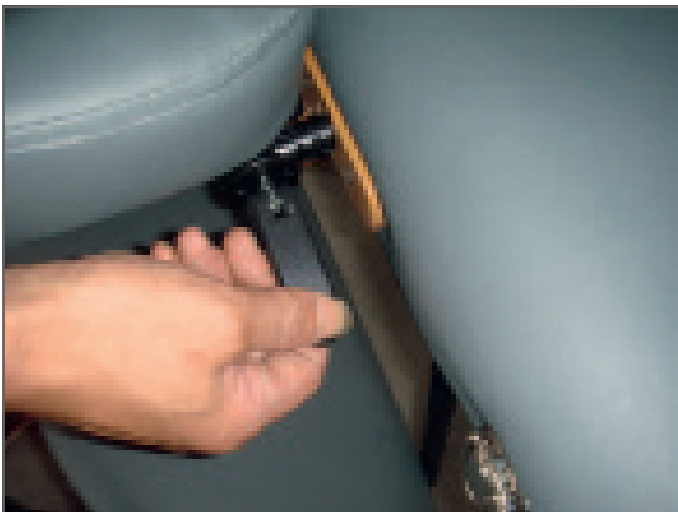
#### ADJUSTING THE ANGLE OF THE HEADREST

PRESS UP ON THE LOCKING HANDLES AND THEN ADJUST THE HEADREST TO THE DESIRED ANGLE. PRESS DOWN ON THE LOCKING HANDLES AFTER ADJUSTABLE TO LOCK INTO PLACE.

## ADJUSTING THE TABLE HEIGHT

1. LAY THE TABLE ON THE RUBBER FOOTED SIDE AND THEN REMOVE THE LEG KNOBS.

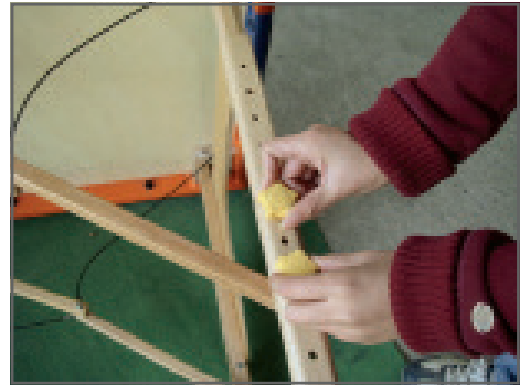2. REMOVE THE ADJUST LEG AND THEN PLACE THE ADJUST LEG BACK ONTO THE BASE LEG MATCHED UP TO THE DESIRED HOLES FOR ADJUSTING THE HEIGHT OF THE TABLE.

3. TIGHTEN THE LEG KNOBS AFTER ADJUSTMENT.

### ATTENTION:

ADJUST ALL OF THE FOUR LEGS TO THE SAME HEIGHT.

## ADJUSTING THE BACKREST

RELEASE THE BUCKLE AND LIFT THE BACKREST TO THE DESIRED LEVEL.


10 LEVELS

## STORAGE

1. LAY THE TABLE ON ITS SIDE AND FOLD AT AN ANGLE.
2. COMPLETELY FOLD THE LEGS FLAT UNDER THE TABLE.
3. MAKE SURE THE CABLES, U-SHAPE HEADREST, FACE CUSHION, ARMREST AND SIDE ARMRESTS ARE INSIDE THE TABLE BEFORE CLOSING IT.
4. STORE THE TABLE IN DRY PLACE TO PROTECT IT FROM MUSTINESS.
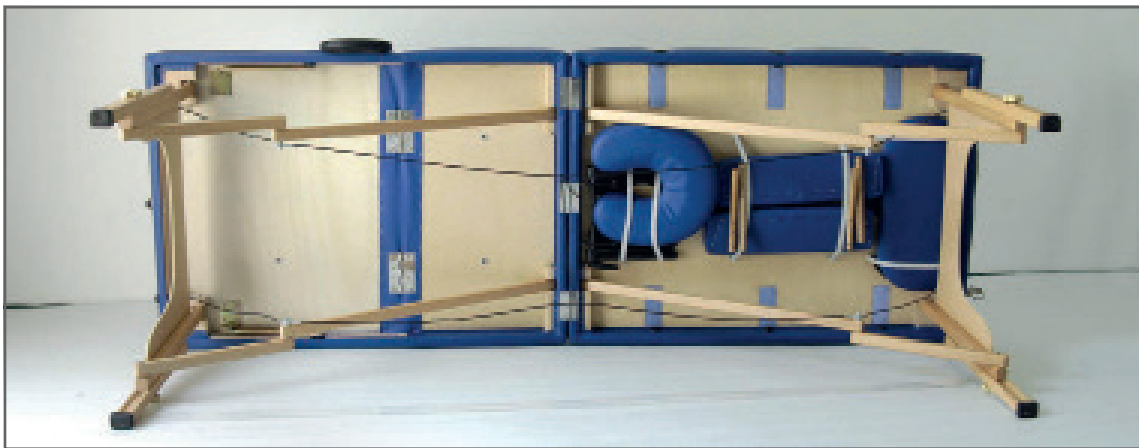
## TRANSPORT

1. THE TABLE IS EASY TO CARRY WITH THE SHOULDER STRAP OF CARRY BAG.
2. PLEASE PUT THE TABLE INTO THE TRUNK OF THE CAR DURING TRANSPORT.

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