VIBRATONE® MAX is a registered EU/CTM trade mark.
International patents pending. All imitations will be prosecuted.

Vibra tone™ MAX

Instruction manual

Manual de instrucciones

Manuale d'istruzioni

Manuel d'instructions

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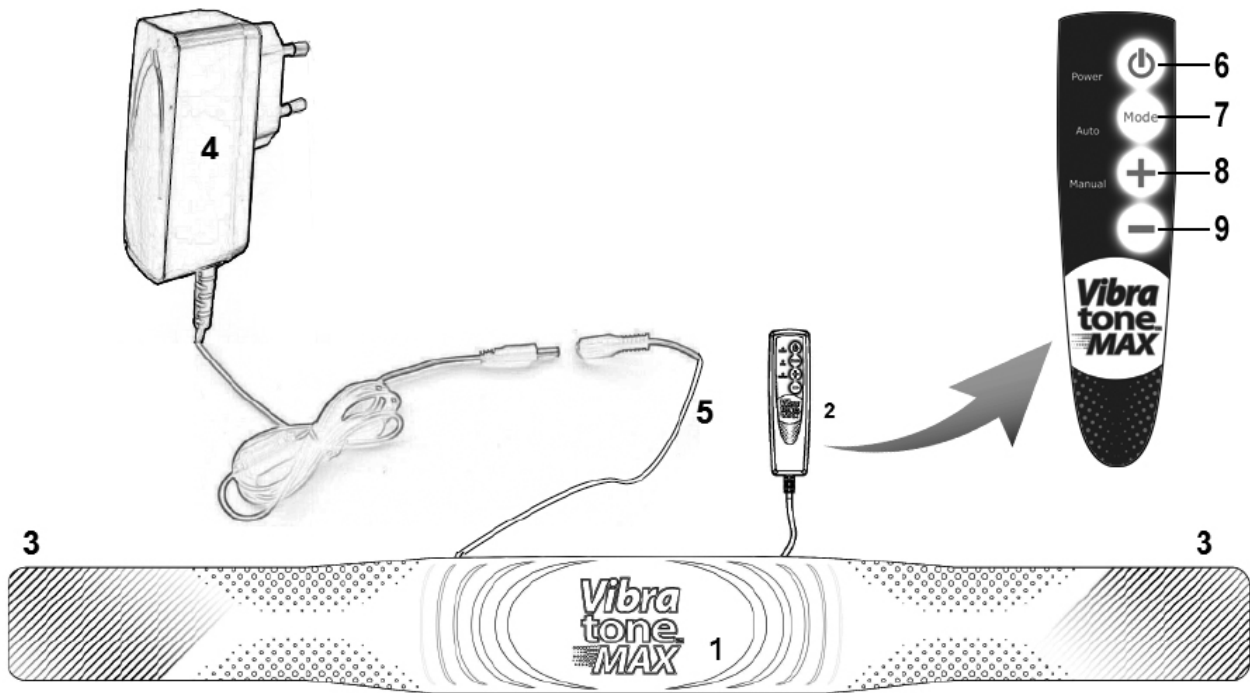
Congratulations on your purchase of the Vibra Tone™ Max, the most advanced vibrating belt on the market today. Muscle-vibration technology provides major benefits, so start enjoying them now! The Vibra Tone™ Max will become your closest ally. Take it with you wherever you go. Be sure to read all the instructions carefully before using the Vibra Tone™ Max for the first time.

IMPORTANT

- Do not use the Vibra Tone™ Max if you are pregnant.
- Keep the Vibra Tone™ Max out of reach of children.
- Do not use the Vibra Tone™ Max if you suffer from thrombosis or are receiving medical treatment.
- Do not use the Vibra Tone™ Max 30 minutes before or after meals.
- If you have suffered from a serious illness, consult your doctor before using the Vibra Tone™ Max.
- Do not use the Vibra Tone™ Max near water or with wet hands.
- Store the Vibra Tone™ Max in the carry bag provided and keep it away from sharp objects.
- Never place heavy objects on top of the Vibra Tone™ Max.
- Stop using the Vibra Tone™ Max immediately if you notice any kind of discomfort.
- Use only the cables and adapter provided with the device.
- Unplug the device from the mains after each use.
- Do not try to disassemble, repair or modify the Vibra Tone™ Max.
- Do not use the Vibra Tone™ Max at high temperatures or in very humid environments.
- The Vibra Tone™ Max is not a medical or therapeutic device.
- This device may be used by children aged 8 and over, by people with reduced physical, sensory or mental abilities, and by people lacking in experience and knowledge if they are appropriately supervised or given suitable instructions on the safe use of the device and understand the dangers its use entails. Children should not be allowed to play with the device. The cleaning and maintenance to be performed by the user should not be done by children without supervision. Children should be supervised to ensure that they do not play with the appliance.

CONTENTS

- 1 Belt with control pad
- 1 Adapter
- 1 Carry bag
- 1 Instruction manual
- 1 Diet plan
- 1 Measuring tape

PRODUCT DESCRIPTION

- 1. Belt
- 2. Control pad
- 3. Rear fastening system
- 4. Adapter
- 5. Belt connector cable

- 6.  button
- 7. Mode button
- 8. + Button
- 9. - Button

PLUGGING ON THE BELT

1. Separate the two parts of the rear fastening system (3).
2. Wrap the Vibra Tone™ Max (1) belt around the part of your body you want to work on and fasten the ends at the rear so the belt fits snugly.
3. Connect the adapter (4) to the connector cable (5) coming out of the belt.
4. Then plug the adapter (4) into the mains.

USING THE CONTROL PAD

When the device is plugged into the mains, press the  button (6) on the control pad (2). You will hear a beep and the red Power light on the control pad will switch on to indicate that the belt is ready for use. You can now choose between two operating modes: Automatic or Manual.

Automatic Mode

Automatic mode consists of a combination of vibration rhythms specially designed to make each use of the Vibra Tone™ Max a pleasurable experience.

1. Press the Mode button (7) on the control pad once to activate Automatic Mode. The red Auto light on the control pad will light up to indicate that Automatic Mode has been activated. Automatic Mode consists of a selection and combination of vibration rhythms that provide all the benefits of a vibratory massage.
2. Adjust the intensity by pressing the "+" button (8) to increase the intensity and the "-" button (9) to lower it.

Manual Mode

Manual mode consists of one constant vibration rhythm specially designed for easy selection of the right intensity for each moment and need.

1. Press the Mode button (7) twice to activate Manual Mode. The red Manual light on the control pad will light up to indicate that Manual Mode has been activated. Adjust the intensity by pressing the "+" button (8) to increase the intensity and the "-" button (9) to lower it.

You can go from one mode to the other by pressing the Mode button (7).

The Vibra Tone™ Max will automatically shut off after 30 minutes of use.

When you are finished, press the  button and the Power light will go out. You can now unplug the device from the mains.

CLEANING AND MAINTENANCE

The Vibra Tone™ Max requires practically no cleaning or maintenance. If you need to clean your Vibra Tone™ Max,

all you have to do is wipe it down with a soft damp cloth. Make sure the device is completely unplugged from the mains before cleaning it.

TECHNICAL SPECIFICATIONS

Belt Technical Specifications

Input: 12V 
Power: 12 W

Adapter Technical Specifications

Input: 100V-240V
Output: 12V  1000 mA
Frequency: 50/60 Hz

Vibra Tone™ Max complies with Directive 2014/35/EU (LVD).
Vibra Tone™ Max also complies with Directive 2014/30/EU (EMC).

Quality Warranty:

This product is covered by a warranty against manufacturing defects subject to the time periods stipulated by the legislation in force in each country. This warranty does not cover damages resulting from inadequate use, negligent commercial use, abnormal wear and tear, accidents or improper handling.

DOUBLE INSULATION

This appliance is double insulated; therefore, no earth wire is required. Always check that the mains voltage corresponds to the voltage on the product rating plate.

CAUTION: DISPOSAL OF MATERIALS

The symbol of the crossed-out wheeled bin means you should find out about and follow local regulations about disposing of this kind of product. Do not dispose of this product as you would other household waste. Dispose of this device in accordance with the corresponding local regulations. Electrical and electronic devices contain hazardous substances that can have harmful effects on the environment and/or human health and should be recycled properly.

Australia & New Zealand only

This product has the benefit of certain Consumer Guarantees. These are prescribed by the Australian Consumer Law & the New Zealand Consumer Guarantees Act 1993 both of which provide protection for consumers. There is no express warranty for this product in Australia or New Zealand. The above paragraph refers to other countries.

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Industax GmbH
ISL ITALY S.r.l.
Best Direct International Ltd.