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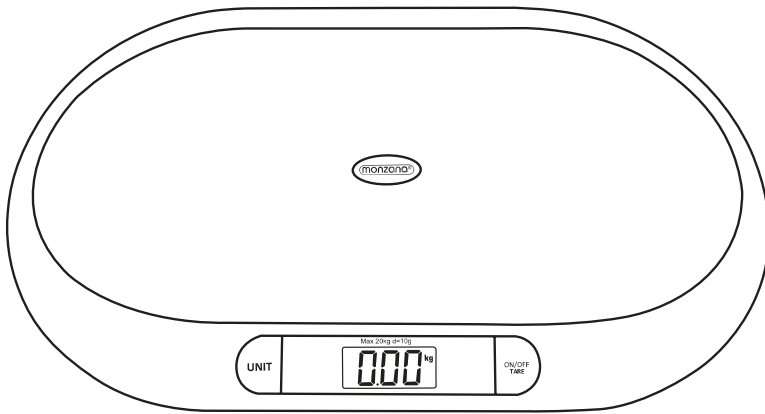
BDA-105021-001 / 21.06.2023

## Babywaage

Artikel-Nr: 105021

Modell-Nr: DBBW001

monzana

[www.DEUBAXXL.de](http://www.DEUBAXXL.de)

# ANLEITUNG

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# MANUAL

Dear customer,

Thank you for choosing our product, you have made the right choice by purchasing one of our brand products. Our products are regularly inspected to ensure that they meet the desired high quality standards and, of course, always comply with the stringent requirements of the European Union.

## **READ THE USER INSTRUCTIONS. IMPORTANT INFORMATION. KEEP IN A SAFE PLACE FOR FUTURE USE!**

Read the entire manual carefully before using the product for the first time. Keep this manual for future reference. If you pass on the product one day, be sure to hand over this manual as well.

### **WARNING! POTENTIAL DANGER!**

Observe the safety information and setup instructions to avoid the risk of injury or damage to the product.

### **DANGER OF SUFFOCATION!**

Keep small parts and packaging material out of the reach of children!


### **WARNING!**

The appliance is not a toy and should not be handled by children.

### **PRODUCT MODIFICATION**

Never make modifications to the product! Modifications will invalidate the warranty and the product may become unsafe or, in the worst-case scenario, dangerous.

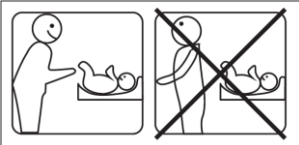**IMPORTANT!** Babies placed on the scales should be kept under constant watch to prevent them from falling out.

## **GENERAL SAFETY INFORMATION**

- Ensure that water cannot get inside the appliance. Never immerse the appliance in water or other liquids.
- Do not use the appliance if it has previously been dropped.
- Do not use the appliance if it exhibits visible signs of damage.
- You must protect the scales from shocks, moisture, dust, chemicals, strong temperature fluctuations and heat sources that are located too close to it (ovens, radiators).
- Do not place any objects on the scales when they are not in use.
- The precision of the scales can be affected by strong electromagnetic fields (e.g. induction hob, mobile phones, etc.).
- Never overload the scales! The scales have a max. load capacity of 20 kg (44 lb). This includes the weight prior to using the TARE function. The results of the weight measurement are displayed in 10 g steps (0.02 lb).
- Remove any transport protection devices that may be present.

## **SPECIAL SAFETY INSTRUCTIONS FOR APPLIANCES RUNNING ON BATTERIES AND RECHARGEABLE BATTERY PACKS**

- Do not use damaged or modified rechargeable battery packs/batteries or appliances. Defective or modified batteries can be unpredictable and may cause fire, explosion and injury.
- Do not expose a battery/rechargeable battery pack or appliance to fire or excessive temperatures. Fire and high temperatures exceeding 130°C (266°F) can result in an explosion.
- Never make modifications to the appliance! Do not attempt to modify or repair the appliance or the rechargeable battery pack/battery.
- In the event of battery leakage, remove the batteries from the product immediately to avoid damage.
- If a battery has leaked, avoid contact with skin, eyes and areas with mucous membranes. Should there be such contact, rinse the affected area with water and seek medical attention immediately.
- Keep batteries out of the reach of children. Seek medical attention immediately in the event of a swallowed battery!
- Only use the type of battery specified.
- Make sure you insert the batteries the right way round, i.e. with correct polarity! This is shown in the battery compartment.

### **INTENDED USE**

The scales are designed exclusively for weighing babies. Use the product only for the purpose it is intended for. The manufacturer accepts no liability for damage caused by unintended use. Any modification to the product may negatively impact safety, cause hazards and invalidate the warranty.

**Only suitable for domestic use; not for commercial use.**

### **UNPACKING INSTRUCTIONS**

Open the box and remove the product. After unpacking, please check the product for damage and to ensure it is complete. Subsequent complaints cannot be accepted.

### **FIRST USE**

#### **Inserting and replacing batteries**

(2 x 1.5 V AAA batteries delivered with the product)

1. Open the battery compartment at the bottom of the scales.
2. Remove the old batteries if applicable.
3. Insert the new batteries into the battery compartment. Make sure the batteries are the right way round, i.e. that the polarity is correct.
4. Close the battery compartment again.

#### **Choosing the right place for it**

Place the scales on a clean, flat surface. If you place the scales on a soft or uneven surface, there is a risk that the baby could fall off them. It is essential to place them on a firm surface to ensure correct measurement.

## OPERATION

### Button functions

ON/OFF/TARE button:	Switch on, switch off, tare
UNIT button:	Set the desired measurement unit in kg or lb.
	Hold function to freeze the weight measurement

### How to weigh

1. To switch on, press the 'ON/OFF/TARE' key. The display indicates '0.00' after 2 seconds.
2. Press the UNIT button to switch the unit of measurement (kg, lb).
3. Place your baby on the scales. Be very careful when doing so and mind the baby's head.
4. The LCD display will indicate the baby's weight.  
**Important!** Do not touch the baby to avoid interfering with the weight measurement.
5. Press the UNIT button to freeze the displayed weight. After checking the weight, press the UNIT button again to reset the display to '0'. You can then take another measurement.
6. Press the 'ON/OFF/TARE' key to set the display to '0'. Empty the tray. Then place your baby in the tray – the weight will be shown on the display.
7. Press the 'ON/OFF/TARE' key to switch off the scales.

### SUGGESTIONS:

If you do not want to place the baby directly on the weighing surface but, for example, prefer placing a towel under the baby, you must lay the towel on the scales BEFORE switching them on. If the towel and other loads you wish to remove during the weighing process weigh more than 50 g, use the TARE function.

### Error displays

As soon as the scales are overloaded, the message 'O-LD' appears on the display. The message "LO" indicates that the batteries need to be replaced.

### MAINTENANCE AND CARE INSTRUCTIONS

Check all the parts for signs of wear. Do not continue to use the product if you find any defects.

Clean the surface of the scales with a soft cloth without using chemical cleaning agents. Make sure that no water enters the scales. Never rinse them under running water.

Regular cleaning, checking and maintenance increase the item's safety and service life.

### TECHNICAL DATA

Item number	105021
Model	DBBW001
Max. weight recommendation	20 kg / 44 lb
Measurement unit	kg / lb
Measurement accuracy	+/- 10 g / 0.02 lb
Battery	2x 1.5V type AAA

### DISPOSAL

At the end of your product's long service life, please dispose of the valuable raw materials properly so that proper recycling can take place. If you are not sure how best to proceed, local waste disposal companies or recycling centres will be happy to help.

## **CORRECT DISPOSAL OF THIS PRODUCT**

Never dispose of electronic appliances as household waste! Take defective or discarded appliances to recycling centres.

Within the EU, this symbol indicates that this product must not be disposed of as household waste. Old appliances contain valuable recyclable materials that should be recycled so that there is no harm to the environment or human health as a result of uncontrolled waste disposal. Therefore, please dispose of old appliances through suitable collection systems or return the appliance for disposal at the place where you bought it. They will then send the unit for recycling.

## **BATTERIES**

As the end user, you are legally obliged to return used batteries. Once used, you can return batteries to us or to designated collection points (e.g. local authority collection points or retailers) free of charge.

You can also send the batteries back to us. We will refund the postage for returning the used batteries.

The symbols shown on the batteries mean the following:

B: Battery must not be disposed of with household waste

Pb: Battery contains more than 0.004 % of lead by weight

Cd: Battery contains more than 0.002 % of cadmium by weight

Hg: Battery contains more than 0.0005 % mercury by weight.

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Copyright by
Deuba GmbH & Co. KG
Deuba LTD