
**NEO<sup>1</sup>**

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**AM - 065**

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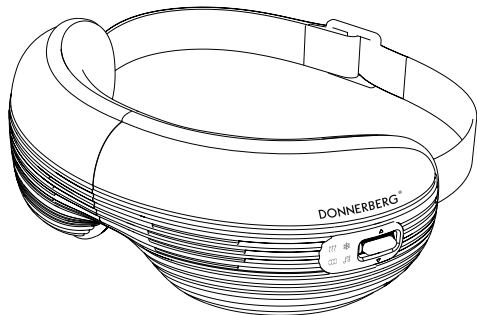
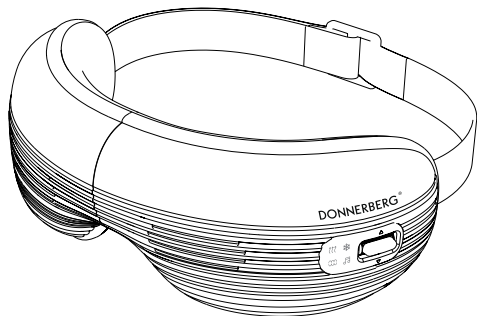
**USER MANUAL**

## Contents

Product features .....	18
Mode and functions.....	18
Instructions .....	19
Important information.....	19
Troubleshooting .....	20
Storage and cleaning.....	21
Safety precautions .....	21
Technical specifications - NEO1 AM-065 .....	23
Warranty policy (2 years).....	24
Bluetooth .....	24
Electric and electronic devices .....	25
Customer service .....	26
Battery removal guide.....	27

## Product features
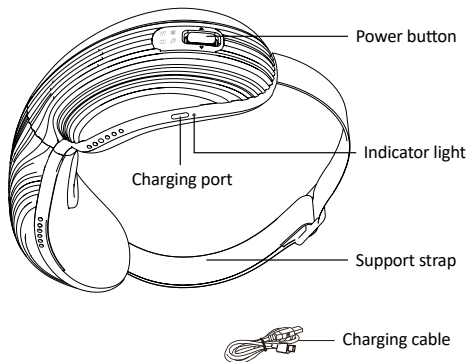
## Mode and functions

- **Dynamic:** Air pressure + Music \*\*air pressure can not be turned off in this mode
- **Heat:** Heat + Air pressure + Music
- **Cool:** Cool compress + Air pressure + Music
- **Auto:** Air pressure + Heat or Cool + Music

## Instructions

- Sit or lie down in a relaxed position.
- Adjust the length of the support strap to comfortably fit your head.
- Position the massager over your closed eyes, adjust strap if necessary.
- Hold the power button for 3 seconds to turn on/off.
- Press the power button twice to pause/resume music.
- Toggle power button UP to rotate modes (Dynamic, Heat, Cool, Auto). The default mode is Dynamic.
- Toggle power button DOWN to select Air Pressure mode (mode 1, mode 2, or air pressure OFF).

## Important information

- Voice Prompt “Low Battery” will sound when device requires charging.
- This device can not be used while charging.
- Default session time on this device is 10 minutes, the device will automatically turn off after each session. However COOL mode will run for 6 minutes by default.
- Charging time is 4 hours. Once completed, the indicator light will change from RED to BLUE. Unplug the charger from the device once the charging cycle is complete.
- Cooling fan will continue to run once device has been turned OFF as a safety precaution to cool the device for storage.

## Troubleshooting

This page summarises the most common problems you could encounter with the product. If you are unable to solve the problem with the information below, please contact Donnerberg customer service.

Problem	Possible Causes	Solution
The product does not turn on or the air pressure is weak	The massager has low battery	Charge the massager
The massager turns off	The massager has finished its 10-minute auto timer session and has stopped	If you want another massage, press the power button again
The product is abnormally warm	The product is likely to be overused	Let the product cool down before using it again
Noise can be heard from the product	Normal sound of air pump inflating or air pressure	No need to do anything

## Storage and cleaning

### Storage

- Ensure device is turned OFF before storing.
- Store in its box away from direct sunlight and charge device if it hasn't been used for more than 3 months.

### Cleaning

- Turn the device OFF and disconnect from mains power.
- Use a soft dry cloth to clean the massager shell.
- A slightly damp cloth may be used to clean the fabric massage area if needed, however avoid getting moisture into the product.
- Do not use harsh or corrosive chemicals to clean the product.

## Safety precautions

**In order to operate this product correctly and safely, please be sure to read the "User Manual" carefully before use.**

1. Use this massager only for its intended purpose as described in this manual.
2. This Massager is not intended to be a substitute for medical attention. If any symptoms persist or worsen after use, contact your medical practitioner.
3. When you are pregnant or suffer from any kind of medical condition, skin condition or illness or have had surgery and are unsure whether you should use this massager, please seek medical advice before use.
4. Stop using if you feel discomfort or pain during use.

5. Do not use if suffering from any kind of physical ailment that would limit operation of the controls.
6. The product is not intended for use by children or people with reduced physical, sensory or mental capabilities or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the product by a person responsible for their safety.
7. Some redness of the massaged area is normal and will disappear soon after massage.
8. Do not use after taking painkillers, sedatives or alcohol. Do not use while sleeping, driving, exercising or during a fever.
9. Remove jewellery and hair accessories that may impact the massage area.
10. Do not apply medical or cosmetic cream on the product or on your face before using.
11. Keep the product away from water, fire and corrosive materials.
12. Avoid the massager or cord getting damaged and do not use if damaged.
13. If damaged or faulty or has an unusual noise during operation, turn off Eye Massager with Heat & Cooling.
14. In order to prevent the risk of electric shock, do not open the product outer cover.
15. Do not place heavy objects on the massager.
16. Keep away from sharp objects.
17. To charge the battery, only use charger cable provided.
18. Check the power supply. If the power supply voltage exceeds the specified range, do not use.

## Technical specifications - NEO1 AM-065

Name	Eye Massager with Heat & Cooling
Model	NEO1 AM-065
Battery capacity	2500 mAh
Charging interface	TYPE-C 5V/1A
Charging time	Approx. 4 - 4.5 hours
Battery run time	Approx. 60-80 minutes
Input Voltage	5V
Rated voltage	3.7V
Rated power	3.5W
Net Weight	0.46 kg
Gross Weight	0.75 kg
Product size	210 x 115 x 80 mm
Heating temperature	38 °C +/- 2 °C
Cooling temperature	19.5 °C +/- 1 °C
Automatic switch off	After 10 minutes

## Warranty policy (2 years)

The company Donnerberg S.G., Von-Werthern-Str.4a 96487 Dörfles-Esbach, Germany grants a valid manufacturer's warranty for the territory of the European Union and the civil-law defect liability. The warranty starts from the time the device was delivered to the customer. Battery is excluded from this warranty. It is not extended by a previous warranty case. The warranty is granted only to the buyer and is not transferable to third parties.

In addition to the statutory liability for defects, we also repair or replace equipment whose defects are not initial defects unless they are caused by improper use or are insignificant. The choice between replacement and repair is ours. The warranty is only granted if the defect is notified immediately upon its occurrence and is reported. You must include the claim form and the copy of your purchase receipt as the proof of warranty. Please note that we will not process warranty claims without the claim form. Upon notifying our customer service team, we will send you a claim form.

## Bluetooth

Open the list of Bluetooth connections on your device and select "Neo 1". A beep will sound once the connection has been established and a voice prompt indicates the connection/disconnection of Bluetooth devices.

## Electric and electronic devices

### Information for private households

The Waste Electrical and Electronic Equipment Regulations 2018 (WEEE) in the United Kingdom contains several requirements concerning the handling of electric and electronic devices as waste. In the following paragraphs, we list the most important requirements for you.

#### 1. Separated recycling of „old devices“

Electric and electronic devices that have become waste are called old devices. Owners of such old devices have to bring these to recycling facilities, where these are handled and recorded separately from urban solid waste. Old devices must not end up in the household refuse. In fact they have to be collected in special collecting and recycling systems.

#### 2. Possibility for disposal of „old devices“

Owners of old devices from private households can dispose these at the public waste disposal authorities or bring them to producer or seller related disposal locations. You can find these disposal locations online: <http://www.recycle-more.co.uk>

#### 3. Privacy information

Old devices often contain private data, especially appliances from data processing and communication equipment, as smartphones or hard discs. In your own interest, please arrange to eliminate all personal and private data from your disposed equipment.

#### 4. Meaning of the symbol „crossed out waste container“

The labelled symbol showing the „crossed out waste container“ on electric and electronic equipment indicates that the concerning equipment has to be disposed separately from urban solid waste.

#### 5. The Waste Electrical and Electronic Equipment Regulations 2018 (WEEE)

Following the implementation of the Directive 2012/19/EU and its British transportation WEEE Regulations 2018, we as a producer of EEE are registered in the Public register of Electrical and Electronic Equipment Producers. The register is held in England and Wales by the Environment Agency.

### Customer service

We have a reliable customer service team that is committed to providing solutions to any issue you might have. If you have any questions or concerns about the products, warranty or exchange do not hesitate to contact us at:

E-mail	de.info@donnerberg.net
Telephone	+ 49 89 215 400 980 + 49 95 619 769 199
Donnerberg S.G. Von-Werthern-Str.4a, 96487 Dörfles-Esbach, Germany	

### Battery removal guide

This electrical device contains the following battery:

Battery type	Chemical composition
YT 123555PL 2500mAh 3.7V	Li-ion Polymer battery

Please do not remove the battery yourself within the warranty period. This invalidates the manufacturer's warranty. Please contact our customer service at [de.info@donnerberg.net](mailto:de.info@donnerberg.net).

**At the end of the product's service life, you can dispose of the battery yourself for recycling as described:**

**Warning: Make sure that the battery is completely discharged and that the device is switched off and the charging cable unplugged.**

1. Remove the one screws at the bottom cover of the device, the cover with Type-C port.
2. Use a flat-head screwdriver to insert it along the gap at the middle part of the cover, and then forcefully pry open the cover.
3. Remove the four screws at PCBA.
4. Pull out all terminals.
5. Carefully remove the battery
6. The battery and the device can now be disposed of separately.

Donnerberg S.G.  
Von-Werthern-Str.4a,  
96487 Dörfles-Esbach,  
Deutschland

Tel: +49 89 215 400 980  
+49 95 619 769 199  
Mail: [de.info@donnerberg.net](mailto:de.info@donnerberg.net)  
Web: [www.donnerberg.net](http://www.donnerberg.net)