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
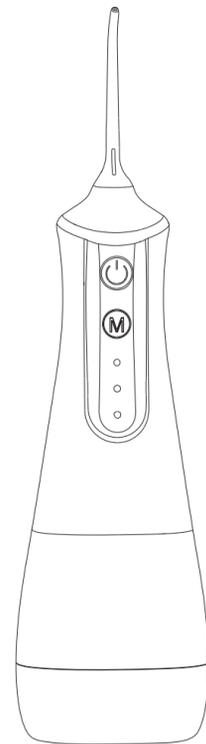
Intelligent Water Flosser Product Manual

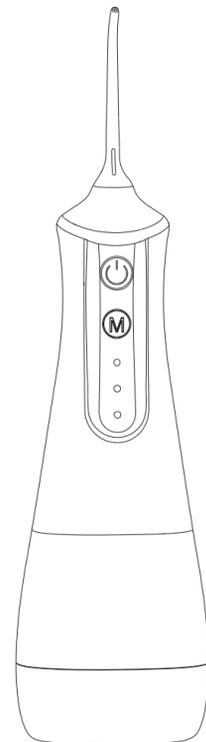Before using the product, please read this product manual carefully and keep it for future reference..

Safety Instruction

Danger

- If you do not understand the operating procedures and safety instructions, please do not use this product to avoid injury caused by improper operation.
- Please do not operate the product by means other than this manual.
- Do not disassemble this product or change the components of this product to avoid damage to this product or high-voltage electric shock, and the warranty is invalid.
- Do not use water with a temperature above 40°C in the water tank and clean the mouth to avoid burns.
- If the tooth cleaner falls into the water while charging, please unplug the power plug of the wall power adapter before removing the product.
- If the power cord or plug is damaged, replace it before using it.
- Do not put the product into fire or heat, and do not charge, use or place at high temperature.


Warning

- Do not use this product for other purposes than oral cleaning.
- Do not press the nozzle firmly on the gums or teeth to avoid damaging the gums.
- Please place this product in a safe place to avoid falling or falling into the water.
- Do not immerse this product in water to avoid failure or cause fire, electric shock, explosion and other hazards.
- Do not plug or unplug the power plug and charging plug with wet hands.
- Before connecting the power plug, please check whether the voltage marked by the power adapter is in accordance with the local voltage.
- Do not damage or modify the power cord, do not stretch, bend, or twist the power cord with force, nor can you place heavy objects on the power cord or jam the power cord between objects.
- Do not place the product, especially the power cord, near hot objects.
- Unplug the power plug and charging plug before maintenance.
- Do not place the product in high temperature or humidity, and avoid direct sunlight.
- Do not charge this product in hot or humid places, such as bathrooms and toilets.
- If not used for a long time, please unplug the power plug and charging plug to avoid fire or danger.

Attention

- This product has no replaceable or repairable parts. If you need repair, please return it to the manufacturer.
- Do not use if there is no nozzle.
- Some mouthwashes may corrode this product, please do not use mouthwashes in this product.
- This product is not intended for use by children, infants and persons who do not know how to handle it, unless they have been supervised or instructed in the use of this product by their safety officer. Keep out of reach of children and infants and do not use as a toy!
- Patients with periodontal disease or those who have had oral surgery in the last 2 months, please consult a doctor before use.
- When plugging or unplugging the power supply, be sure to hold the plug instead of the power cord.



Structure Introduction

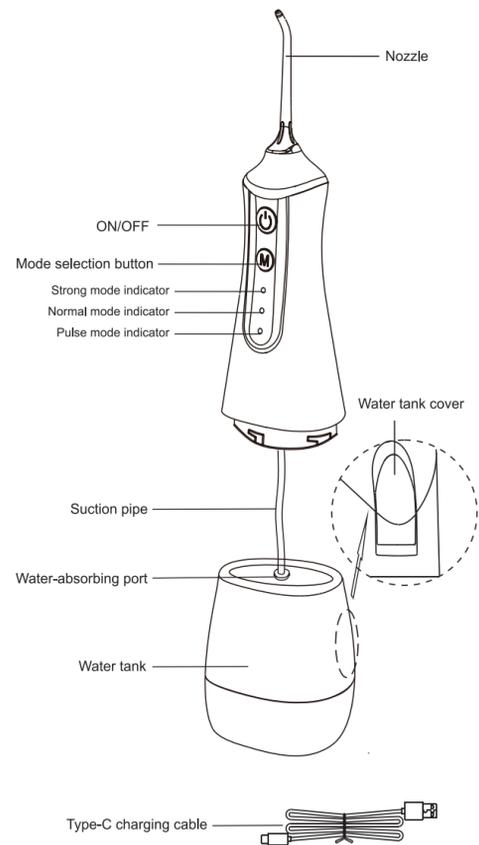Operation Introduction

Before use, please charge the product as follows:

1. Check to make sure that the product and hands are dry, and then connect the USB plug to the 5V adapter (the adapter is prepared by the user) and plug it into the socket (Figure 1).
2. Take off the charging cover, insert the other end of the USB cable into the charging socket of the punch, please pay attention to the charging plug and make sure it is inserted into place (Figure 2).


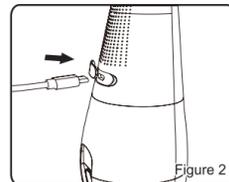
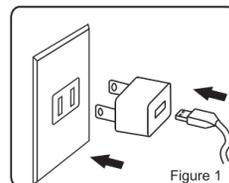
Warm tips:


If the indicator light flashes during or after use, the battery is low and needs to be recharged. It takes 3-4 hours to fully charge. When charging, the indicator light flashes to indicate that it is charging. The indicator light stop flashing to indicate that it is fully charged.

When using, please follow the steps below:

1. Insert the nozzle into the hole on the top of the water flosser and rotate it 90° to fix it (Figure 3-1). If you need to replace the nozzle, first shut down, then rotate the nozzle 90° and then pull out to replace it. (Figure 3-2).
2. Fill the water tank with water. There are two methods for filling water, as follows:
 - a. Open the water tank cover, hold the water flosser horizontally, put water into the water tank and then close the cover. (Figure 4).
 - b. Rotate the water tank in the unlocking direction to remove the water tank, hold the water tank vertically to fill it with water. Then

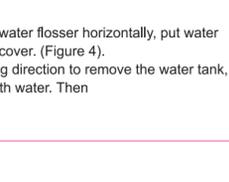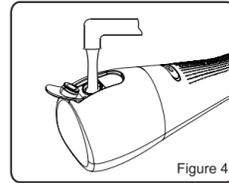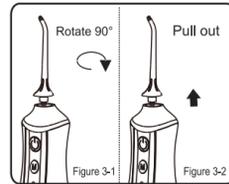push it up and rotate it in the direction of the lock to lock. (Figure 5).


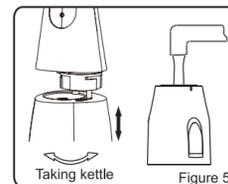
3. Select the desired mode by pressing the mode button. You can also switch to the desired mode by pressing the mode button during use. This product has three modes: Strong, Normal and Pulse.

The specific functions are as follows:


Strong mode: Deep cleaning, quick cleaning.
Normal mode: Standard water pressure to meet daily cleaning.
Pulse mode: Massage gums to promote blood circulation, improve gum health, suitable for sensitive teeth and first-time users.

Tips:

Dental professionals recommend that first-time users use the Pulse mode for 1 week and then use the other modes after getting used to it.

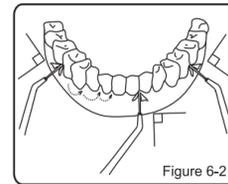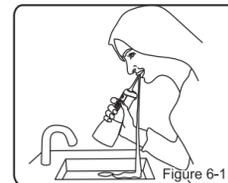
4. Before use, you need to hold the water flosser vertically, the nozzle extends into the cavity, aligns the teeth or gums, and the mouth opens slightly so that the water can flow out smoothly (Figure 6-1 & 2).

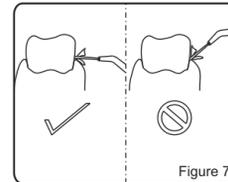5. Press the switch button to use. When using, please pay attention to the following points: (Figure 7)
 - a. The direction of water flow is perpendicular to the teeth or gums;
 - b. Move slowly along the teeth and clean the teeth in order;
 - c. The nozzle is aligned with the gum line and perpendicular to the gum.



Tips:

Do not spray water directly along the periodontal pocket. After using, please press the switch button to shut down. This product has a 2-minute timer function. It will automatically shut down after 2 minutes. If you need to continue to use it, please press the switch button.

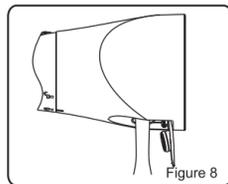After using, please follow the steps below:

1. Open the water tank cover and pour out the water (Figure 8).
2. Press the switch button to flush the teeth and automatically drain the accumulated water (Figure 9).
3. Press the switch button to shut down the water flosser and then close the water tank cover.
4. Dry the water flosser with a cloth (Figure 10).

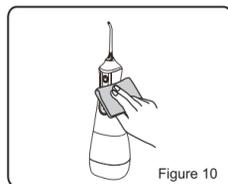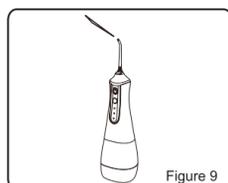Tips:

In order to avoid the spread of bacteria, please ensure that there is no residual water in the dental flusher. If you plan not to use it for a long time, please clean the product and store it after drying.

Routine Maintenance

Please use clean water or neutral detergent to clean the product, do not use corrosive or abrasive detergent, because these will damage the product. Do not wash with hot water above 50°C.

Cleaning of the body:

- Please clean the body with a dry cloth or a damp cloth.
- Do not immerse the body in water for cleaning.

Cleaning of water tank:

- It can be cleaned with a dry or wet cloth, or removed and rinsed with clean water.

- If you plan not to use it for a long time, please clean the water inside.

Cleaning of nozzle:

- It can be cleaned with a dry cloth or a damp cloth, or removed and rinsed with clean water, then wipe dry with a cloth.
- Do not bend, pull or twist the nozzle hose.
- Dental experts recommend replacing the nozzle every 6 months.

Troubleshooting

Problems	Possible causes	Solution
Stop working	New purchase or long time no charge leads to low battery	Please charge for more than 3 hours
Low water pressure during use or no water spray	There is something blocked in the nozzle	Replace the nozzle
	The suction inlet is blocked	Clean the suction inlet
	It's been sitting idle for a while	Pull out the nozzle to flush the nozzle hole on the main unit When the tank is full of water, pick it up horizontally and switch it on
Only works for a short time after charging	The battery has reached the end of its useful life	Please contact the distributor
	Insufficient charging time	Please charge for more than 3 hours
The charging indicator does not light up during charging	Power connector not connected properly	Please connect the power supply correctly
	Product charging circuit is faulty	Please contact the distributor

Waste disposal

This product uses rechargeable batteries. At the end of product life, before discarding the product, be sure to remove the battery and recycle or dispose of it in accordance with local waste management regulations. If in doubt, please consult your local waste management department.