**LEBEN
LANG**

USER MANUAL

KÖRPERANALYSE- WAAGE

BODY COMPOSITION SCALE

LBL702

Download on the
App Store

GET IT ON
Google Play

Deutsch	3
English	12
Français	19
Español	26
Italiano	33
Nederlands	40

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Thank you for trusting LEBENLANG. We pursue the highest standards of quality and design and wish you great pleasure with your new product.

Before installing your OKOK International app, carefully check whether your tablet or mobile phone meets these minimum requirements:

- Compatible with iOS 7.0 or higher and Bluetooth 4.0 or higher.
- Compatible with Android 4.3 or higher and Bluetooth 4 or higher.

Scan the QR code to
download the APP

Install the OKOK International App:

1. Search for “**OKOK International**” in the App Store or Google Play or scan the above QR code.
2. Download and install the app on your device.
3. After installation is complete, the app icon will appear on your phone or tablet.

GENERAL SAFETY INSTRUCTIONS

When using electrical appliances, basic safety precautions should always be followed, including the following:

- Please read this user manual to ensure effective use. If you pass the product on to someone else, also hand over this manual.
- LEBENLANG devices comply with applicable safety regulations. Repairs and replacement of the connection cable may only be carried out by specialists authorized by the manufacturer. Improper repairs can pose significant dangers to the user.
- The scale is intended solely for private and not for commercial use.
- Children must not play with the scale. Cleaning and maintenance must not be carried out by children.
- The scale may be used by children aged 8 years and older and persons with reduced physical, sensory, or mental capabilities or lack of experience/knowledge if they are supervised or instructed in the safe use of the scale and understand the hazards involved.
- Keep the scale away from heat, direct sunlight, moisture, sharp edges, etc.
- Regularly check the scale for damage. Do not use the device if it is damaged.
- Do not use the scale outdoors.
- Do not repair the scale yourself; contact an authorized expert.
- If the scale is damaged, it must be replaced by the manufacturer, its service personnel, or similarly qualified persons to avoid hazards.
- Always place the scale on a hard, dry, and flat surface before performing a

measurement.

- When inserting batteries, pay attention to the correct polarity (+/-). Incorrect insertion may damage the device.

SPECIAL SAFETY INSTRUCTIONS

- Do not stand on the edge of the scale or jump on it.
- Do not overload the scale (Max. 180kg).
- Do not drop the scale or drop objects onto it, as this may damage the sensors.
- Do not immerse the scale in water and do not use organic or chemical cleaning agents. A lightly dampened cloth is sufficient to clean the scale.
- Any measurement taken with this scale is for reference only and should not be considered as medical advice.
- Consult your doctor or physician before making changes to your diet, exercise plan, or physical activities.
- Remove the batteries if the scale will not be used for an extended period.

PACKAGE CONTENTS & TECHNICAL DETAILS

- Intelligent Bluetooth Body Scale (max. 180kg)
- 2 x AAA batteries
- 1 x User manual
- Weight unit: kg/lb/st (units switchable)
- Optimal room temperature: 5~35°C
- Optimal storage temperature: -10~60°C
- Optimal humidity: <90%RH

OKOK INTERNATIONAL APP – FEATURES

Body weight	Body fat	Visceral fat	Muscle mass
Bone mass	BMR	Biological age	BMI

1. **Data Analysis:** Track your body metrics such as weight, body fat percentage, muscle mass, and water content over time.
2. **Graphical Representations:** Visualize your progress with charts and statistics to better understand changes in body composition.
3. **Set Goals:** Set personal health and fitness goals and monitor your progress.
4. **Multiple User Profiles:** Create profiles for different users, allowing several people to use the scale and save their data separately.
5. **Integration with Fitness Apps:** Sync your data with other popular fitness or health apps (e.g., Apple Health, Google Fit) for a comprehensive overview of your health.


- 6. **Reminders and Notifications:** Receive reminders for regular weigh-ins or set goals.
- 7. **Data Backup:** Store your data in the cloud for easy restoration when needed.

USAGE/ OPERATION

Step 1: Open the OKOK App.



Pair the device, register, and fill in user information

- Register as a new user in the app
- Set how you want your information to be displayed
- Enter your personal details


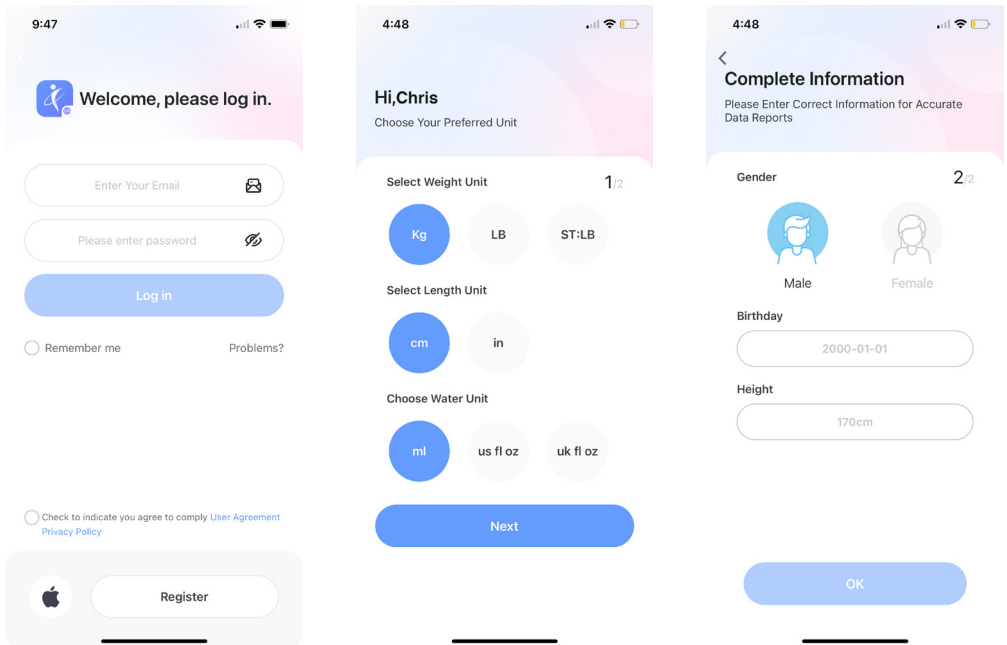
Step 2: Weight Calibration

- Add the scale to the app
- Step on the scale barefoot and always in the correct direction without syncing the data with your phone.

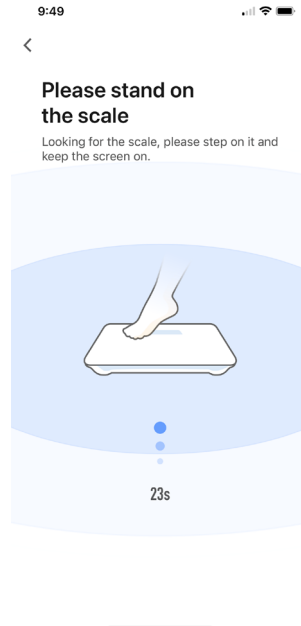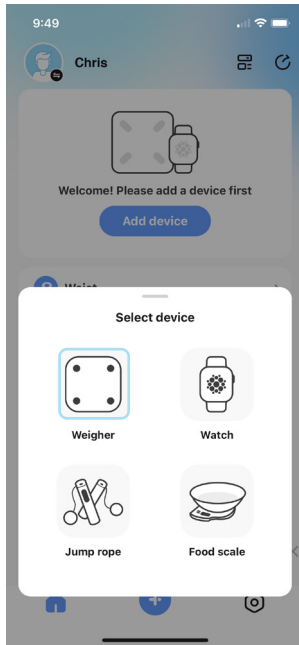Step 3: Results

- After 10-30 seconds, your results will be displayed.
- The app shows the user's weight information and other users' data. It also provides a history report and gives health-related advice based on the measurements.
- The app provides a weekly, monthly, and quarterly report. These changes/progress are displayed through corresponding trend curves.

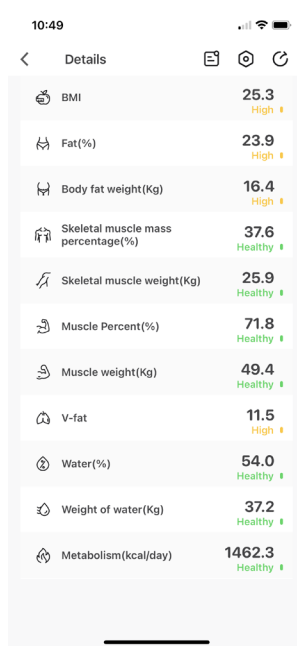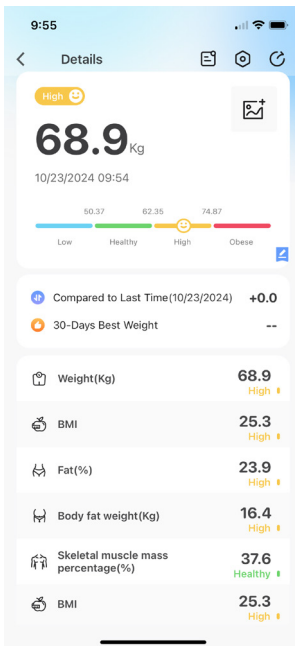
Daily Weighing

- Open the OKOK App.
- You will be prompted to turn on Bluetooth.
- Once Bluetooth is on, the pairing process starts automatically.
- Begin weighing as usual
 - » The weight data will be instantly synced with the app.
- If the Bluetooth connection is lost during weighing, the weight data will be stored in the scale's memory and synced with the app once the Bluetooth connection is re-established.

Note:

- You can change the weight units displayed on the scale (KG/LB/ST) via a button on the back of the scale.
- This change can also be made in the app.

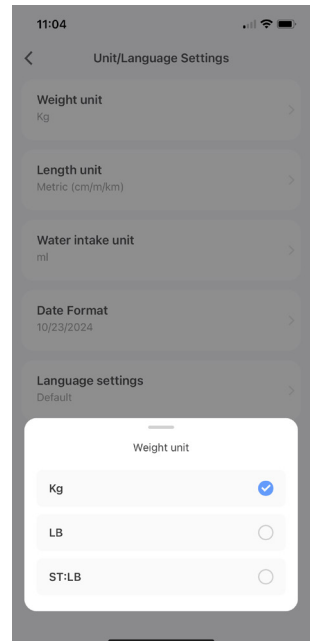TROUBLESHOOTING

1. No display on the screen

- Cause: Batteries are dead or not correctly inserted.
- Solution: Check the batteries and replace them if necessary. Ensure the batteries are correctly inserted (+/- poles).

2. Weight display is inaccurate

- Cause: The scale is on an uneven surface.
- Solution: Ensure the scale is on a solid, level surface.

3. Bluetooth connection does not work

- Cause: The Bluetooth signal is disrupted or the device is out of range.
- Solution: Ensure your smartphone or tablet is near the scale (max. range 10 meters).
- Cause: Bluetooth is not activated on the smartphone/tablet.
- Solution: Check if Bluetooth is enabled on your smartphone or tablet, then search for the scale in Bluetooth settings again.
- Cause: The scale is connected to another device.
- Solution: Ensure the scale is only connected to one device at a time. Disconnect from the previous device if necessary.

4. Scale turns off during use

- Cause: Batteries are low.
- Solution: Replace the batteries with new ones.

5. Weight freezes on display

- Cause: Software error or scale overload.
- Solution: Turn the scale off and back on. If the problem persists, remove the batteries for a few minutes and reinsert them.

6. No synchronization with the app

- Cause: Outdated app version.
- Solution: Update the app to the latest version in the app store.
- Cause: Outdated smartphone operating system.
- Solution: Ensure your smartphone's operating system is up to date.

7. Display shows nothing even though the batteries are full

- Cause: Possible contact issue or corrosion.
- Solution: Check the battery contacts and clean them carefully with a dry cloth.

BATTERY INSTRUCTIONS

Please observe the following guidelines when handling batteries to ensure your safety and protect the environment:

- Disposal: Never dispose of batteries in household waste. Take them to

designated collection points or recycling centers.

- **Replacement:** Only use the recommended battery types. Do not mix old and new batteries or different types of batteries.
- **Empty Batteries:** Remove empty batteries from the device immediately to avoid damage from leaking batteries.
- **Storage:** Store batteries in a cool, dry place. Avoid extreme temperatures and direct sunlight.
- **Safety:** Keep batteries out of reach of children and pets. Swallowing batteries can be life-threatening. In such cases, seek medical help immediately.
- **Short Contact:** Avoid contact between battery terminals and metal objects to prevent short circuits.

By following these guidelines, you contribute to a safe and environmentally friendly handling of your batteries.

- This symbol indicates that this product should not be disposed of with other household waste in the EU. Disposal can be done through locally available collection and return points.

- The batteries should be removed before disposing of the device and disposed of separately. To protect the environment, non-rechargeable and rechargeable batteries must not be disposed of with regular household waste but should be taken to appropriate collection points. Please also observe the corresponding legal regulations for battery disposal.

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LEBENLANG GmbH