For FAQs and more information, please visit:

**[www.eufylife.com](http://www.eufylife.com)**

Health & Wellness

## Owner's Manual

Smart Scale P1

Smart Scale P1

Báscula Inteligente P1

Balance Intelligente P1

Bilancia Smart P1

01

English

07

Deutsch

13

Español


19

Français


25

Italiano

## What's In the Box

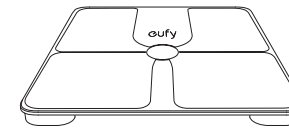Smart Scale P1

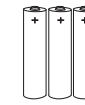
AAA Battery (x3)

Owner's Manual &  
Other Documents

## Inserting the Batteries


- 1 Turn the scale over and place it on a flat surface.
- 2 Open the battery door and insert three AAA batteries. Make sure the positive and negative ends are facing the correct polarity directions marked in the battery compartment.  
- The digits "0000" will show on the LED display.
- 3 Close the battery door and wait until the digits "00" are shown on the LED display.

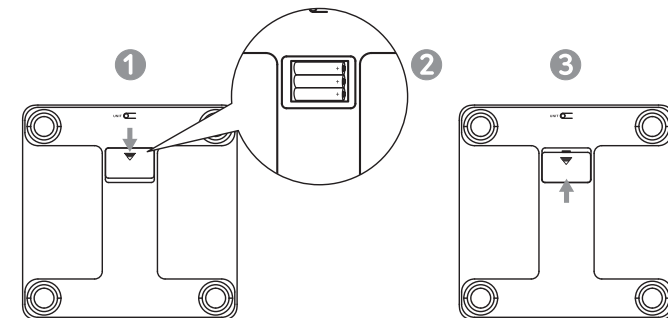**LO** will appear on the LED display when the batteries are low. Replace all three batteries at the same time. Do not mix new and used batteries.


## Selecting the Measuring Unit

Press the **UNIT** button at the bottom of the scale to cycle through and select the measuring unit: **lb** (Pound) / **kg** (Kilogram). The measuring unit can be selected from within the EufyLife app as well.
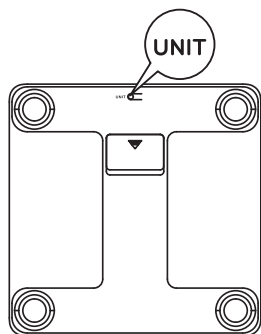
## Installing the EufyLife App and Pairing

- 1 Download the EufyLife app from the App Store (iOS devices) or Google Play (Android devices).
- 2 Enable Bluetooth on your smartphone or tablet.
- 3 Open the EufyLife app and follow the on-screen instructions to pair the scale with your device.
  - Once paired, the app and scale will automatically connect the next time the app is opened.

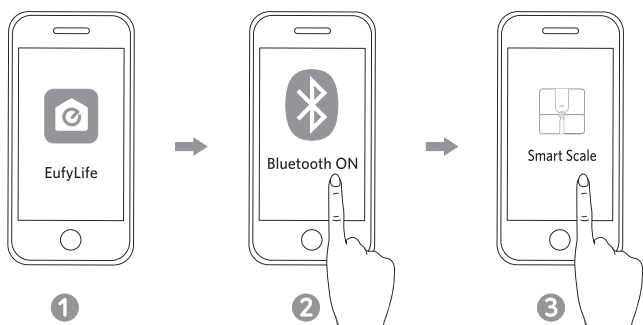	Status Icon	Indicates
	Flashing	Pairing
	Solid	Connected

- This scale is compatible with devices running iOS 9.0 and above or Android 4.4 and above.
- Do not pair the scale from your smartphone's Bluetooth menu. Use the EufyLife app to pair the scale with your smartphone.


## Starting a Measurement

Before taking your first measurement, you will need to create and set up a user ID in the EufyLife app.

- 1 Place the scale on a hard and flat surface.
- 2 Open the EufyLife app and select Smart Scale on the device homepage.
- 3 Step barefoot onto the scale. An analysis of your body's composition will begin automatically.

- Step barefoot onto the center of the scale. Stand still and keep full contact with the electrodes until the measurement is complete.
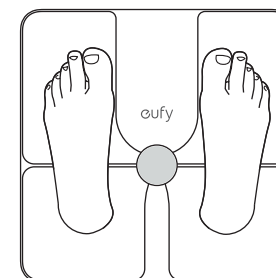
## Weight-Only Measurement

This scale can also be used as a conventional weight-reading scale.

To measure your current weight, you can

- step on the scale wearing your shoes; or
- select the weight-only measuring mode in the EufyLife app.

## Maintenance

Regularly clean the scale with a dry or slightly dampened soft cloth.

- Do not wash the scale with water or immerse it in water.
- Never use detergents to clean the scale as this may cause it to malfunction or discolor.

## Troubleshooting

### Error Message

Error Display	Description	Solution
OL	Overload. The device will power off.	Stop using this scale for measurement.
Lo	Low battery. The device will power off.	Replace all three batteries at the same time.
Err	Measurement failure.	<ul style="list-style-type: none"> <li>• Step barefoot and keep full contact with the electrodes.</li> <li>• Stand still until the measurement on the LED display locks.</li> </ul>

## When Measuring

Issue	Cause	Solution
Abnormal measurement results: <ul style="list-style-type: none"> <li>• Too high</li> <li>• Too low</li> <li>• Huge difference between two recent measurements</li> </ul>	<ul style="list-style-type: none"> <li>• Incorrect posture.</li> <li>• Device is placed on soft carpeting or uneven surface.</li> </ul>	<ul style="list-style-type: none"> <li>• Step onto the scale barefoot and stand still.</li> <li>• Place the device on a hard and flat surface.</li> </ul>
LED display does not turn on after stepping onto the scale.	<ul style="list-style-type: none"> <li>• Batteries not installed.</li> <li>• Worn batteries.</li> </ul>	<ul style="list-style-type: none"> <li>• Install three new AAA batteries</li> <li>• Replace all three batteries at the same time.</li> </ul>
Scale does not proceed to analyse body composition.	<ul style="list-style-type: none"> <li>• Taking a measurement while wearing shoes or socks.</li> <li>• Using the weight-only measuring mode in the EufyLife app.</li> </ul>	<ul style="list-style-type: none"> <li>• Stand barefoot on the scale and keep feet in full contact with the electrodes when taking a measurement.</li> <li>• Change the measuring mode in the EufyLife app.</li> </ul>

## Specifications

Power Source	4.5 V (3 x AAA batteries)
LED Display	2.6 x 1.2 in / 67 x 29 mm
Measurement Unit	Pound / Kilogram
Unit Increments	0.1 lb / 0.1 kg
Body Fat Unit Increments	0.1 %
Body Water Unit Increments	0.1 %
Muscle Mass Unit Increments	0.1 %
Measurement Range	11-397 lb / 5-180 kg
Working Environment	Temperature: 5 - 40 °C Humidity: 20 - 90% RH
Storage Environment	Temperature: -20 - 50 °C Humidity: <= 90% RH
Operating Frequency Range	2402 - 2480 MHz
Radio-Frequency Power	4 dBm
Dimension	11.0 x 11.0 x 0.9 in / 280 x 280 x 24 mm
Net Weight	3.31 lb / 1.5 kg

## Customer Service

For FAQs and more information, please visit [www.eufylife.com](http://www.eufylife.com).

✉ Email: [support@eufylife.com](mailto:support@eufylife.com)

💬 Live chat: [www.eufylife.com](http://www.eufylife.com)

☎ Customer service hotline:

(US) +1 (800) 988 7973 Mon-Fri 9:00am - 5:00pm (PT)

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(UK) +44 (0) 1604 936200 Mon-Fri 6:00 - 11:00 (GMT)

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(DE) +49 (0) 69 9579 7960 Mo-Fr 6:00 - 11:00

 @EufyOfficial
  @EufyOfficial
  @EufyOfficial
  @EufyOfficial

Anker Innovations Limited

Room 1318-19, Hollywood Plaza, 610 Nathan Road, Mongkok, Kowloon, Hong Kong

## AUSTRALIA & NEW ZEALAND WARRANTY

This document details the terms and conditions of the warranty provided with this product.

Please disregard any other representations contained on this product's packaging or provided with this product, which constitute or appear to constitute a warranty.

Subject to the exclusions and limitations set out below, we warrant to you that if you purchase a product directly from us, or, from any of our authorized resellers, that it will not contain defects which are a result of faulty manufacturer workmanship or materials. This warranty subsists for 15 months from the date of purchase. This warranty is provided in addition to your rights under the Australian & New Zealand Consumer Law.

Subject to any rights you have under the Australian & New Zealand Consumer Law, we will not be liable for defects which are a result of damage or loss caused by factors beyond our control, any product which has not been maintained according to any care or maintenance instructions provided with the product, any alterations to a product which are not performed by us, or damage or defects caused to the product due to unusual or non-recommended use. We will not be liable for any special, indirect, consequential or economic loss or damage or loss of profits whatsoever suffered by you or any other person resulting from any act or omission by us. Our total liability in respect of this Warranty against Defects is limited to the replacement of the goods.


This Warranty against Defects is provided in addition to other rights and remedies that you may have at law. Our goods come with guarantees that cannot be excluded under the Australian Consumer Law. You are entitled to a replacement or refund for a major failure and for compensation for any reasonably foreseeable loss or damage. You are also entitled to have the goods repaired or replaced if the goods fail to be of acceptable quality and the failure does not amount to a major failure.

In order to claim the warranty, you must return the product to the retailer from which it was purchased or if that retailer is part of a National network, a store within that chain, along with satisfactory proof of purchase. The retailer will then return the goods to Directed Electronics Australia authorized by Anker Innovations Limited. Directed Electronics Australia will repair, replace or refurbish the product at its discretion. The retailer will contact you when the product is ready for collection. All costs involved in claiming this warranty, including the cost of the retailer sending the product to Directed Electronics Australia, will be borne by you.


This warranty is offered by Directed Electronics Australia.

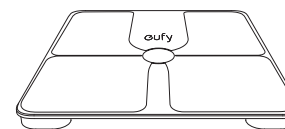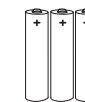Customer service support number (+61) 3 8331 4800 Mon-Fri 8:30am - 5:00pm (AEST) and customer service email address [support@myanker.com.au](mailto:support@myanker.com.au).

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