GYM form[®]



GLUTES



 ***Instruction Manual***
 ***Manual de instrucciones***
 ***Mode d'emploi***

 ***Bedienungsanleitung***
 ***Istruzioni per l'uso***
 ***Manual de instruções***
 ***Gebruiksaanwijzing***

International patents pending. All imitations will be prosecuted.
Gymform® Glutes is a registered EU/CTM trade mark.

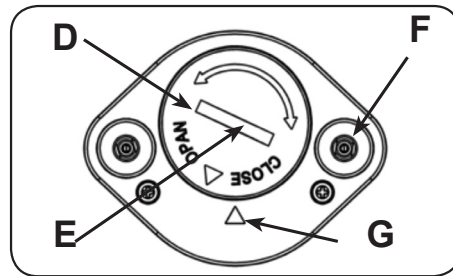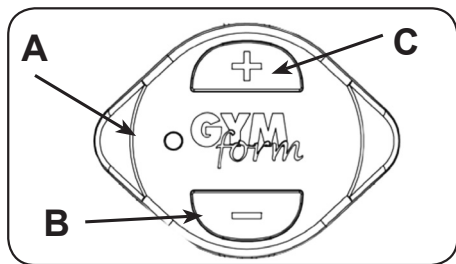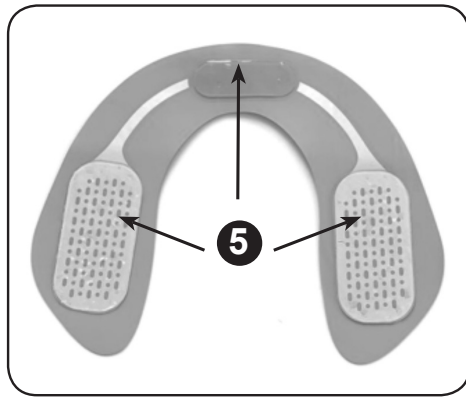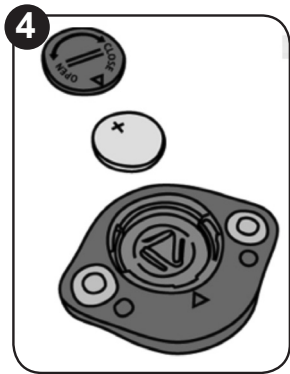Gymform® Glutes is the system of Electro Muscle Stimulation (EMS) designed to work the gluteal muscles, helping to shape, lift and firm up the buttocks.

WARNINGS

Before using Gymform® Glutes, consult your GP and physiotherapist.

Be sure to read the instructions, warnings and precautions carefully before use.

Gymform® Glutes is designed for adults in good health. However, it needs to be used with care. Read the following warnings and instruction manual before using this product.

- Application of electrodes close to the chest may increase the risk of cardiac fibrillation.
- Simultaneous connection of a

user to high frequency surgical equipment may cause burns at the contact points of the Gymform® Glutes and the Gymform® Glutes may be damaged.

- The operation of shortwave therapeutic equipment or a microwave in the vicinity (1m for example) can cause instabilities in the power output of the Gymform® Glutes.

- This is a powerful muscle stimulation device, injuries may occur if you use the most advanced exercise programme directly. Start at lower levels and work your way up.

- Do not use on children or animals.
- Do not use while driving or sleeping.
- Do not use in the bathroom or other places with high humidity.
- Batteries must be removed from

the appliance before disposal.

- Batteries must be removed for safe disposal.

- Non-rechargeable battery is not to be recharged.

- New and used battery cannot be mixed.

- The supply terminals should not be short-circuited.

- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.

You should NOT use Gymform® Glutes:

- If you have an electronic implant (e.g. a pacemaker) or suffer from any other heart problem.

- If you are pregnant or think you may be pregnant.

- If you suffer from cancer, epilepsy or are under medical supervision for cognitive dysfunction.

- During breast-feeding.

- If you have a high fever.

Ladies, to use Gymform® Glutes you must wait until:

- At least six months after the birth of your child (in consultation with your doctor).

- One month after the insertion of an intrauterine contraceptive.

- At least three months after a caesarean section (in consultation with your doctor).

- The first few days of your period have passed.

To use the Gymform® Glutes unit you must obtain permission from your doctor or physiotherapist:

- If you have a serious illness or injury not mentioned in this section.
- If you have recently undergone surgery.

- If you take insulin because of diabetes.

Always remember when applying the product:

- Do not apply to: the neck, genitals, scars, moles, varicose veins, or head.

- Do not apply to recent scars, inflamed skin, infected or acne prone areas, thrombosis, other vascular problems (varicose veins), or any part of the body with limited sensitivity.

- Do not place on areas with injuries or limited mobility (fractures or sprains).

- Do not place the product directly on metal implants. They should be placed over the nearest muscle.

Possible adverse reactions:

There are isolated reports of skin reactions to the use of electro stimulation units, including allergies, prolonged skin redness and acne. If you experience any of these symptoms, turn off the unit and consult your doctor.

If you have never done an exercise programme before, or if it has been a long time since you last did it:

Consult your family doctor.

Cautions:

• Do not use Gymform® Glutes while sleeping or drowsy.

• Do not use Gymform® Glutes

when your body is in contact with water.

• Keep out of the reach of children
• Do not use Gymform® Glutes in conjunction with any other device that transfers electrical current to the body or with any other machine.

• Do not connect the metal brackets or patches to any other object.

• Do not modify the unit.

• If the gel pad/patch has not been applied correctly, you will feel a stinging sensation.

• Never use damaged, deteriorated or dirty gel pads/pads, otherwise you may experience pain or burns.

• During use, make sure that the gel pads are in contact with the electrodes. If the electrode touches the skin directly, injury or burns may occur.

• Attach the product to the area

in question and then switch the appliance on.

• When you have finished, switch off the appliance before detaching it from the body.

• If you have any abnormal sensation on your skin or body, stop use immediately and consult a doctor. Failure to do so may result in an accident or injury.

• If you use the product at full power or for a long time, it may cause muscle fatigue, tingling or redness of the skin.

• Remove any jewellery or accessories such as necklaces, belts or earrings before using the device. Failure to do so may result in an accident or injury.

• Do not use the device with electronic devices such as watches, pedometers, mobile phones, etc., as their operation

may be affected.

- Do not use products such as body cream or lotions with the device as this may reduce adhesion.

- The gel patch adheres better to the skin if the hair has been removed from the area. But do not use it on the same day you have epilated.

- Do not heat or cool the product before use, as it may become deformed. Also do not use it in a way that is subject to sudden temperature changes, such as pouring hot water over it or in a sauna, etc.

- Clean the area where you want to use it.

- The product is not waterproof. Do not use it in a place with high humidity or where it may come into contact with water.

- When you have finished using it, turn it off before removing it from your body.

PRODUCT DESCRIPTION

1. Gymform® Glutes console (figure 1)
2. self-adhesive patch (figure 2)

- A. Ignition light
- B. Off/Lo key
- C. On/Hi key
- D. Battery cover
- E. Slot
- F. Press button
- G. Arrow mark

On / +: To turn the unit on or to gradually increase the intensity up to level 10.

Off / -: To turn the unit off or to gradually reduce the intensity down to level 1.

BATTERY INSTALLATION.

Gymform® Glutes is powered by a CR2032 battery (not included).

1. Turn the device upside down to open the battery cover by inserting a coin into the battery cover slot on the back of the microprocessor device and turn

counterclockwise towards the arrow mark. (see figure 3)

2. Insert the battery so that the + pole is facing upwards, making sure that the battery is correctly positioned. (see figure 4)

3. Put on the battery cover and using a coin turn the battery cover clockwise to close it.

USE Gymform® Glutes

1. Attach the console to the Gymform® Glutes silicone unit using its brackets. Make sure it is perfectly fixed.

- 2- Remove the transparent film protecting the self-adhesive patches (keep the protective film for storage after use). (figure 5)

- 3- Next, attach the pad to your buttocks, the patch must be securely attached to the area to be trained. (see figure 6)

- 4- Press the “+” button once to turn the device on, and the indicator light will turn green and go to intensity level 1.

- 5- Press the “+” button again to increase the intensity of the massage level. Press “-” to decrease the intensity of the massage level.

- 6- We advise you to start with low intensity levels and gradually increase

the intensity as you become familiar with the device.

You can switch off the device at any time during use by pressing the “-” key for 2 seconds.

- 7- Gymform® Glutes has only one automatic programme. The programme will automatically switch off after 20 minutes.

10 minutes a day is enough to get great results, if you want more, use it for the programmed 20 minutes a day.

- 8- Gymform® Glutes will automatically switch off if the self-adhesive patches are not in contact with the skin.

- 9- When the battery needs to be changed, and the voltage is below 2.4V, the red light will flash for 5 seconds and the Gymform® Glutes will automatically switch off.

ATTENTION:

The weaker the battery becomes, the weaker the electrostimulation will become. If the battery runs out, the keys may not work properly. Replace the battery with a new one

MAINTENANCE

Main device

- Always remove the battery when the product is not in use for a long period of time to prevent damage from battery leakage.
- Always keep the main device clean by using a soft, dry cloth to wipe the surface of the unit.
- Do not spill liquid on the device or immerse it in water.
- Keep the product in a cool, dark place, but NOT hot or humid.
- Keep out of reach of children and pets.

Buttock Pad

- The buttock pad is consumable. When it loses its adhesiveness, the adhesive surface (hydrogel) dries out or becomes trapped with dirt that cannot be removed. It must then be replaced with a new gel pad. Use only Gymform® Glutes pads. Do not use pads from other similar devices.
- Use your fingers to remove water droplets and wipe the adhesive surface of the hydrogel (adhesive surface) to remove dirt.
- After draining the water thoroughly, place the adhesive side (hydrogel side) face down in a well-ventilated

shade to dry.

- Dry it naturally until there is no moisture on the surface, then replace the hydrogel on the protective film.
- Keep the hydrogel pads clean in a dust-free and dry place, and keep them away from oily or sticky places.
- Do not apply to damaged skin. If a skin rash occurs, discontinue use and contact your doctor.

PRODUCT SPECIFICATIONS

Initial curved shape	Two-phase square-wave pulse
Pulse duration	320µs
Pulse rate	3.3 ~ 33.3Hz
Output voltage	36V+/-5V Load resistance 1000 ohm
Input voltage	3V
Duration of treatment	Automatic shutdown after 20 minutes
Intensity	Adjustable from 1 to 10
Mode	1 frequency mode
Operating conditions:	10°C~40°C (50°F~104°F); 30~85% RH
Storage conditions	-10°C~50°C (14°F~122°F); 10~95% RH

This symbol indicates that the device is protected against electrostatic discharge with double insulation type BF.

ISO 7000-1641 User's Guide

BATTERY PRECAUTIONS

Be sure to install it correctly before turning on the unit.
To change the battery, switch off the unit.
Remove the battery if the unit will not be used for a long period of time.
Take the necessary precautions when disposing of the battery and do not dispose of it in a fire..

HOW TO DISPOSE OF MATERIALS

The crossed-out wheeled bin symbol indicates that you should inform yourself and follow local regulations for the disposal of this type of product.
Do not dispose of this product in the same way as you would dispose of general household waste. You must dispose of it according to your local

regulations.

Electrical and electronic products contain hazardous substances that have harmful effects on the environment or human health and must be properly recycled.
The battery must be disposed of separately. To protect the environment, batteries, including rechargeable batteries, should not be disposed of in the same way as normal household waste. Take them to a collection point and check that you comply with the applicable legal requirements for battery disposal.

Quality assurance:

This product is covered by a warranty against manufacturing defects subject to the time periods stipulated by the legislation in force in each country.
This warranty does not cover damage resulting from misuse, negligent business use, abnormal wear and tear, accidents or mishandling.

Australia & New Zealand only
This product has the benefit of certain Consumer Guarantees. These are pres-

PROBLEM SOLVING

Problem	Cause	Solution
THE LIGHT DOES NOT ILLUMINATE	<ul style="list-style-type: none">- Battery not installed- Battery installed incorrectly- dead battery	<ul style="list-style-type: none">- Installing the battery- Check the polarity “+” towards you- Replace battery
THE PILOT LIGHT TURNS ON BUT THE UNIT DOES NOT WORK	<ul style="list-style-type: none">- No skin contact- The console is not correctly attached to the Gymform® Glutes.- Intensity level in 0- The microprocessor has crashed- Battery charge is low	<ul style="list-style-type: none">- Make sure the gel patch is in direct contact with the skin.- Check that the console is properly connected with the brackets to the Gymform® Glutes unit.- Press the “+” button to increase the intensity.- Reset by removing the battery and reinserting it after 1 minute.- Replace the battery.
NONE OF THE BUTTONS WORK WELL	<ul style="list-style-type: none">The microprocessor crashedBattery charge is low	<ul style="list-style-type: none">- Reset the unit by removing the battery and reinserting it after 1 minute- Replace the battery
The gel pad does not attach properly on the glutes	The adhesive surface (hydrogel surface) is dirty or damaged	Replace the gel pad with a new one
Redness or prickling on the skin	The hydrogel pads don't adhere properly to the skin	Adjust the placement of the hydrogel pad. If the pad is attached on a surface that is too small you will feel some prickling sensation.

cribed by the Australian Consumer Law & the New Zealand Consumer Guarantees Act 1993 both of which provide protection for consumers. There is no express warranty for this product in Australia or New Zealand. The above paragraph refers to other countries.

Made in China

*This product is not medical equipment.
*Do not disassemble or maintain this product by yourself.

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ISL GmbH
Best Direct (International) Benelux B.V.
Best Direct International Ltd,
ISL ITALY S.r.l.