
STARKBOOTS®

RECOVERY & HEALTH

STARKFELD® - 110% or nothing!

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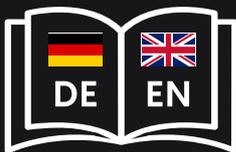CONTENTS

Safety information	02
FAQs-(Frequently Asked Questions)	05
Specifications	06
Legal information	07
Customer service – STARKSURANCE®	08
Troubleshooting	08
Product Description	09
Operating instructions	10
After use	12
Cleaning	13
Storage & maintenance	13
Recycling	13
Scope of delivery	14

SAFETY INSTRUCTIONS

Important note:

Please read these safety instructions carefully before first use and keep them in a safe place for future reference.

Warnings:

1. No modifications to the device:

Do not modify or change the device under any circumstances. Any unauthorised modification can lead to malfunctions and represents a safety risk.

2. Pain or unusual symptoms:

If severe pain or unusual symptoms occur during use, switch off the appliance immediately by pressing the main switch. Remove the boots and consult a doctor immediately.

3. Consult your doctor before use if you suffer from the following conditions or symptoms:

Inflammation or open wounds

Broken bones or recently healed fractures

Current medical treatments or ongoing therapies

Osteoporosis

Heart disease or the use of a pacemaker

cardiac pacemaker

Peripheral neuropathy or other nerve disorders

Physical injuries that make massage or compression unsuitable (e.g. broken legs)

Malignant tumours or cancer

4. Keep out of the reach of children and persons with reduced capabilities:

Keep the appliance out of the reach of children and persons who cannot operate it independently and safely.

5. Use the original power supply unit:

Only use the original power supply unit supplied to operate the device. A defective or incompatible power supply unit can cause damage to the device or pose a fire hazard.

6. Do not damage the mains cable:

Do not damage, modify, excessively bend, pull or twist the power cable of the adapter. This can lead to a short circuit or electrical damage

7. Do not use a defective power supply unit:

Do not use the device if the power supply unit is defective or visibly damaged.

8. Avoid wet hands:

Do not plug in the mains adapter with wet hands to avoid the risk of electric shock.

9. Do not use in ceilings or at high temperatures:

Do not place the control unit under ceilings or use the device in environments with high temperatures. This can cause the device to overheat.

10. No independent repairs or conversions:

It is prohibited to modify, disassemble or repair the appliance without authorisation. Contact customer service in the event of technical problems.

CAUTION

1. Only use the mains adapter supplied:

only use the original power supply unit supplied with the device.

A non-compatible power supply unit can damage the device or lead to safety risks.

2. Do not use in damp locations:

do not use the appliance in damp environments or outdoors to avoid the risk of electric shock.

3. Keep out of the reach of children:

keep the appliance, power cord and all accessories out of the reach of children and pets to avoid accidents or damage.

4. Unplug the appliance when not in use:

unplug the appliance when not in use to save energy and protect the appliance from possible damage.

5. Avoid movement during use:

Only use the device in a sitting or lying position. Avoid moving during use to ensure that the device functions correctly and to prevent injuries.


6. Protect the system from damage:

Store the appliance, the mains cable and all accessories in such a way that they cannot be damaged by children, pets or liquids

7. Designed for use on persons only:


The device is designed exclusively for use on humans and must not be used on animals.

MODE 1: SEQUENCE MODE

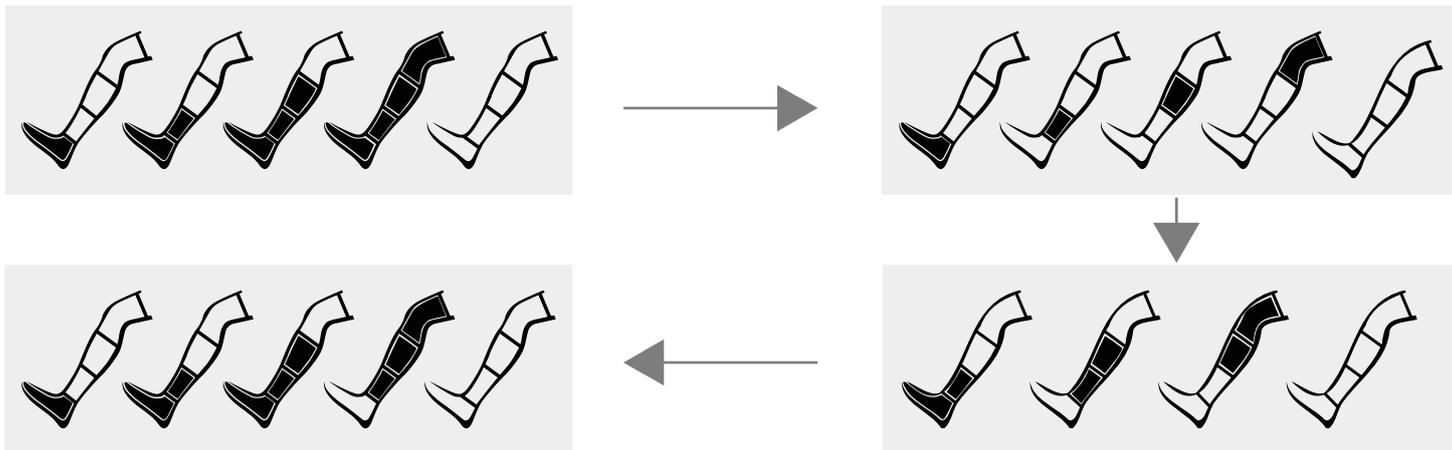
- Sequential pressure build-up from bottom to top.
- Each chamber is activated in turn.
- Ideal for even pressure without stressing areas that have already been treated.

MODE 2: CIRCULATION MODE

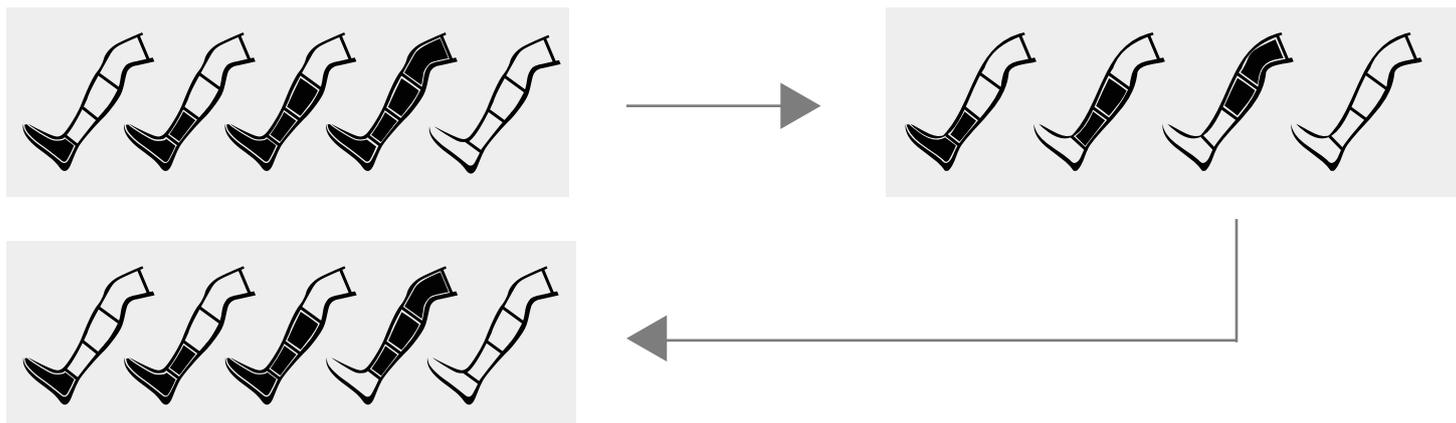- Cyclical pressure on all areas simultaneously.
- Effectively promotes lymphatic drainage.
- Perfect for swelling and quick recovery.

MODE 3: HYBRID MODE

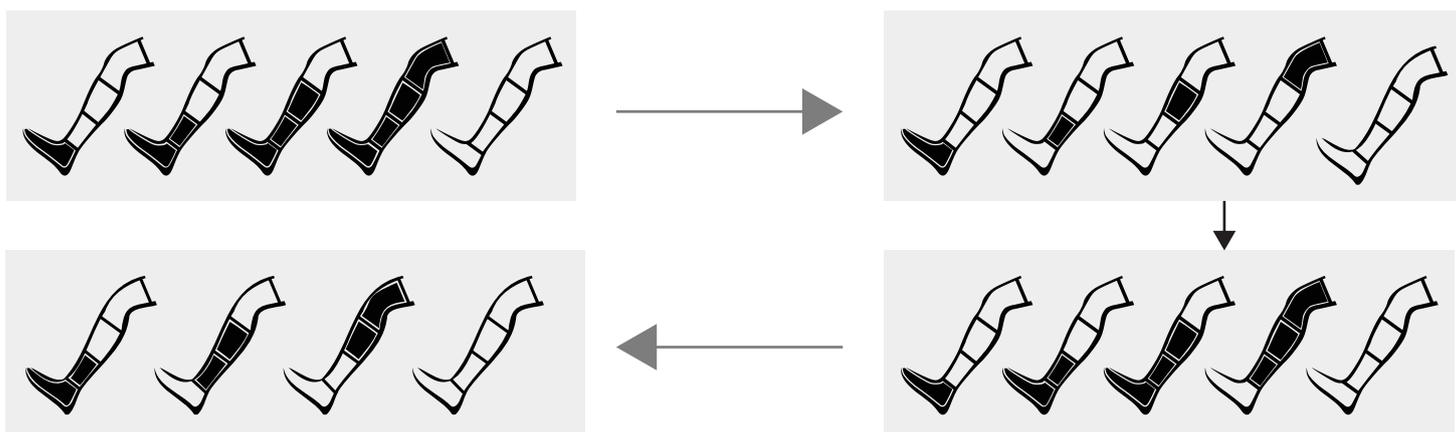FAQS (FREQUENTLY ASKED QUESTIONS)

Q1: How often and for how long should I use the Recovery Boots?

A1: It is recommended to use the Boots 3-4 times per week for 20-30 minutes per application.

The exact duration and frequency will depend on your individual recovery needs. If you are unsure please consult a doctor.

Q2: Can I use the Boots even if I have just exercised?

A2: Yes, the Recovery Boots are ideal after training to promote regeneration and relieve muscle tension.

Q3: Are the boots suitable for all body sizes?

A3: The STARKBOOTS - Recovery Boots are available in different sizes and can be customised to cover different leg circumferences.

Q4: How do I clean the boots after use?

A4: Use a damp cloth to gently wipe the outside of the boots.

Do not immerse the appliance in water and do not use aggressive cleaning agents.

Q5: Can I use the Boots while sleeping?

A5: No, it is not recommended to use the Boots while sleeping. Only use them when you are awake and stick to the recommended duration of use.

Q6: What should I do if the Boots are uncomfortable during use?

A6: If you feel discomfort during use, switch off the device immediately and remove the boots. If necessary, consult a doctor before continuing use.

Q7: What exactly does STARKSURANCE mean for me and how can I convince myself of the quality of the Recovery Boots?

A7: With our STARKSURANCE you have full service for at least 2 years.

If you are not satisfied during this time, we will arrange a replacement, repair or refund.

We want you to be completely satisfied with the STARKBOOTS.

Q8: For which applications are the STARKBOOTS - Recovery Boots suitable?

A8: STARKBOOTS are ideal for runners, strength athletes, cyclists and people with restless legs or other complaints. They promote regeneration, alleviate discomfort and support slim, toned legs through lymphatic drainage.

Q9: How powerful are the STARKBOOTS - Recovery Boots from STARKFELD?

A9: The STARKBOOTS offer adjustable air pressure from 40 to 200 mmHg. This range is ideal for effective treatment. A pressure of up to 160 mmHg is often perceived as optimal and is above average compared to the competition.

Q10: Which mode is suitable for which area of application?

A10: The choice of mode basically depends on the area of application. For complaints such as swelling or to support lymphatic drainage, we recommend the circulation mode, as this is the best way to drain the lymphatic fluid. Sequence mode can also be used, especially if it is uncomfortable to maintain pressure on areas that have already been treated. Hybrid mode is particularly recommended for athletes who want to perform a longer treatment.

Q11: What should I do if a leg no longer inflates?

A11: Check that the hoses are correctly inserted into the control unit. Make sure that no hose is kinked or damaged.

Q12: What should I do if the air supply is weak?

A12: Check whether the hoses are possibly blocked. Clean them carefully or replace them if necessary.

Q13: What should I do if the boots build up pressure unevenly?

A13: Check all connections and hoses for leaks. If the problem persists, contact customer service.

SPECIFICATIONS

Model	SF-002
Name	STARKBOOTS® - Recovery Boots
Adapter	AC Input: 100-240 Volts, AC 50/60Hz, DC Output: 12V 2A
Rated power	24W
Compression modes	3 modes: Sequence, Circulation, Hybrid
Intensity levels / pressure ranges	5332.88 - 26664.4 Pa (40-200 mmHg)
Timer duration	15-60 minutes
Net weight	3.5kg
Gross weight	14.8kg
Boot type	4-chamber boot - foot / lower shin / upper shin / thigh
Available sizes	S, M, L, XL

LEGAL INFORMATION

1. Exclusion of liability:

Starkfeld accepts no liability for damage caused by improper use or failure to observe the operating instructions. Use of the product is at your own risk.

3 Warranty:

We offer you the statutory warranty for all products that you purchase from us. As a rule, the warranty period is 2 years from receipt of the goods. It covers all defects that are attributable to material or manufacturing faults.

3. Starkfeld® trademark rights:

STARKFELD® and STARKSURANCE® are registered trademarks of Starkfeld. Unauthorised use of these trademarks or the associated product names is prohibited and may have legal consequences.

4 Copyright:

All contents of these operating instructions, including texts, images and graphics, are protected by copyright. Reproduction, distribution or modification of the contents is not permitted without the prior written consent of Starkfeld.

5. Data protection:

We collect and process personal data exclusively in accordance with the applicable data protection regulations. For more information, please visit our privacy policy on our website.

6. declaration of conformity:

This product complies with the applicable health and safety requirements in accordance with European standards and bears the CE mark.

7. disposal:

This device contains electronic components and must be disposed of in accordance with the regulations

Electrical and Electronic Equipment Act (ElektroG). Please dispose of the device at the appropriate collection centres or contact us to enable an environmentally friendly return. This product has been certified by TÜV Süd and complies with the applicable European standards, including CE certification. It fulfils all relevant safety and quality standards

CUSTOMER SERVICE - STARKSURANCE® SATISFACTION GUARANTEE

At **Starkfeld**, we want you to be 100% satisfied with your purchase. That's why we offer you the **STARKSURANCE®** satisfaction guarantee. You no longer need a separate warranty card - your satisfaction is our top priority. If you have a problem with your product within the warranty period, simply contact our customer service:

Contact:

E-mail: kundenservice@starkfeld.com

Information: Please provide your order number and a description of the problem. We will take care of it and provide a quick solution - be it through repair, replacement or refund.


With **STARKSURANCE®** we want to offer you security without complicated procedures. Your

TROUBLESHOOTING

Problem	Possible cause	Solution
Device does not work	<ol style="list-style-type: none">1. device is not switched on.2. power supply unit not connected.3. power supply unit defective.	<ol style="list-style-type: none">1. Switch on the device.2. Check the connection of the power supply unit.3. Test the power supply unit or replace it.
Screen on, but boots do not work	<ol style="list-style-type: none">1. a boot is not connected.2. boots are damaged.	<ol style="list-style-type: none">1. connect both boots correctly.2. check the boots for damage.
Air pressure decreases during operation or does not continue to build up	<ol style="list-style-type: none">1. hoses not connected correctly.2. power supply unit defective.3. timer expired.	<ol style="list-style-type: none">1. check the hoses.2. make sure that the power supply is working properly3. start a new session if the timer has expired.


PRODUCT OVERVIEW

Smart Monitor

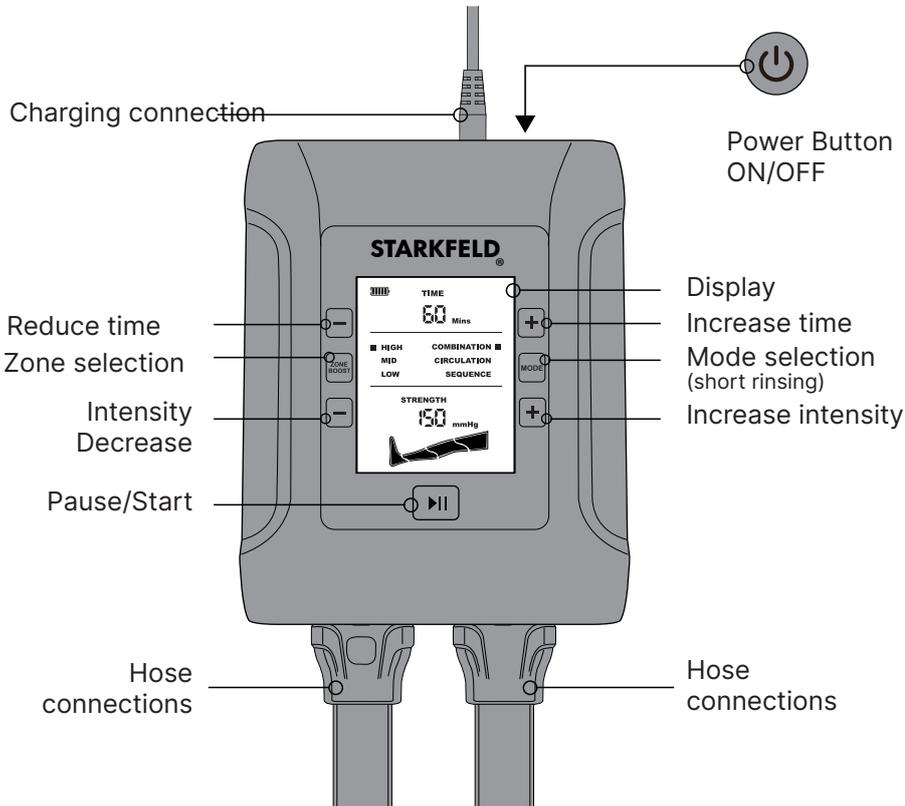
Display

Displays the current settings settings:

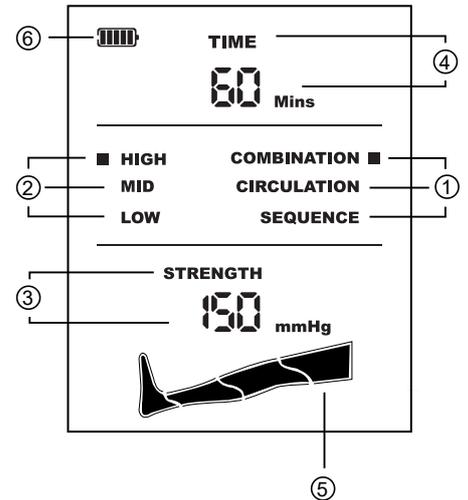
- ① Current mode
- ② Current intensity
- ③ Current pressure
- ④ Remaining time
- ⑤ Pressure chamber
- ⑥ Battery charge status

Boots

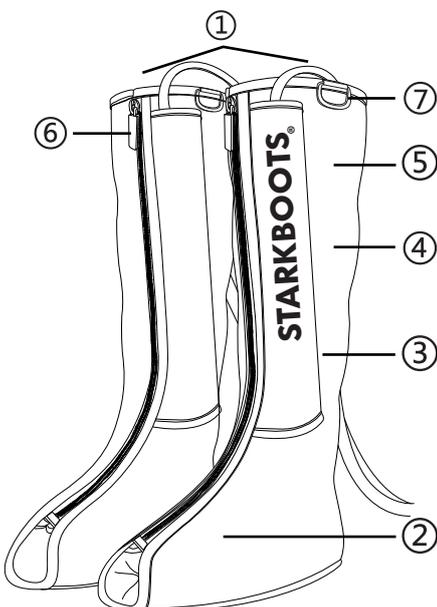- ① Hose connection of the boat
- ② Foot air chamber
- ③ Lower calf chamber
- ④ Upper calf chamber

Power supply unit

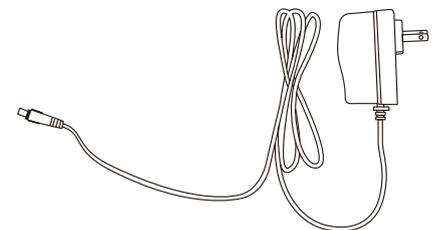- ⑤ Thigh chamber
- ⑥ zip fastener
- ⑦ D-buckle


OPERATING INSTRUCTIONS

Please read the safety instructions carefully before using the product.

1. Connect the connections:


Attach the boat connections firmly to the hoses.

Note: Both boots must be connected.

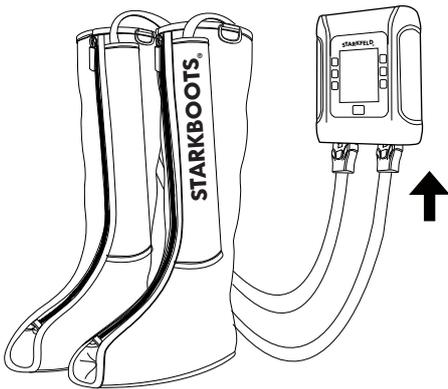
2. Connect the hoses to the Smart Monitor:

Attach the boat connections firmly to the hoses.

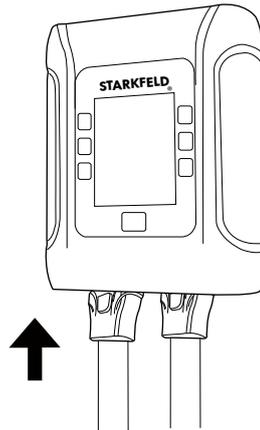
3. Connect the power supply unit:

Connect the adapter to the socket and the Smart Monitor. The display shows that the battery is charging.

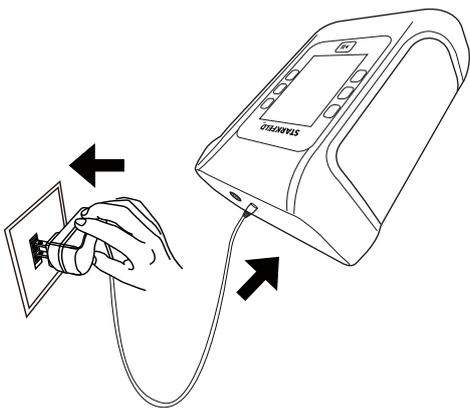4. Put on a boot:

Put the boots on correctly and make sure they fit snugly and comfortably.

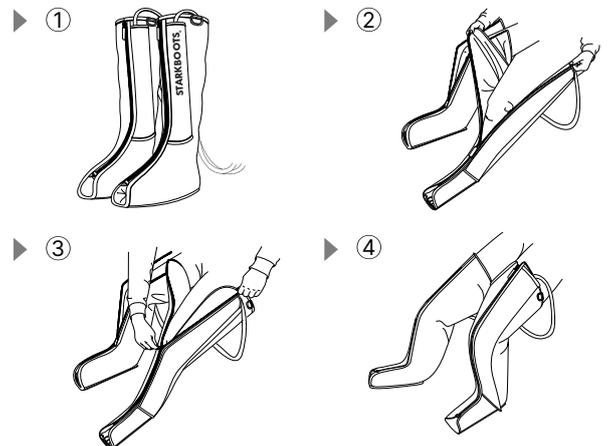TIPS FOR USE:

• Close the zip:


Make sure that the boot zip is fully closed before starting the system. Never use the system if the zip is partially or fully open.

- **Use of the D buckle:**

Use the D-buckle to make it easier to put on the boots. This allows you to put the boots on more quickly and easily.


5. Get into position and start:

Find a comfortable position in a sitting, lying or half-sitting position. Press the on/off switch on the Smart Monitor to start the session.

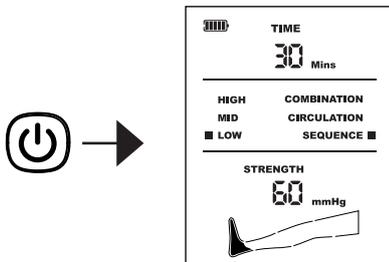6. Current settings: Selection.

After starting, a signal tone sounds and the screen lights up to display the current settings. By default, the device starts in sequence mode with an intensity of 60 mmHg and a duration of 30 minutes.

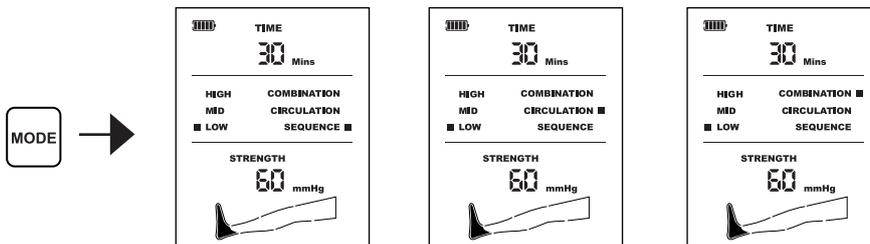7. Switch mode:

Press the mode button to switch between the three modes (sequence, circulation, hybrid).

8. Select massage zone:


Press the zone button to select a specific massage zone if you only want to treat one area.

9. Adjust intensity:

Use the intensity buttons (+/-) to adjust the pressure between 40 mmHg and 200 mmHg. Start with a lower intensity for the first application and increase it gradually.

10. Set time:


Use the time buttons (+/-) to select the treatment duration between 15 and 60 minutes. The time can be adjusted at any time during the session.

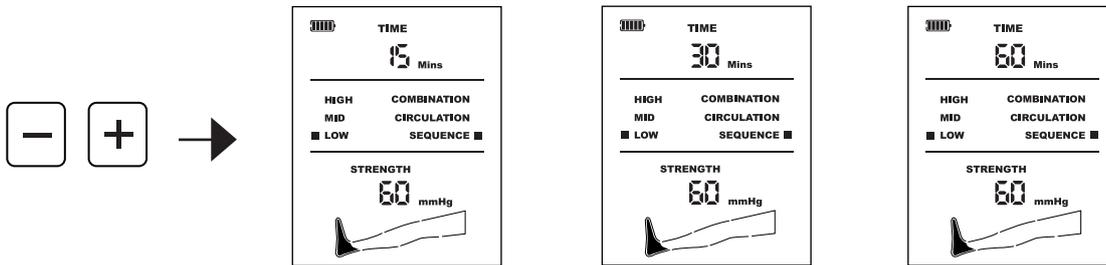11. Automatic switch-off:

The device switches off automatically after the set time has elapsed. If you want to extend the treatment, you can restart the session manually.

NOTES AFTER USE

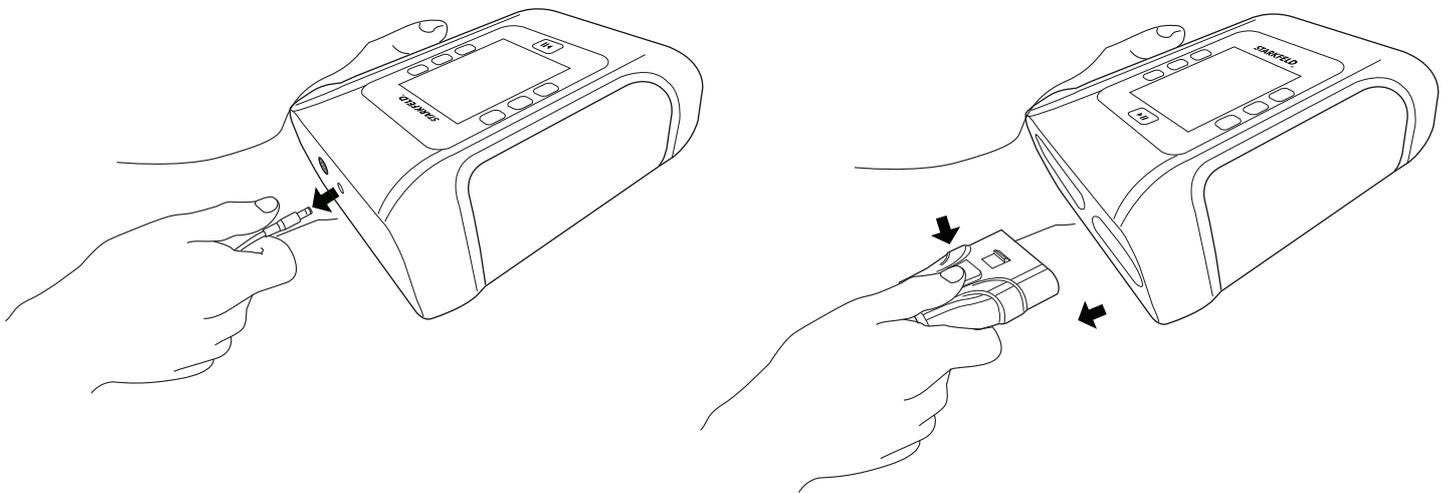1. Switch off the Smart Monitor:

After use, switch off the Smart Monitor by pressing the on/off switch. Wait until the display has completely gone out before continuing. Disconnect the hose from the control unit by pressing the button on top of it and pushing the button back.

3. Disconnect the power supply unit:

Carefully unplug the power adapter from the Smart Monitor and then from the socket. Take care to pull the plug out directly so as not to damage the cable.

4. Remove the hoses:

Press the hose connection button on the Smart Monitor and carefully pull off the hoses. Check the hoses for any kinks or blockages before stowing them away.

4. Packing and storage:

Store all accessories, including the boots, hoses and Smart Monitor, in the original packaging or storage bag. Store the device in a dry and clean place to protect it from dust and damage.

5. Maintenance and care:

After each use, check that the Smart Monitor, boots and hoses are in perfect condition. Look for cracks, dirt or loose parts and clean the device according to the cleaning instructions if necessary.

CLEANING

1. Switch off the appliance before cleaning.
2. Wipe the control unit, boots and hoses with a damp, clean cloth. 3 Dry all parts thoroughly with a clean cloth.
3. dry all parts thoroughly with a clean cloth.
4. Do not use aggressive cleaning agents (such as petrol, alcohol or thinner) to avoid damage.
5. prevent foreign bodies from entering the hoses.
6. do not use chemical cleaning agents.

STORAGE & MAINTENANCE

1. Do not disassemble the appliance yourself. If you have any problems, contact customer service.
2. keep the appliance out of the reach of children.
- 3 Do not store the device in extreme temperatures or humid conditions.
4. avoid direct sunlight.
- 5 Do not insert needles into the boots or hoses.
6. do not place any heavy objects on the device.
7. keep the connections of the control unit and hoses clean.
8. regularly remove dust or fluff from the control unit and hoses.

WASTE DISPOSAL

The STARKBOOTS® - Recovery Boots contain electronic components and batteries that must be recycled. Do not dispose of the boots in household waste. In accordance with the Electrical and Electronic Equipment Act (ElektroG) and the Battery Act (BattG), such devices must be disposed of in an environmentally friendly manner in order to recycle valuable materials and keep harmful substances out of the environment.

1. Recycling:

The boots contain valuable materials such as copper and rare earths as well as batteries. These must be disposed of separately and recycled. This helps to conserve resources and safely dispose of hazardous substances such as lead and mercury.


2. Return to collection centres:

Electrical appliances and batteries must not be disposed of with household waste. Return the boots and the batteries they contain to municipal collection points (e.g. recycling centres) or to shops that take back electrical appliances and batteries (e.g. specialist electrical shops or supermarkets).

3. Return to the manufacturer:

You can also return the boots and batteries directly to the manufacturer, who will take care of their proper disposal.

Please observe the local disposal regulations and help to protect the environment by recycling correctly.

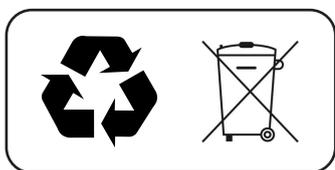SCOPE OF DELIVERY

- 2x recovery boots
- 1x Smart Monitor
- 1x Power adapter (AC100-240V, 50/60Hz; DC12V 2A)
- 1x User manual
- 1x Portable storage bags