
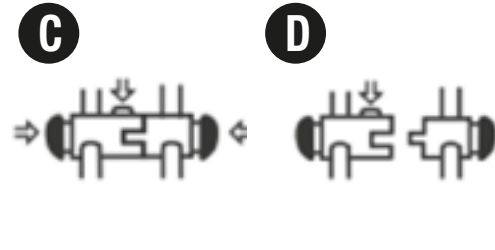
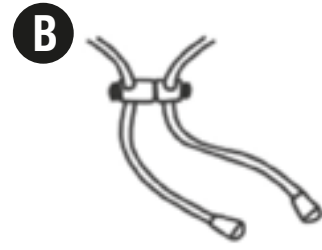
# WELLPRIME® Miracle Massager
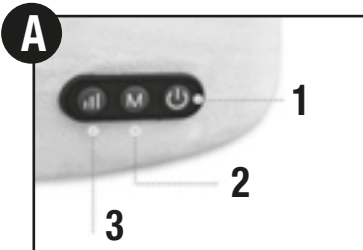
## ENGLISH

Thank you for purchasing Wellprime® Miracle Massager, the ergonomically designed neck massager that adapts to your neck and keeps your head in a comfortable position. Please read the operating instructions and all safety instructions carefully before use. Keep the manual in an accessible place for future use.

### SAFETY WARNINGS

- Do not place the Wellprime® Miracle Massager in the bathroom or in a damp place.
- Always keep the Wellprime® Miracle Massager dry. Do not store it in damp places or near containers of water into which it may fall.
- It is not waterproof; if the unit should fall into water, do not touch or reuse it. Do not submerge the unit under water.
- Children under 3 years of age should be kept away unless continuously supervised.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children must not play with the appliance. Cleaning and maintenance shall not be made by children without supervision.
- Do not leave the massager exposed to direct sunlight or high temperatures for a long period of time.
- This product has a built-in rechargeable battery. Please dispose of this product in accordance with the applicable environmental protection regulations.
- Never charge the battery if the temperature is below 0°C. Use, charge and store the device at a temperature between 5°C and 35°C.
- Avoid charging the battery by exposing it to direct sunlight or near a heat source.
- Never charge the battery continuously for more than 48 hours.
- To recharge the battery, use only the detachable power supply unit supplied with this apparatus.
- The appliance may only be operated with the power supply unit supplied.
- This product is not a toy.
- Do not use this product if it is not working properly and the massage device is damaged or defective.
- In case of severe or intense pain during the massage, stop the massage immediately and seek medical advice.
- Do not use this product directly on swollen, inflamed or inflamed areas or rashes.
- Do not charge the product when the power supply is unstable.
- Do not touch the power supply and plug with wet hands to avoid malfunction or electric shock.

### COMPONENTS

- FIG. A.  
1- On/Off button  
2- M: Modes (change mode/function)  
3- Speed

### INSTRUCTIONS FOR USE

1. Press the on/off button for 2 seconds and the massager will start working. Pressing the power button by default will set the massager to the function: kneading + heating + vibration.

2. Press the "M" button to change mode: By pressing M 1 time: "click 1" you will get the mode: heat +

### Vibration

If you press it 2 times "click 2": it will switch to heat mode. If you press it 3 times "click 3"- it will change to: kneading massage + warmth  
And if you press it 4 times "click 4-: you will get only the kneading massage function.

3. Press the "III" button speed:  to change speed. Wellprime® Miracle Massager has 3 speeds: low, medium and high. (FIG. A)  
4. When you stop using it, press and hold the power button for 2 seconds to turn it off. (FIG. A)

Each programme has a preset duration of 14 minutes; after 14 minutes the massager will switch off.

The massager is powered by a built-in lithium battery. If the battery is about to run out of power or has run out. Charge the massager before use. When charging, use the charging cable of this product to check if the power connection is normal and plug it into the DC interface to charge. Do not use the product while it is charging.

The Wellprime® Miracle Massager takes about 2-3 hours to charge. This product includes accessories only for use with this machine. Other accessories on the market are not a substitute for this product.

(FIG. B). The Wellprime® Miracle Massager has a fastener to fit around the neck.  
(FIG. C). Press the unlock button  
(FIG. D). Press up and down to loosen

Composition:  
Fabric: 100% polyester  
Filling: Viscoelastic foam  
Pillowcase Cleaning and Maintenance:  
The pillowcase, which is soft and plush, can be easily removed for washing. Unzip the pillowcase. Remove the foam filling and wash the pillowcase according to the care and washing instructions on the Wellprime® Miracle Massager composition and care label. NOTE: The memory foam cannot be washed. Only the pillowcase can be washed.

TECHNICAL SPECIFICATIONS  
Power rating: 3.7w  
Input: 5V=1A  
Battery: 3.7V= 1200 mAh 4.44Wh

QUALITY WARRANTY:  
This product is covered by a warranty against manufacturing defects subject to the time periods stipulated by the legislation in force in each country. This warranty does not cover damage resulting from inadequate use, negligent commercial use, abnormal wear and tear, accidents or improper handling.

Australia & New Zealand only  
This product has the benefit of certain Consumer Guarantees. These are prescribed by the Australian Consumer Law & the New Zealand Consumer Guarantees Act 1993 both of which provide protection for consumers. There is no express warranty for this product in Australia or New Zealand. The above paragraph refers to other countries.

 Disposal of old electric appliances  
The European directive 2002/96/EC on Waste Electrical and Electronic Equipment (WEEE), requires that old household electrical appliances must not be disposed of in the normal unsorted municipal waste stream.

Old appliances must be collected separately in order to optimize the recovery and recycling of the materials they contain, and reduce the impact on human health and the environment. The crossed out "wheeled bin" symbol on the product reminds you of your obligation, that when you dispose of the appliance, it must be separately collected. Consumers should contact their local authority or retailer for information concerning the correct disposal of their old appliance.

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