
A I L O R I A

DEPLAQUED

DE Bedienungsanleitung **EN** User manual **FR** Manuel de l'utilisateur **IT** Manuale di utilizzo **ES** Manual de Usuario **NL** Gebruiksaanwijzing **PL** Instrukcja obsługi **RO** Manual de utilizare **PT** Manual do utilizador **SV** Bruksanvisning **DA** Brugervejledning **CS** Uživatelská příručka **NO** Bruksanvisning **SK** Používateľská príručka **FI** Käyttöohje **LT** Naudotojo vadovas **BG** Ръководство за потребителя **HR** Korisnički priručnik

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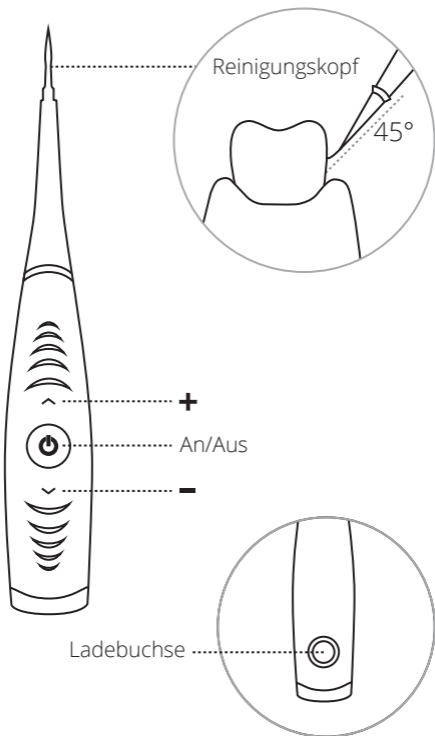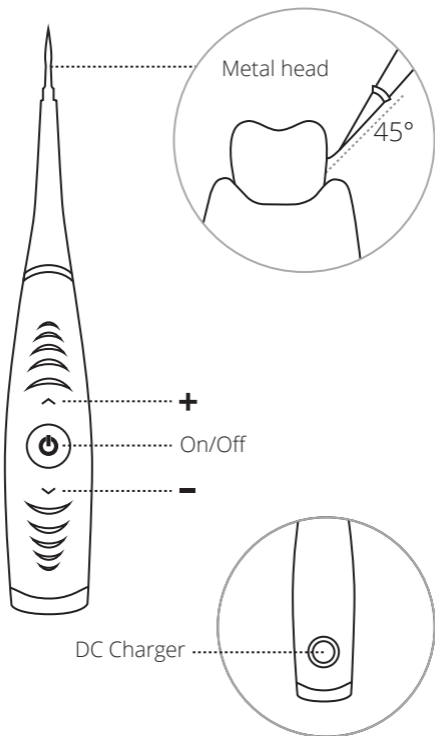
PRODUCT FEATURES

The AILORIA DEPLAQUED sonic scaler is ideal for cleaning teeth at home. It removes tartar, plaque and discolouration, even from hard-to-reach areas.

- Intelligent high-frequency vibration
- Ergonomic design for comfortable handling
- 5 intensity levels
- Cordless battery operation
- Charging time approx. 2 hours
- Battery life approx. 6-8 hours when fully charged
- Charging of the lithium-ion battery via USB

COMMISSIONING

Charge the sonic scaler before the first use. Charging time is approximately 2 hours to charge the sonic scaler before first use or fully discharged battery.

CHARGING PROCEDURE

The connector for the charging cable is located at the bottom of the back of the unit. Open the silicone cap and connect the charging cable to the sonic scaler. Connect the other end to a power source (wall socket) using the USB charging adapter provided. Alternatively, you can also connect the USB cable to a USB interface (e.g. PC or notebook).

The sonic scaler cannot be used during the charging process. A red light shines in the transparent base for the duration of the charging process. When the charging process is completed, this is indicated by a green light.

APPLICATION

Press the ON/OFF button to switch on the sonic scaler. Select the desired intensity using the arrow keys.

The following 5 intensity levels are available:

Intensity level 1: 8321 vibrations per minute

Intensity level 2: 9787 vibrations per minute

Intensity level 3: 11195 vibrations per minute

Intensity level 4: 12423 vibrations per minute

Intensity level 5: 13320 vibrations per minute

Start cleaning the teeth by gently placing the metal attachment on the areas with plaque and discoloration. Now slowly move the metal tip over the areas to be cleaned. This process can be repeated several times until the desired result is achieved. The sonic scaler can also be used on the edges of the teeth and in the spaces between the teeth. Move the sonic scaler carefully to prevent possible injuries to the gums. Repeat the application 1 to 2 times per week. The DEPLAQUED can also be used on dental crowns and bridges. Proceed carefully.

Clean the metal head with water after each use and wipe the housing with a damp cloth. Sterilize the sonic scaler with alcohol if necessary.

SAFETY INSTRUCTIONS

- Always keep the USB cable and USB adapter away from water and other liquids.
- Only use the USB cable provided with this product. Using alternative cables may lead to damage.

- Avoid excessive bending of the USB cable to maintain its integrity and charging capability.
- Ensure the power plug used with the product's USB cable is appropriate and compatible.
- Discontinue use immediately if you observe any product malfunctions or signs of potential fire risk.
- Remove the USB charging cable before cleaning the ultrasonic scaler.
- Do not use the product if it is damaged. The sonic scaler can only be repaired by the manufacturer.
- Children and persons with reduced physical or mental capabilities or lack of experience and knowledge should only use the product under adult supervision and instruction. Always keep DEPLAQUED out of the reach of children.
- DEPLAQUED is only suitable for cleaning teeth and should not be used for any other purpose.
- If excessive bleeding of the gums or pain occurs and persists while using the sonic scaler, please consult your dentist.
- The product should not affect pacemakers in any way. However, as a maximum precaution, we recommend that you consult your doctor before use.
- If you have undergone jaw or gum surgery in the last two months or if you have any health problems, consult your doctor before using DEPLAQUED.
- Clean the DEPLAQUED metal attachment with

water after each use and wipe the housing with a damp cloth. Sterilize the tooth cleaner with alcohol if necessary.

- The sonic scaler is protected against water splashing according to IPX6 protection standard. However, do not keep the product under water for a long period of time.
- Avoid overcharging the battery; disconnect the device from the charger once it is fully charged.
- Store the device in a cool, dry place and steer clear of extreme temperatures.
- Do not expose the battery to direct sunlight and do not store the device in damp environments.
- Protect the device from impacts, drops, and excessive pressure to prevent damage to the battery.
- Do not attempt to remove or replace the battery yourself.
- If you notice any signs of damage or leakage from the battery, stop using the device immediately and contact customer service.
- Avoid contact with leaking fluids; if skin contact occurs, wash thoroughly with water right away.
- Dispose of the device at the end of its life cycle according to local regulations for electronic waste.
- Do not throw the battery into fire or water.

- If the device overheats, emits smoke, or shows any other signs of malfunction, turn it off immediately and disconnect it from any power source.
- Keep the device away from flammable materials until the issue is resolved.

TECHNICAL DATA

Rated voltage: 3.7V

Rated current: 200mA

Battery capacity: 300mAh

Charging voltage: 5V


DISPOSAL INFORMATION

Waste electrical and electronic equipment must not be disposed of with household waste. You must therefore not simply throw away your old appliances. They must be recycled. In order to improve waste separation, the Electrical and Electronic Equipment Act therefore stipulates that all electrical and electronic equipment must be permanently marked with the corresponding symbol, the crossed-out dustbin. Free collection points for old electrical appliances are available for you to return them. You can obtain the addresses from your city or municipal administration. Batteries and accumulators that are not enclosed in the old appliance, as well as lamps that can be removed from the old appliance without being destroyed, must be removed from the old appliances before disposal. As a consumer, you are generally obliged to hand in used batteries and accumulators as well as lamps separately at the special collection points. Make sure that the batteries and accu-

mulators are completely empty. Then carefully remove the batteries and accumulators from the old appliance. The following batteries are included with this item: 1 rechargeable battery, type Li-Ion, 3.7V, 300mAh.

Subject to alterations and errors

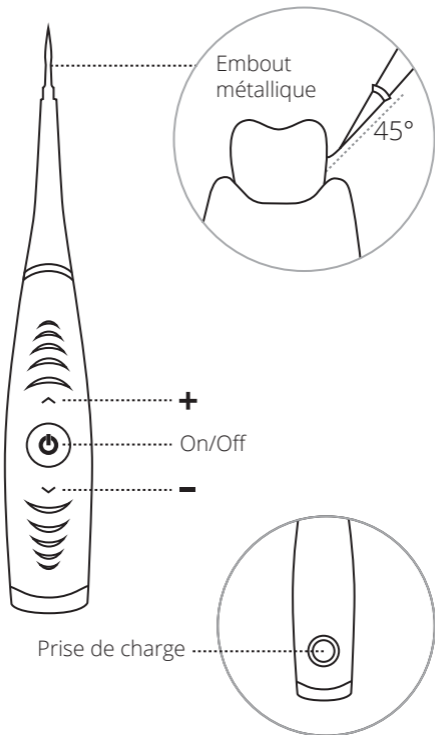The manufacturer accepts no liability for damage caused by incorrect, improper or incompatible use of the product.

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© VEHNS GROUP GmbH
DONG GUAN YI YUAN DA ENERGY TECHNOLOGY CO. , LTD
VEHNS GROUP GmbH Theatinerstr. 40-42, 80333 Munich,
VEHNS GROUP GmbH, Theatinerstrasse 40-42