


A302-9 MASSAGE CHAIR

Operation Instructions

SAVE THESE INSTRUCTIONS

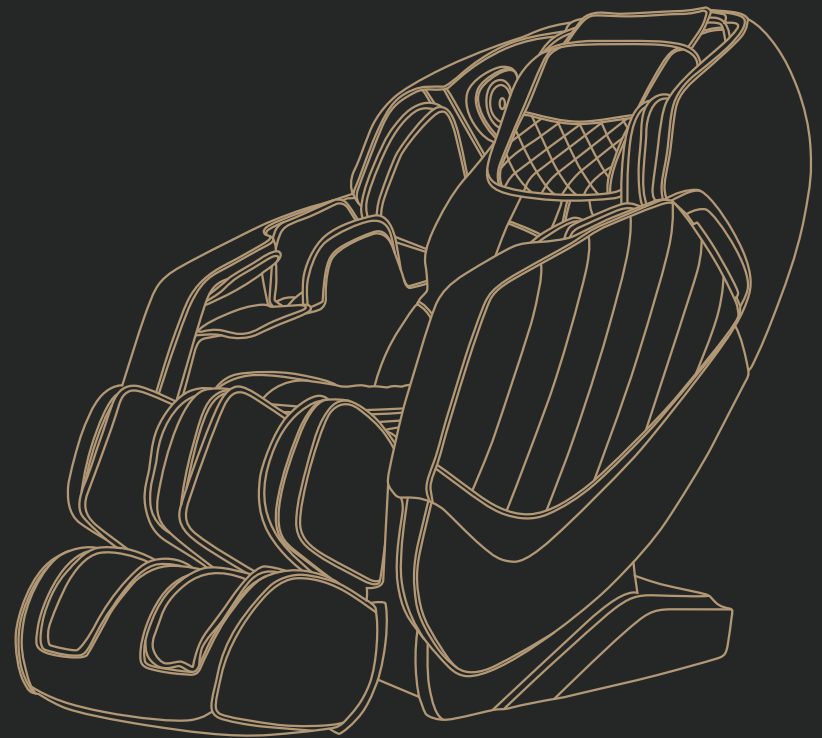

Thank for your purchasing this product. Before your use it ,Please read the Operation instructions carefully for the correct usage of this product , especially the important safety warnings. Keep the operation instructions properly for future reference.

Note: Our company reserves the right to change and explain the design of the product. If there is some difference between the product function in market circulation and the description in the instructions, it is a normal update of the product without prior notice. The specific product function and color are subject to the actual product.

CONTENTS

Before Use

Safety Precautions ----- 02-04
Installation Parts ----- 05
Place Massage Chair ----- 06-07
Before Use ----- 08-09
Product Structure ----- 10-13

Usage Method

Operation Instructions ----- 14-17

Maintain



Cleaning and Maintenance ----- 18

Important

Troubleshooting ----- 19
Product Specifications ----- 20

Safety Precautions must be followed


- Please read all relevant instructions carefully before use.
- The precautions indicated here will help you use this product safely and correctly.
- In order to clearly indicate the degree of harm and damage, the precautions are divided into "Warning" and "Caution" that may be caused by incorrect use. Both of them are important for safety and must be followed.

 Warning	Indicates that serious injury or death may result from improper use.
 Caution	Indicates that misuse may cause minor injury or property damage.

- Please pay attention to the symbols in front of each safety precaution.


 Forbidden	Absolutely Prohibition.
--	-------------------------

Following people should not use massage chair


 Forbidden	<ul style="list-style-type: none"> • The device is not intended for use by persons (including children) with physical disabilities, sensory or neurological impairments or who lack experience and common sense, unless they are supervised and instructed in the use of the device by those responsible for their safety. Children must be supervised to ensure they do not play with instruments. • Osteoporosis patients. • Patients with heart disease or wearing electronic medical devices such as a pacemaker. • People with fever. • Pregnant or menstruating women. • People who are injured or have a skin disease. • When unsupervised, any children are under the age of 14 and those unconscious are not allowed to use this product. • People who are told by their physicians to have rest or who feel unwell. • People whose bodies are wet are forbidden to use this product. • This product has a hot surface. Persons insensitive to heat must pay attention when using this product. • If the power cord is damaged, in order to avoid danger, it must be replaced by a professional from the manufacturer, its maintenance department or similar department.
--	--

Safety Precautions must be followed


Important safety warnings

 Warnings	<ul style="list-style-type: none"> • Keep children away from movable parts of this product. • Use the well-grounded power supply unit suitable for this product. • Pull out the plug after use or before cleaning to avoid any physical injury or damage to the product. • Use this appliance as instructed in the manual. • Do not use any parts or accessories which are not recommended. • Do not use this product outdoors. • Please read the Operation Instructions carefully before operation. • Never use this product for any other purpose not listed herein. • 20 minutes of usage each time is recommended. • Do not use this product if the cover or leather is broken or damaged. • Do not drop anything into this product. • Never fall asleep while using this product. • Do not use this product if you are drunk or feel unwell. • Do not use this product within one hour after a meal. • Do not use this product with excessive force to avoid any injury.
---	--

Safety precautions

 Caution	<ul style="list-style-type: none"> • Check the voltage if meets the specifications for this appliance. • Do not plug or unplug the power plug with wet hands. • Keep dry, do not let water flow into this appliance to avoid electric shock or cause damage. • Do not pull the power cord when plugging or unplugging it. Rough handling is forbidden. • Do not damage the wire or modify the circuit of this product. • Never use a wet cloth to clean live parts such as the power switch, plug, etc. • Be away from this product in case of power failure to avoid any injury if power is restored suddenly. • Stop using this product at once while it works abnormal, and consult your local supplier. • Stop using this product if you feel unwell, and consult your health care practitioner.
--	---

Storage conditions

 Caution	<ul style="list-style-type: none"> • Storage room temperature: 5 °C to 40 °C. • Relative humidity (RH) between 20% and 80%; free of corrosive gases, and a well-ventilated room.
--	--

Safety Precautions must be followed

Product Care and Maintenance



Caution

- This appliance may only be repaired by a service center designated by the manufacturer. Users cannot disassemble or repair it by themselves.
- Remember to turn off the main power switch after use.
- Do not use this product if the power socket becomes loose.
- If this product will be left unused for a long period of time, it is necessary to roll up the power cord of this product and store this product in a dry and dust-free environment.
- Do not store this product at a high temperature or near open flame. Avoid prolonged exposure to direct sunshine.
- Please clean this product with dry cloth. Never use thinner, benzene or alcohol.
- The mechanism of this appliance were specially designed and manufactured and don't require special maintenance.
- Do not use a sharp object to stab at this product.
- Do not roll or pull this product on any uneven ground. Instead, it should be lifted before moving.
- Please use this product intermittently. Do not keep it working continuously for a long period of time.

Operating environment



Caution

- Do not use this product in a very moist environment such as the bathroom.
- Do not use this product in any heavily-dusty or caustic environment.
- Do not use this product in a small space or a place without good ventilation.

Solutions to ordinary malfunctions







Caution

- If the motor sounds during operation of this product, this is the normal operation sound.
- If the controller cannot work normally, check whether the plug and the power socket are connected reliably, and whether the power switch is turned on.
- If the rated working time is over, this appliance will automatically activate the power-off button; If this product keeps operating for a long period of time, the temperature protector will stop this product automatically. The product can be used again after half an hour.

Installation parts

Calf installation instructions

1. Insert the wire of the calf and the air tube into the front socket of the steel frame assembly (Figure 1).
2. Then push the connected wire and air tube into the hole in the steel frame (Figure 2)
3. After the wire docking is completed, open the calf fixing covers on both sides of the steel frame assembly (Figure 3).
4. Lift the calf into the slot, close the fixing covers, and secure it in place with 2 pieces of screws (Figure 4).

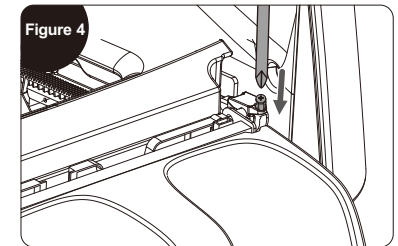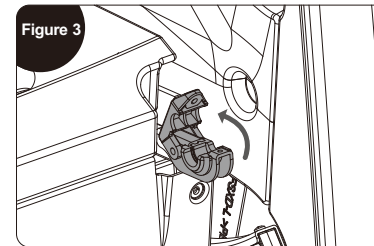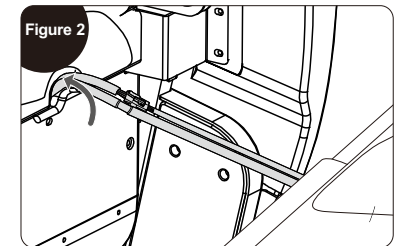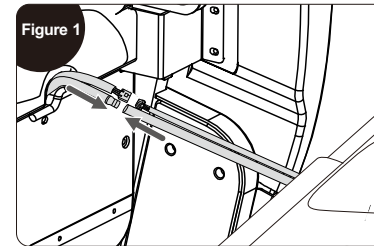


Place massage chair

Installation placement

Make sure there is enough space for the massage chair can be tilted at least 15 cm.



Do not expose the massage chair directly to high temperature environment due to sunlight (such as: in front of heating equipment), and place it to where can avoid high temperature and direct sunlight. It is recommended to lay a soft cushion for long-term use in a designated location. When laying the cushion, the size of the cushion should be enough to cover the range where the machine touches the floor and the range where your feet touch the floor.

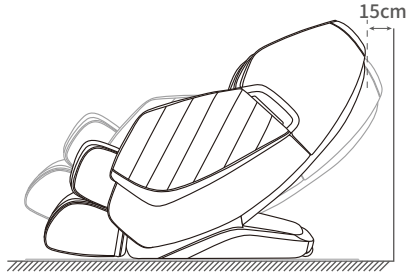


Product lift instruction

Require team work to lift and move the chair to avoid injury.



1. Please lift the massage chair to move it on vulnerable ground such as wooden floor.
2. When moving on a similar wooden floor, it is recommended to lift it. To avoid damage to the casters. When the floor is lifted and moved to the designated position, be careful to keep your hands out of the way to avoid injury, and be careful of your feet being crushed. Do not let go of the unit until the unit is completely flat, the floor may be damaged, therefore, it is recommended to place mats, etc. on the floor.

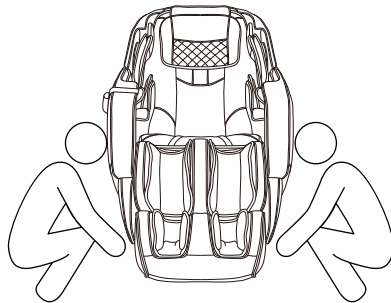

Place massage chair


Move method

Use the rear casters to move the massage chair. Move the chair by tilting it down and backward about 45° as shown in the diagram.

NOTE: Excessive force to title will cause the chair to tip over completely and product damage.






1. Before moving, please shut down the machine properly to reset the whole machine, and then cut off the power.
2. Please keep the chair vacant when moving. (People, pets or clutter will not be able to stay on the chair.)
3. When moving, do not cross a drop barrier higher than 2mm or a gap greater than 5mm.
4. When using casters to move this product, single continuous movement is not supported; if the distance exceeds 50 meters, please make this product move slowly and evenly.

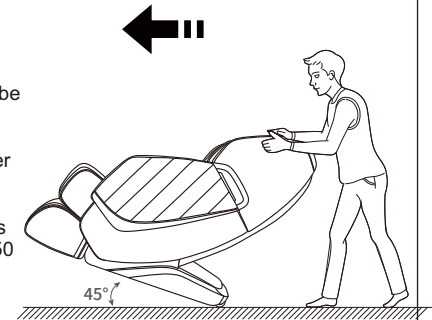

How to disassemble and use headrest cushion and backrest cushion


Using a headrest cushion can reduce the intensity of kneading massage on the neck and shoulders, and you can decide whether to use a headrest cushion according to your needs (recommended to use a headrest cushion). The backrest cushion and the backrest are connected by a zipper (1), and the headrest cushion and the backrest cushion are connected by a zipper (2).

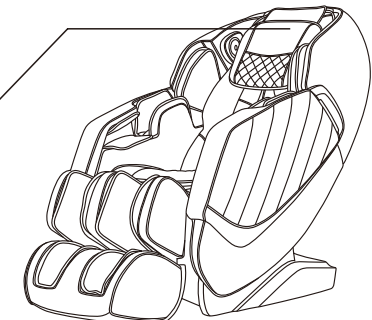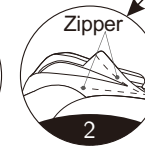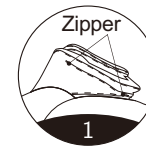

Before use


Check around

Make sure that there are no people, pets and other objects around the unit.

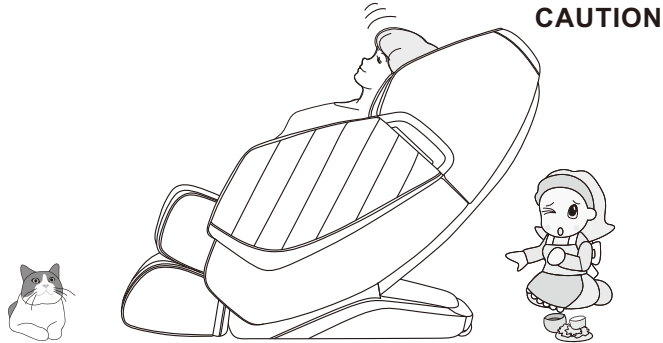

Make sure there are no foreign objects in the gaps of the machine


Before massaging, make sure that no foreign objects are caught in the legs, feet, etc.
Under no circumstances allow hands, feet, and head to penetrate into the area between the lower leg mechanism and main body, or the area between the armrest and the capsule.
Under no circumstances should hands, feet, or head sink under the fabric or leather covering of the massage area.

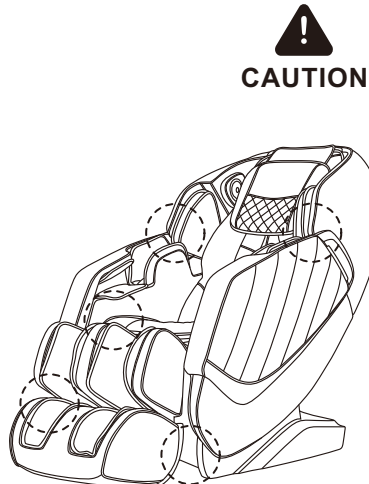


Before Use

Turn on the power


CAUTION

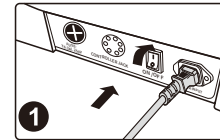

Plug in power cord

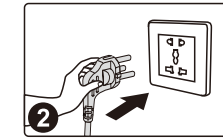

Schematic diagram of power on

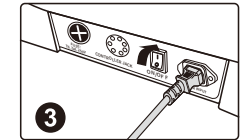

Turn on the power switch ("I" position)




Warn

- Please check whether the power cord and power plug are damaged before use.
- Please supply power to this product in accordance with the working power requirements of this product's calibration.
- It is strictly forbidden to use a power supply other than the nominal power supply of this product.
- It is strictly forbidden to connect various power conversion devices to supply power to this product without authorization.


Product Structure

Description of Components

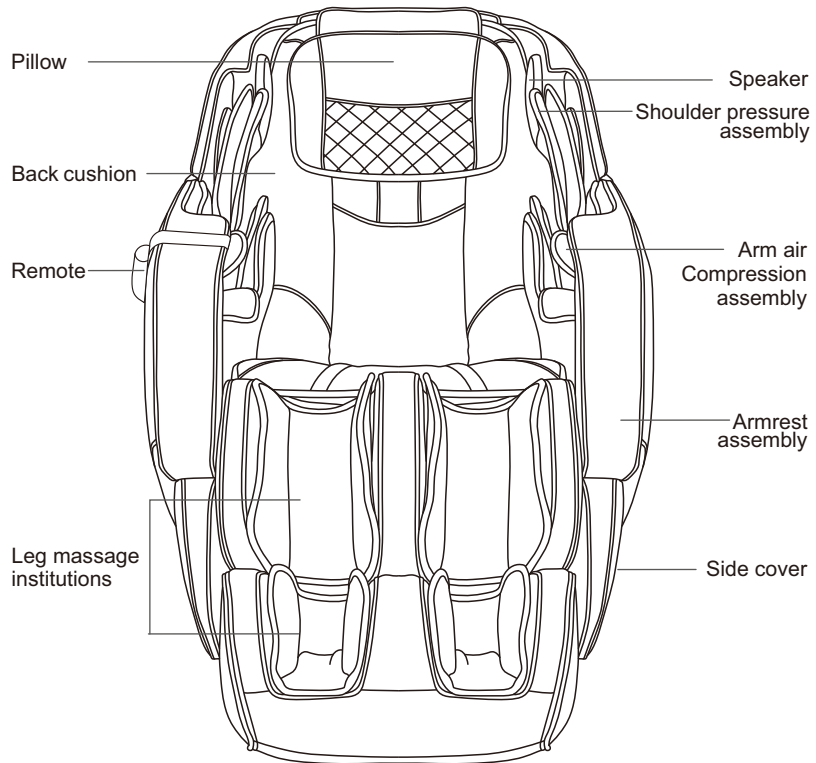


Product Structure

Function distribution

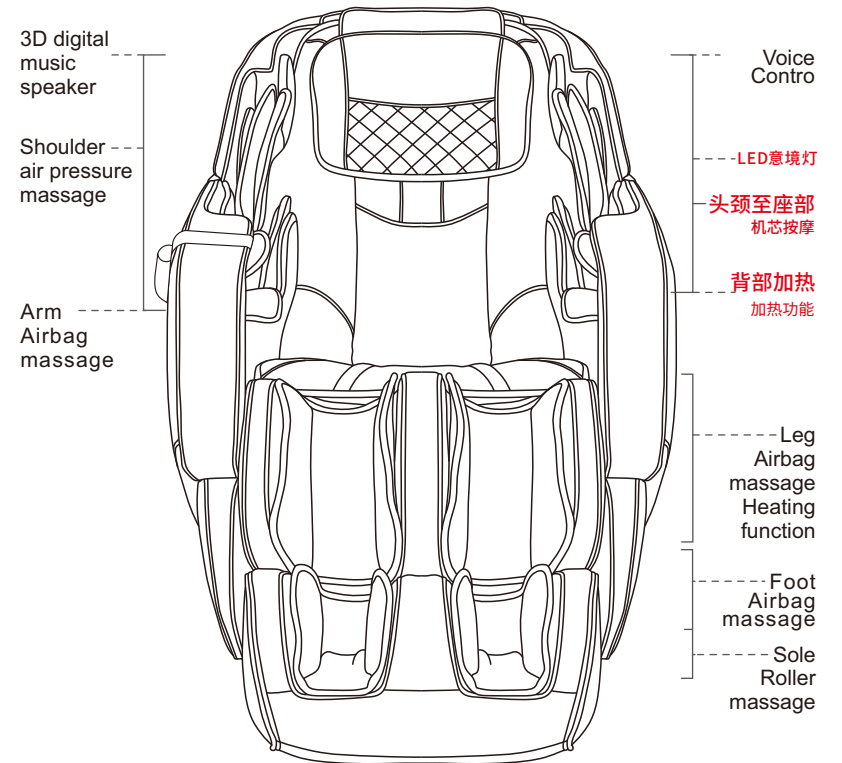

Product Structure

I .Introduction of the remote controller

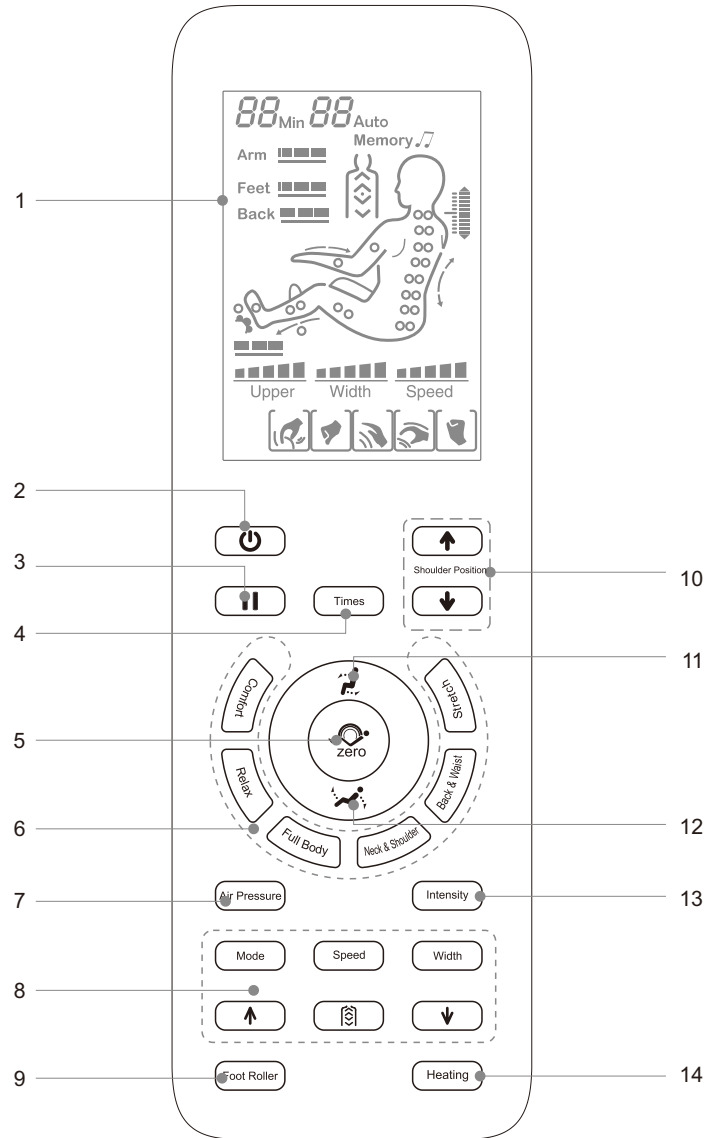

Product Structure





- 1.LCD display:
 - Display all massage functions
- 2.Power switch:
 - Turn on or Turn off all massage functions.
- 3.Pause button:
 - Stop or start all massage functions..
- 4.Timing button
 - Time setting
- 5.Auto lying Button
 - There are three automatic lying position, choose your favorite one.
- 6.Auto massage Button
 - Please refer to page 10 of this manual for details.
- 7.Air pressure mode Button
 - Turn on/off the air pressure massage function (refer to page 11 of this manual for details).
- 8.Manual massage button
 - Please refer to page 10 of this manual for details.
- 9.Sole roller button
 - Turn on/off sole roller massage function
- 10.Shoulder adjust button
 - Make the adjustment on shoulder acupressure points.
- 11.Lifting button
 - Rise up the backrest and lower down the footrest.
- 12.Lying button
 - Lower down the backrest and rise up the footrest.
- 13.Air pressure intensity button
 - Make the adjustment on air pressure intensity.
- 14.Heating
 - Turn on/off heating massage function


Operation Instructions

II. Start massage

- ①.Press the power button to start the massage, the default is relax massage.
- ②.Automatic angle adjustment (you can press any key to skip)
- ③.Automatic human shoulder position detection, keep the sitting posture fit with the massage chair (you can press any key to skip).
- ④.Shoulder position fine-tuning, if the detected position deviates from the actual shoulder position, adjust it to the appropriate position through the up and down key of the shoulder position (or you can press any key to skip), there are 11 shoulder positions in total. If there is no adjustment within 10 seconds, it will enter into the massage function.

Note: Automatic function can be selected during startup. Repeated selection of the same automatic function can skip the startup process and directly massage.



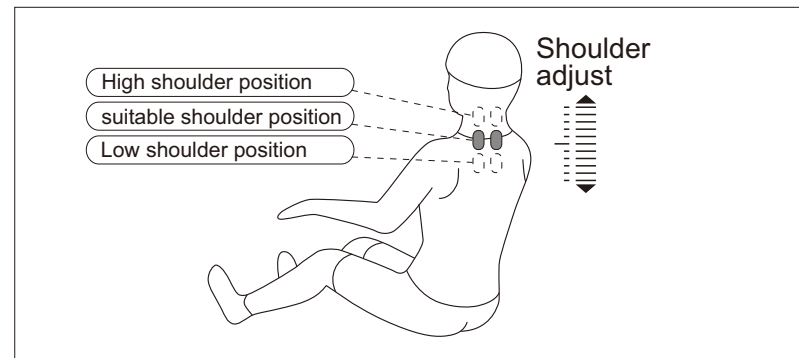



Button	Description	Display
	Move the massage balls upward to a suitable shoulder position.	
	Move the massage balls downward to a suitable shoulder position.	







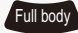













During the detection of shoulder position and massage finger pressure point (the hand controller screen displays the shoulder position adjustment icon and hears the "Didi" of the hand controller. At the same time, press the up and down buttons of the shoulder position to adjust the shoulder position).

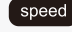Operation Instructions

· Features automatics massage




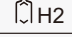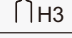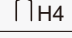




Function	Button	Description	Display
Features automatics massage		Loosen the muscles.	F1 Auto
		Deep-tissue massage throughout key acupressure points to relax the sore area.	F2 Auto
		Full body massage, relaxing the sore areas.	F3 Auto
		Focus on the neck and shoulder area	F4 Auto
		Focus on the back waist area.	F5 Auto
		The footrest and backrest lifting together, combined with air pressure massage to stretch the legs, remove fatigue, and recover your body.	F6 Auto

· Upper body manual massage

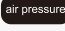

Function	Button	Massage ball adjustment	Display
Massage		Shiatsu Speed, width	
		Kneading: Speed	
		Tapping: Speed, width	
		Kneading and tapping: speed	
		Knocking: Speed, width	
Speed		3 adjustable level	 Speed
Width		3 adjustable level	 Width
Back stretch		The massage hand stay in fixed position	 H0
		The massage hand will move up and down in a partial position	 H1
		The massage hand will move up and down on full back	 H2
		The massage hand will move up and down on back and shoulder area	 H3
		The massage hand will move up and down on back and waist area	 H4
		Massage position adjustment (in fixed point, partial mode, the massage position can be adjusted).	Massage hand move up
	Massage hands go down		

Operation Instructions

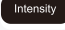

• Air pressure massage

Function	Button	Description	Display
Air pressure mode		Full air pressure: air pressure massage on arm, leg and foot.	Arm  Feet 
		Arm: Air pressure massage focus on arm.	Arm 
		Legs: Air pressure massage focus on leg and foot.	Feet 
Intensity		3 levels adjustable	Arm  Feet 

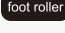

• Foot roller manual massage function

Function	Button	Description	Display
Foot roller		The foot roller massage function turn on or off.	




• Heating

Function	Button	Description	Display
Heating		Press this button to turn on/off the back infrared heating, and it will be warm after turning on 3 minutes.	

• Pause/start massage button


Function	Button	Description	Display
Pause/Start		Press this button during the massage, all massage actions are stop immediately (the timer will continue to count down), and press this button again, all the massage will be resume.	

• Backrest and calf position adjust

Function	Button	Description
Lifting button		Press and hold this button, it will raise the backrest and lower the footrest simultaneously, a "Di.Di sound indicate that it has reach the maximum position and stop. Press this button again while lifting, then immediately stop.
Lying button		Press and hold this button, it will lower the backrest and raise the footrest simultaneously, a "Di.Di sound indicate that it has reach the maximum position and stop. Press this button again while lifting, then immediately stop.
Automatic lying mode		Switch to an auto lying mode when you press key each time, total 3 modes, cycle switching change. Long press the zero gravity button of remote controller for 3 seconds to turn on or off voice control function.

Operation Instructions

• Other manual functions

Function	Button	Description	Display
Timing		Timing adjustment, press one time, add 5 minutes, not more than 30 minutes. If it exceeds 30 minutes, it will return to 5 minutes.	20 min


III.Operation Instructions for Bluetooth Speaker Play

- Enter the smart device main menu. Open "Settings" and select the "Bluetooth" option. Open "Bluetooth", and Bluetooth will automatically "Search for a New Bluetooth Device".
- The smart device will search for the corresponding Bluetooth model (e.g. IMCM-XXXX). Click the "Match" button. If you are asked whether to match with the device, you just select "Match" to complete matching.
- Turn on "Music Player" in the smart device and select the song to play.


IV.End massage

1. During the massage process, press the power key off touch screen controller to immediately turn off all massage functions, the backrest and leg rest will reset, or when the massage time expires, it will automatically turn off all massage functions, but the backrest and leg rest will not reset.


Cut off the power of the whole machine. Figure (switch power off)

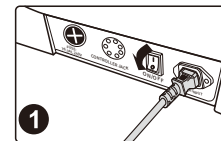

Power off ("0" position)

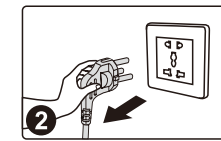

Take off the power plug

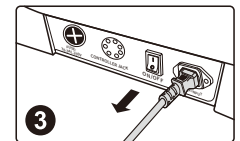

Take off the product power cable.

Cleaning and Maintenance

Clean PU leather and plastic parts

Wipe the leather with a soft, dry cloth.

- Do not use chemicals such as thinner, gasoline, alcohol, etc. to clean the unit.
- When using commercially available leather care products (rags), please follow the instructions for use.
- If the PU leather is particularly dirty, wipe it by the following methods.
 - ① Immerse a soft cloth in water or hot water containing 3% to 5% neutral detergent, then wring it out thoroughly.
 - ② Use the above soft cloth to wipe the leather surface.
 - ③ Wipe off the cleaning solution with a cloth soaked in clean water and thoroughly wrung out.
 - ④ Wipe with a wrung soft cloth.
 - ⑤ Let it air dry naturally.
- When it is difficult to remove stains, please immerse commercially available melamine foam in neutral detergent, and then wipe the machine with it.
- Do not use a hair dryer to dry the surface.
- As the color of the fabric will stain the surface of the PU leather, be careful when wiping the chair with denim and colored fabrics.
- Do not keep these parts in contact with plastic for a long time, as this will cause discoloration.
- PU leather may be discolored, so when using hair dye, please cover the touched part with a towel.



Troubleshooting

Fault	When the calf mechanism or the backrest reaches a certain position, press the corresponding adjustment button and it will no longer change, and a continuous "beep, beep" sound will be emitted.
Solution	Due to product structure design and ergonomic requirements, this product has a limited protection design for the limit state of attitude adjustment. At the same time, a reminder sound will sound when the limit is reached.

Fault	The function does not work or some functions do not work after the machine is turned on.
Solution	After the machine is turned on by pressing the switch key, it is necessary to manually select an automatic massage program or other working state. If it is not operated, the machine will automatically shut down after 20 minutes. If the manual operation mode is selected, the functions that are not selected may be disabled, and the working mode needs to be manually set in the manual mode.

Fault	The machine makes a leather rubbing sound when adjusting its posture.
Solution	This product adopts high-grade anti-corrosion leather that meets environmental protection requirements. Due to the inherent characteristics of leather materials, there will be a normal slight noise when rubbing at the relative movement, which is a normal phenomenon. Note: If you hear fabric or leather tearing sound when the machine is running, please stop using it immediately and seek professional maintenance.

Fault	mechanical hand did not reach the shoulders or neck.
Solution	If the head does not touch the head cushion or the back does not touch the backrest, the shoulder position may be lower than the actual position during the body shape detection process. Sit in the deepest position of the seat, rest your head on the headrest, and start over from the beginning.

Fault	The unit is damaged. The power cord or power plug is abnormally hot.
Solution	To prevent accidents, be sure to contact an authorized service center.

Please feel free to contact the designated service personnel for any questions or doubts about this product.

Product Specifications

Model: A302

Configure: 9

Name: Massage Chair

Rated voltage: 110-120V~ 60Hz

220-240V~ 50Hz/60Hz

Rated power input:120W

Rated time: 20Min

Safety structure: Class I

Material introduction: PVC, PA, steel parts and electric & electron parts

MATERIAL	PU	LEATHER	CLOTH	WOOD