
M.CURO

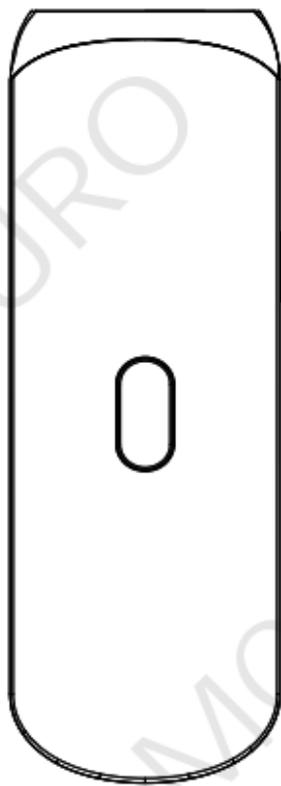IPL

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Deutsch

WILLKOMMEN

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MCURO

IPL Hair Removal Device


Product Model: BHRL-06
Colors: Rosa, Blau
Adapter Input: AC/100V-240V
Product Power: 48 W
Voltage: 12V 4A

Manufacturer: Bestway Plastic & Metal Products.,Ltd.
Address: Building B, Hongtengyuan Science and Technology Park, Mapu'ao Industrial Zone, Zhuweitian Village, Fenggang Town, Dongguan,GuangDong

MADE IN CHINA

M.CURO

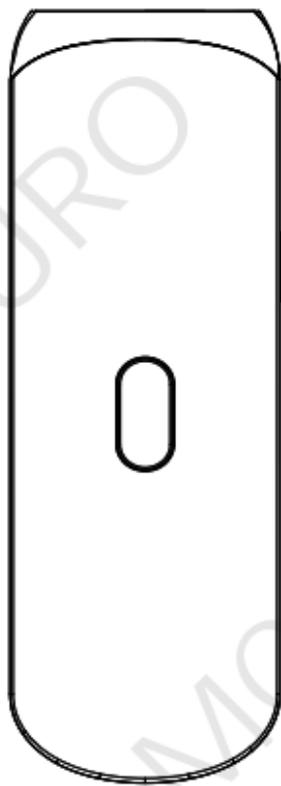IPL

HAIR REMOVAL DEVICE

English

WELCOME


The BHRL-06 hair removal device employs Intense Pulsed Light (IPL) technology, renowned for its effectiveness in reducing and preventing hair regrowth. By integrating the optical hair removal techniques commonly used in professional beauty salons, this device allows you to safely and conveniently achieve effective hair removal results in the comfort of your home. Enjoy soft, smooth skin and confidently showcase your natural beauty.

Important: Please read this user manual thoroughly before use.

CONTENTS

01. Welcome.....	01
02. Table of Contents.....	01
03. Device Structure.....	02
04. Important Safety Information.....	03
05. Contraindications.....	04
06. Common Skin Reactions.....	06
07. Product Overview.....	07
08. How It Works.....	08
09. Before Using Your Device.....	09
10. Steps for Usage.....	09
11. Troubleshooting.....	11
12. Device Maintenance.....	12
13. After-Sales Service.....	12

DEVICE STRUCTURE

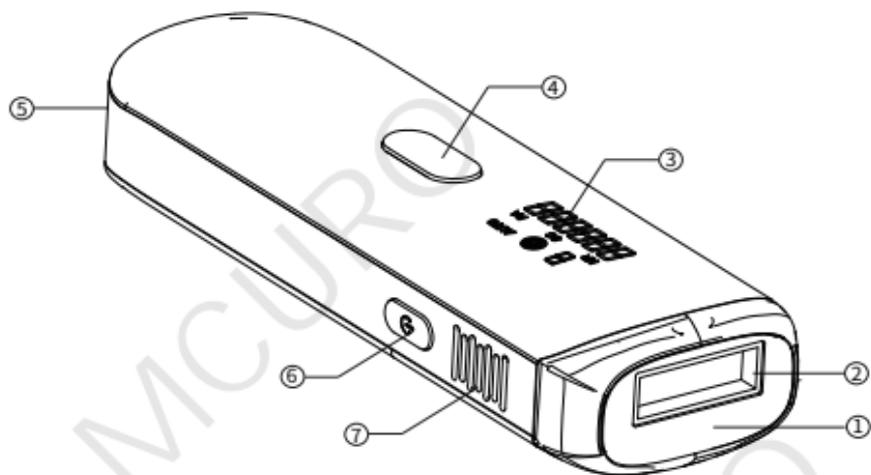- ① Ice panel ② flash window ③ Display screen ④ Flash button
⑤ DC socket interface ⑥ On/off button ⑦ Air outlet

IMPORTANT SAFETY INFORMATION

This device is not intended for use by individuals (including children) with reduced physical, sensory, or mental capabilities unless they are supervised or instructed by a person responsible for their safety.

Children should be supervised to ensure they do not play with the device. This device is not suitable for children under 15 years old. Teenagers aged 15 to 18 may use the device only with the consent and supervision of their parents or legal guardians. Adults aged 18 and above can use the device independently.

Always inspect the device before use. Do not use the device or its adapter if it is damaged. Replace any damaged parts with original components only.

The adapter includes a transformer. Do not replace the adapter with a different plug, as this may create a hazardous situation. Additionally, do not use the device if the UV filter on the light exit window is damaged.

This device is not waterproof. Do not immerse it in water or rinse it under running water.

For hygiene purposes, the device should only be used by a single individual.

Use the device only with settings appropriate for your skin type. Using settings higher than recommended may increase the risk of skin irritation or adverse reactions.

Only operate the device as intended and described in this user manual.

Do not clean the device with compressed air, abrasive pads, harsh cleaning agents, or aggressive liquids such as petrol or acetone.

This symbol indicates that the product must not be disposed of with regular household waste. Please follow your local regulations for the separate collection of electrical and electronic products. Proper disposal helps to prevent potential harm to the environment and human health.

Contraindications

Prohibition for Minors: This device must not be used by individuals under the age of 18.

Skin Type Restrictions: If your skin rarely or never sunburns, or if your skin tone is very dark, do not use this device.

Pregnancy and Breastfeeding: Pregnant or breastfeeding women should not use this device, as it has not been tested for safety during pregnancy or lactation.

Medical Conditions: Do not use this device if you have skin diseases, are undergoing cancer treatment, have early-stage cancer lesions, or have several atypical moles.

Sensitive Skin: Do not use this device if your skin is sensitive to light, prone to rashes, or susceptible to allergic reactions.

Recent Operations: Avoid using this device on areas that have undergone surgery within the past three weeks.

Health Conditions: This device should not be used by individuals with diabetes, lupus, porphyria, or congestive heart disease.

Skin Abnormalities: Do not use this device on moles, freckles, large blood vessels, darker areas, scars, or any skin abnormalities without consulting a doctor.

Sunburn and Tanning: Avoid using this device on sunburned skin, recently sun-exposed areas, or artificially tanned skin.

Affected Areas: Do not use the device on areas with tattoos, infections, eczema, psoriasis, sunburn, inflammation of hair follicles, open wounds, abrasions, Herpes Simplex, trauma, or other skin diseases.

SKIN COLOR CHECK

Applicable Skin Types: This device is suitable for Fitzpatrick skin type IV only.

The Fitzpatrick Scale is a classification system that divides skin tones into six levels based on their reaction to sunlight, ranging from very fair (Type I) to very dark (Type VI). Before using this device, it is essential to identify your skin tone to ensure safe and effective use.

Skin color Hair color	White	Beige	Light brown	Medium brown	Dark brown	Brown black and darker
	White/Grey	✗	✗	✗	✗	✗
Red	✗	✗	✗	✗	✗	✗
Light blond	✗	✗	✗	✗	✗	✗
Dark blond/light brown	regular mode power mode	regular mode power mode	regular mode power mode	soft mode regular mode	✗	✗
Brown	regular mode power mode	regular mode power mode	regular mode power mode	soft mode regular mode	✗	✗
Dark brown	regular mode power mode	regular mode power mode	regular mode power mode	soft mode regular mode	✗	✗
Black	regular mode power mode	regular mode power mode	regular mode power mode	soft mode regular mode	✗	✗

SUGGESTED TREATMENT INTERVALS

Initial Use:

It is generally recommended to use the device 1-2 times per week for 4-6 weeks, until significant hair reduction is visible.

Make sure your skin is clean and free of any skincare products before each use.

Maintenance Phase:

After completing the initial treatment, use the device 1-2 times per month to maintain the hair removal effect. As hair growth decreases, the frequency of use can be reduced accordingly.

Usage Duration:

Each session typically lasts around 10-20 minutes, depending on the size of the area being treated and the power setting of the device.

COMMON SKIN REACTIONS


Your skin may experience mild redness, tingling, or a warm sensation. These reactions are harmless and should subside shortly. After using pulsed light treatments for hair removal, your skin may feel dry or itchy. You can relieve the itching by applying ice or increasing moisture. If the skin remains dry, wait 24 hours after treatment before applying a fragrance-free moisturizer to the affected area.

PRODUCT OVERVIEW

The BHRL-06 epilator uses advanced technology with a combination of height and width, hyperspectral, non-continuous light sources to achieve both hair removal and skin rejuvenation effects.

IPL Technology Principles:

The device emits gentle pulses of light to the skin, which are absorbed by the hair roots. The darker the hair, the more effectively it absorbs the pulsed light. This light stimulates the hair follicles, causing them to enter the telogen phase, leading to natural hair shedding and new hair regeneration. As a result, the skin becomes softer and smoother.

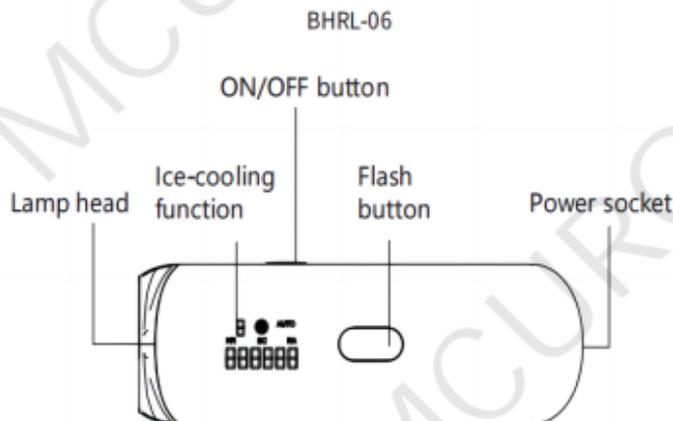Rated voltage: 12VDC/4A

Input power: 48W


Output light energy: 4.5-12J/ cm²

Working temperature: 5-35°C

Working humidity: 30%-85%

Storage temperature: 0 -40 °C

HOW IT WORKS

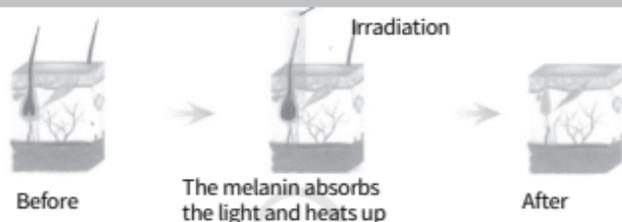Hair removal

The IPL emitted by the device is absorbed by the melanin (pigment) in the hair follicles. As the pigment absorbs the light energy, it heats up, damaging the hair follicle without affecting the surrounding skin. This heat leads to the destruction of the follicle, preventing future hair growth. The damaged cells are then naturally removed by the body's metabolic processes, resulting in effective hair removal.

Skin Rejuvenation

The IPL emitted by the device generates a thermal effect on the skin, which accelerates blood circulation and stimulates fibroblasts to produce collagen and elastin fibers. This process enhances skin elasticity, helping to reduce wrinkles and achieve a rejuvenated, smoother appearance.

Acne Removal

Through photochemical action, porphyrin produced by *Propionibacterium acnes* interacts with triplet-state oxygen to form a structure. The unstable singlet oxygen produced in this process effectively kills *Propionibacterium acnes*. IPL technology is also effective in treating acne and reducing scars, promoting clearer skin.

BEFORE USING YOUR DEVICE

Before using the BHRL-06, it is important to prepare your skin by removing any surface hair. This allows the light to be absorbed by the hair beneath the skin's surface, ensuring effective treatment. You can shave, trim, epilate, or wax the area.

Avoid using depilatory creams, as the chemicals may cause skin reactions. If you choose to wax, wait 24 hours before using the BHRL-06 to allow your skin to rest.

We also recommend taking a shower before the treatment to ensure any remaining wax residue is completely removed from your skin.

Pretreat the areas you intend to treat with the BHRL-06.

Clean your skin thoroughly, ensuring it is hair-free, completely dry, and free of any oily substances.

STEPS FOR USAGE:

1. Pre-Usage Instructions:

Before using the device, shave the hair off using a razor, then dry your skin. Keep your skin moisturized to prevent burns caused by dry skin. Wear the provided goggles to protect your eyes. Do not directly expose the device to black clothing, as this may damage the device. The darker the skin, the hotter it may feel during the treatment. Start with the lowest setting (1st gear) and gradually increase the intensity.

2. Powering On and Initial Setup:

Press and hold the switch button for 3 seconds. The display will light up with "HR" and the 1st gear indicator. You will also hear sounds from the suction and exhaust ports, indicating the device is operating normally.

3. Selecting Functions and Adjusting Intensity:

Double-click the switch button to select the desired function based on your needs. Press the switch button to adjust the light intensity. Gear 1 is the lowest, and Gear 9 is the highest. You can cycle through the intensity levels.

4. Using the Device for Hair Removal:

After setting the intensity, place the light outlet of the device on the area you want to treat. The selected function indicator will flash on your skin to confirm the light outlet is in contact with the skin. You can then press the light button to activate the flash.

5. Flash Mode Options:

Press and hold the flash button for 3 seconds to toggle between automatic and manual flash modes. When the "AuTo" indicator light is on, the automatic flash mode is active. In this mode, the device will flash automatically every 2-3 seconds. When the "AuTo" indicator is off, manual flash mode is active, and you need to press the flash button to trigger each flash. To turn the ice cooling mode on or off, press both the switch button and the flash button simultaneously.

For effective hair removal, each treatment area should be flashed once, then repeat the process twice for a total of three flashes. Do not flash the same area continuously to avoid overheating and potential burns.

6. Powering Off:

After use, press and hold the switch button for 3 seconds to turn off the device. Disconnect the power supply immediately. The device will automatically shut down and display "oFF" after 5 minutes of inactivity.

TROUBLESHOOTING

Problem	Possible Causes	Solutions
The device and power converter may generate varying temperatures during use.	This is a normal phenomenon	Does not indicate a malfunction
The buzzer doesn't sound, but the cooling fan is running when the device is placed on the skin.	Overheating protection has been activated.	Allow the device to cool for about 15 minutes without powering it off, then resume use.
The device does not emit light when you press the Flash button, but everything else appears normal.	The device requires a reset.	Turn off the device, unplug it, and wait a moment before powering it back on.
Unable to increase or decrease the light intensity setting.	The device requires a reset.	Turn off the device, unplug it, and wait a moment before powering it back on.
There is an odor when using the device, or there is intolerable pain on the skin during treatment.	<ul style="list-style-type: none">- Accessories or filters are too dirty.- Hair in the treatment area has not been fully removed.- Light intensity is set too high.- UV filter on the light outlet is damaged.	<ul style="list-style-type: none">- Carefully clean all parts of the device to remove dirt or debris.- Ensure hair is removed from the skin's surface before use.- Adjust the light intensity to an appropriate level for your skin.- If the UV filter is damaged, stop using the device and contact customer support.
Skin reactions last longer than usual after treatment.	<ul style="list-style-type: none">- Light intensity was set too high.- Light intensity was set too low for effective treatment.- The treatment area was not fully covered during the session.- IPL treatment response may vary.	<ul style="list-style-type: none">- Next time, use a lower light intensity to avoid prolonged skin reactions.- If results are not effective, try a higher light setting.- Ensure the light covers the full treatment area during use.- Continue regular treatments for at least 6 months, as hair regrowth will decrease over time.

DEVICE MAINTENANCE

1. Clean the lamp holder by removing any impurities from the light window before use.
2. Place the device near a wall outlet, ensuring it is kept away from liquids or damp areas.
3. Before turning on the device, confirm that the lamp holder is securely installed.
4. After use, turn off the device, unplug it, and allow it to cool completely.
5. Once cooled, use a clean, slightly damp cloth to wipe the following parts: the light window, the sides of the device, and the reflector in the accessories.
6. Thoroughly dry all parts before storing the device.
7. Store the device in a dust-free, dry place with a temperature range between 0°C and 60°C.

AFTER-SALES SERVICE

We are committed to providing you with high-quality products and excellent service. Within two years of purchase, if any damage occurs due to manufacturing defects or component issues, you are eligible for free repair or replacement.

As this product comes into direct contact with the skin, to ensure the hygiene and safety of all users, we do not resell used products. Therefore, we encourage you to use the product for a period of time to observe its effects. If you are still not satisfied after using it, please contact our customer service to request a refund. You do not need to return the product. We hope this service gives you peace of mind and confidence in using our product.

SAFETY PRECAUTIONS


Ensure both the device and the adapter remain dry at all times.

If the device is damaged, avoid touching any internal components to prevent the risk of electric shock.

Do not insert any objects into the device under any circumstances.

BODY PARTS FOR USE IN TREATMENT

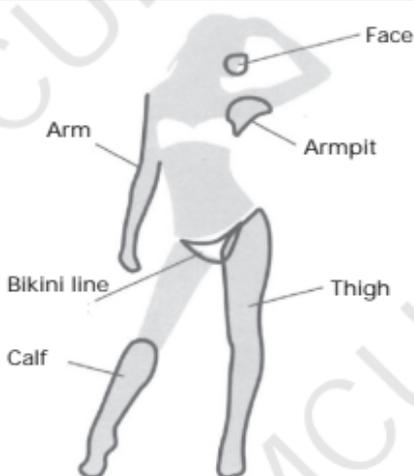MCURO

IPL Hair Removal Device

Product Model: BHRL-06

Colors: Rosa, Blau

Adapter Input: AC/100V-240V

Product Power: 48 W

Voltage: 12V 4A

Manufacturer: Bestway Plastic & Metal Products.,Ltd.

Address: Building B, Hongtengyuan Science and Technology Park, Mapu'ao Industrial Zone, Zhuweitian Village, Fenggang Town, Dongguan,GuangDong

MADE IN CHINA