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**+49 (0) 6431 9 21 96 00**

## Content

Product components.....	15
General information.....	16
Explanation of symbols.....	16
Safety.....	16
Proper use.....	16
Safety and warning instructions.....	16
Application of the accessories.....	19
Before using for the first time.....	19
Charging the battery.....	20
Battery indicator.....	20
Filling the water tank.....	20
Mouthwash.....	20
Use.....	20
Brushing teeth.....	21
Mouthwash.....	21
Teeth brushing + mouthwash at the same time.....	21
Using the tongue cleaner.....	21
30-second pause.....	22
2-minute timer.....	22
Using the appliance when travelling.....	22
Cleaning and storage.....	22
Troubleshooting.....	23
Disposal.....	24
Technical details.....	25
Customer service.....	25

## Product components

- 1 Brush head
- 2 Mouthwash nozzle
- 3 Storage box
- 4 Water tank
- 5 Water hose
- 6 Water hose bracket
- 7 Toothbrush on/off button
- 8 Toothbrush mode indicator
- 9 Mouthwash on/off button
- 10 Mouthwash water pressure button
- 11 Mouthwash water pressure indicator
- 12 Charging indicator
- 13 Charging station
- 14 USB-C charging socket
- 15 Charging contacts
- 16 Charging connector
- 17 Travel case
- 18 USB-C charging cable

BEFORE USING THIS APPLIANCE FIRST READ ALL INSTRUCTIONS AND SAFETY PRECAUTIONS IN THIS INSTRUCTIONS FOR USE.

KEEP THIS MANUAL FOR FUTURE REFERENCE.

ONLY GIVE THE APPLIANCE TO THIRD PARTIES TOGETHER WITH THE INSTRUCTIONS FOR USE.

WHEN USING ELECTRICAL APPLIANCES, THE FOLLOWING BASIC PRECAUTIONS MUST BE OBSERVED AT ALL TIMES, INCLUDING THE FOLLOWING:

These instructions for use apply to Genius Dental.

The illustrations in these instructions for use are intended to schematic representation of the correct use. They may differ from the original product in detail.

These operating instructions cover all available accessories for the appliance. The content of the set that you have purchased is indicated on the packaging or a separate insert. You can also purchase additional accessories which are not included in your set.

## General information

### Explanation of symbols

The following symbols and signal words are used in this manual, on the device or on the packaging.

**WARNING!**

Failure to observe this warning could lead to a serious or potentially lethal accident.

**CAUTION!**

Failure to observe this caution could lead to injury.

**NOTE!**

Information about potential material damage.

## Safety

### Proper use

- Only use the appliance for cleaning teeth, gums and tongue.
- Use the device only as described in this manual. Any other use is regarded as improper and may result in damage to property or even personal injury. The device is not a toy for children.
- The manufacturer and distributor accept no liability for damages resulting from improper or incorrect use. Safety instructions
- Only use the unit with original accessories, otherwise you will damage the appliance and increase the risk of injury.
- Do not use a brush head with kinked or crushed bristles. Replace the brush head every 3 months or sooner if there are signs of wear.
- Stop using the appliance and consult your dentist/doctor if excessive bleeding of the gums occurs or if bleeding of the gums persists for more than one week or if you experience discomfort or pain during use.
- If you have had a surgical procedure on your teeth or gums in the past 2 months, consult your dentist before using the appliance.
- If you have any concerns about your health, consult your doctor before using this appliance.
- This appliance is a device for individual oral care and is not suitable for use by multiple patients in a dental practice or dental clinic.

## Safety and warning instructions

**WARNING!**

**RISK OF ELECTROCUTION!**

Improper electrical installation or excessive mains voltages can result in electric shock.

- Only connect the device if the mains voltage of the socket matches the information on the nameplate of the power supply unit.

- Only connect the power supply to an easily accessible plug socket, in order that it can be quickly disconnected from the mains in case of a fault.
- Do not operate the device if it displays visible damage or if the mains adapter is defective.
- Only use the included power adapter to charge the device.
- The power cable of the power supply cannot be replaced. Dispose of the power supply if it is damaged and replace it with a power supply of the same type.
- Do not under any circumstances open the housing of the device: repairs must be left to specialists. For repairs, visit a specialist workshop.
- Never touch the mains plug with wet hands.
- Do not immerse the unit in water.
- Never pull the power supply from the socket by the power cable, but rather always grip the plug itself.
- Never use the power cable as a carry handle.
- Keep the device and mains adapter away from open flames and hot surfaces.
- Do not kink the power cable and do not run it over sharp edges.
- Do not leave the device unsupervised whilst in operation.
- When you are not using the device, for cleaning or in the event of a fault, always turn the device off and remove the mains plug from the socket.
- Make sure that you can always check from your seat that the mains plug is still disconnected from the socket.
- Store the unit so that it cannot fall into the water at any time.
- Do not use the appliance when it is connected to the adapter or the mains.
- Misuse is a hazard and can lead to serious injury.
- Charge and use the product at temperatures between 5°C and 40°C.
- Keep the appliance and batteries away from fire and do not expose it to direct sunlight or high temperatures.
- If the product becomes unusually hot, develops an unusual odour, changes colour or if charging takes much longer than usual, stop using and charging the product and contact customer service.
- If the batteries are damaged or leaking, avoid contact with the skin or eye. If this is the case, immediately rinse the corresponding area thoroughly with water. Seek medical help.

**WARNING!**

**DANGER FOR CHILDREN AND PEOPLE WITH LIMITED PHYSICAL, SENSORY OR MENTAL CAPACITIES (FOR EXAMPLE PARTIALLY DISABLED PEOPLE, OLDER PEOPLE WITH LIMITED PHYSICAL AND MENTAL CAPACITIES) OR THOSE LACKING IN EXPERIENCE AND KNOWLEDGE (E.G. OLDER CHILDREN).**

- Do not allow children aged under 8 years or people with physical, sensory or mental handicaps or a lack of experience or knowledge use the appliance without supervision or instruction on the safe use of the appliance and potential risks. Never allow children to play with the appliance.
- Keep children under 8 years of age away from the device and mains adapter.
- Do not allow children to play with the packaging film. Children could become entangled or suffocate when playing.
- Make sure that children do not insert any objects into the unit.

**WARNING!****RISK OF EXPLOSION AND CHEMICAL BURNS!**

There is a risk of explosion if batteries and accumulators are not handled correctly. In case of contact with leaking battery acid, chemical burns may occur.

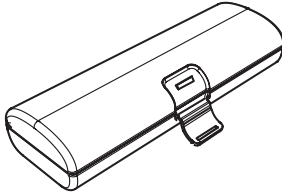
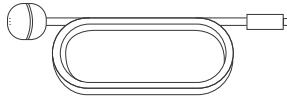
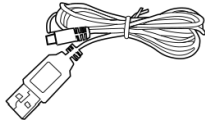
- Do not heat the battery and never throw it into an open fire. Do not destroy or open the battery, do not submerge it in liquid, and do not short-circuit it.
- Do not expose the battery to direct sunlight, strong temperature fluctuations or temperatures above 60 °C.
- Avoid battery acid from coming into contact with the skin, eyes and mucous membranes. In case of contact with battery acid, rinse the affected areas immediately with plenty of clean water and if necessary consult a doctor

**CAUTION!****RISK OF INJURY!**

Failure to use the device correctly can result in injury to you or others.

- Do not aim the device at people or animals.

## Application of the accessories

Accessories	Image	Use
Travel case		For safe transport of the toothbrush including accessories.
Water hose		For using the water flosser on the go.
Charging cable		For direct charging at the toothbrush.
Brush head (sensitive/standard) Brush head (sensitive/large) Brush head (medium/standard) Brush head (medium/large) Brush head (intensive/standard) Brush head (intensive/large)		Interchangeable brush heads in different sizes and degrees of hardness.
Tongue cleaner		For cleaning the tongue.

## Before using for the first time

**TIP!**

### DANGER OF DAMAGE

If you open the packaging carelessly with a sharp knife or other pointed objects, the appliance may be damaged.

- Open the packaging carefully.
  1. Remove the appliance from the packaging.
  2. Remove the packaging material and any labels on the product that are used for protection during transportation or advertising.
  3. Check that everything has been included.
  4. Place the product so that it is on an even surface and cannot fall into a sink or bathtub.
  5. Check whether there is any apparent damage to the appliance or its accessories. If this is the case, do not use the appliance. Contact the customer service address given in the Customer Service section.
  6. Fully charge the appliance before using it for the first time.

### Mounting/Dismounting the brush head and main unit

1. Put the brush head on the main unit. Press firmly until it clicks into place. (See Fig. A)
2. To remove, pull firmly on the brush head to separate it from the main unit.

### Connecting and disconnecting the water hose and the main unit

1. To connect, insert the water hose firmly into the water inlet on the back of the main unit.
1. To disconnect, pull out the water hose, and immediately insert the water hose into the small hole on the bottom of the charging station to prevent water leakage from the tank. (See Fig. B)

## Charging the battery

Make sure that the main unit is switched off

### Option 1: Charging in the charging station with water hose

1. Insert the mains plug into the charging socket on the back of the charging station. (See fig. D).
2. Plug the mains adaptor into a functioning socket.
3. Dry the contact surface of the main unit and the base.
4. Place the unit in the charging station so that charging contact is established.

### Option 2: Direct charging of the main unit

1. Plug the Type-C USB cable into the charging socket on the underside of the main unit. (See fig. E)
2. Plug the mains adaptor into a functioning socket.

## Battery indicator

**TIP!**

When using the appliance for the first time or if it has not been used for more than 6 months, the battery indicator may not light up at the beginning of the charging process. In this case, please charge the appliance for at least 6 hours.

- When the battery is low, the battery indicator flashes red.
- During charging, the battery indicator lights up red.
- When the battery is fully charged, the battery indicator lights up green.

## Filling the water tank

1. Remove the storage box.
2. Fill the water tank with cold or lukewarm water.
3. Put the storage box back on the water tank.

## Mouthwash

1. Hold the appliance so that the handle and nozzle attachment point into the water basin.
2. Press the floss button to get the water flowing.
3. After water comes out of the nozzle, deactivate the water flow by pressing the floss button again.
4. Once you get used to using the function, you can continue using it.

## Use

**TIP!**

In cleaning mode, the hose can be disconnected from the main unit. There is also no need for water in the water tank.

## Brushing teeth

1. Put the brush head in your mouth before pressing the button.
2. Press the brush button to start the main unit.
3. By pressing this button for a longer time (approx. 3 seconds), you can select between the different cleaning speeds: Low, Medium and High. The selected speed is indicated by a green LED light under the button:
  - a. Low: The left LED lights up
  - b. Medium: The middle LED lights up
  - c. High: The right LED lights up
4. By pressing this button again, you can switch off the device.

### Recommendation

For best results, hold the brush head at a 45-degree angle parallel to the gumline. Switch on the toothbrush. Move the toothbrush slowly and with light pressure along the teeth. This ensures each individual tooth gets cleaned.

## Mouthwash

1. Always put the brush head in your mouth before pressing the floss button.
2. Press the floss button to get the water flowing.
3. Bend over the sink and close your mouth so that the water does not splash out, but can run into the sink.
4. Adjust the water pressure by pressing the speed button. Before use, set the water pressure to the lowest level. Gradually increase the pressure until you reach the desired setting or the setting recommended by your dentist. The water pressure settings are shown by LEDs under the speed button (level 1-5).
5. Deactivate the water flow by pressing the floss button again.

### Use mouthwash and other solutions

The device can also be used together with mouthwash and antibacterial solutions. After using a solution, flush the appliance to avoid clogging the tubing. For rinsing, replace the solution-and-water mixture in the water tank with clean tap water. Operate the mouthwash to rinse the hose and brush head with clean water.

### Recommendation

For best results, align the brush head so that the water nozzle points slightly towards the gumline. Move the brush head along the gumline. Pause briefly after each tooth so that the interdental space can be cleaned by the water jet. Proceed from the molars to the incisors until all teeth in the upper and lower jaw have been cleaned. Then empty the tank.

## Teeth brushing + mouthwash at the same time

1. First place the brush head in your mouth.
2. Press the brush key followed by the floss button.

## Using the tongue cleaner

1. Plug the tongue cleaner onto the main unit
2. Turn on the main unit and brush your tongue with a slow and steady forward and backward motion.
3. Rinse your mouth out.

## 30-second pause

### TIP!

The Genius Dental pause means that there is a short pause every 30 seconds to remind you to move the toothbrush on to the next area of the mouth. Brush the outside and then the inside of your teeth for 30 seconds at a time. This process should then be repeated on the lower and upper jaws, respectively, so that all your teeth are cleaned within two minutes.

## 2-minute timer

### HINWEIS!

The brush turns off automatically after 2 minutes of running time.

## Using the appliance when travelling

The Genius Dental can be separated from the base and used when travelling, both as a toothbrush and as a mouthwash.

1. When travelling, transport the main unit and the desired attachment brushes and accessories in the travel case supplied.
2. Get a glass and fill it with water.
3. Insert the end of the water hose with the gravity ball into the glass and connect the other end to the appliance.
4. You can also use the Genius Dental without a tube as a travel toothbrush.

## Cleaning and storage

### WARNING!

#### Risk of electric shock

Switch the appliance off and disconnect it from the mains before cleaning or user maintenance.

### TIP!

The product parts must not be washed in the dishwasher. Store the appliance in a safe place and protect it from excessive heat or direct sunlight. Keep the charger away from water.

#### 1. Clean the attachment brush:

When all the functions are switched off, remove the attachment brush from the main unit after each use and clean the brush head under running water.

#### 2. Cleaning the appliance:

- a. Clean the main unit by wiping it with a slightly damp cloth.

- b. If necessary, clean with a mild, non-abrasive detergent.
- c. Dry the main unit every time you use it/clean it.

### 3. Removing limescale deposits:

Depending on the hardness of the water

- a. Add two tablespoons of vinegar to the water tank with warm water.
- b. Point the main unit and the brush head into the water basin.
- c. Switch on the appliance and let the tank run empty.
- d. Repeat the rinse with a full tank of warm, clean water.
- e. You can repeat this process every one to three months.

## Troubleshooting

### Only for electrical appliances

### WARNING!

Switch off the appliance and disconnect the mains plug before addressing any faults! Failure to do so may result in injury or electric shock.

Problem	Possible causes	Solution
Appliance isn't working	Battery is empty.	Charge the appliance.
Battery will not charge	Appliance or charging station is not correctly connected to a power socket.	Plug the mains adaptor securely into a functioning power socket
	Contact with charging station has not been established.	Position the hand-held appliance correctly in the charging station Check that the main unit and the charging station are dry.
Water leakage from the water hose	Damaged water hose.	Water hose needs to be replaced.
	Water hose not properly connected to the main unit.	Connect the water hose into the main unit firmly and check it is properly connected.
Water leakage from the water tank	Leaks.	Water tank needs to be replaced.
Water remaining in the water tank.	Hose blocked.	Clean the unit (see "Cleaning and storage").
Water pressure low	Not enough water in the water tank.	Fill the water tank.
	Wrong water pressure setting selected.	Adjust the water pressure settings.

### ATTENTION!

If the fault cannot be rectified, please contact our customer service. Never open up the appliance!

## Disposal

(Applicable in the European Union and other European states with systems for the separate collection of waste and recyclables).

### Packaging

The packaging material can be recycled. Dispose of the packaging in an environmentally friendly way for recycling.

### Electrical appliances

The symbol with the crossed out dustbin means that this product cannot be disposed of in the household waste. Electrical waste should be taken to a local acceptance point.

You can obtain further information about the disposal and recycling of this appliance from your municipal utility company, the retailer where you purchased it or the manufacturer.

Before disposing of the appliance remove the batteries because they have to be disposed of separately. To protect the environment, batteries and rechargeable batteries must not be disposed of alongside normal household waste, regardless of whether they contain harmful substances\* or not, and should instead be taken to the relevant collection point. Contact your local authority to find out where you can dispose of batteries in your area.

Even empty batteries have a residual charge. Before you dispose of them, cover the poles with insulating tape to prevent them short circuiting, which could cause fire or an explosion.

\*marked with: Cd = Cadmium, Hg = Mercury, Pb = Lead

## Technical details

<b>Model</b>	<b>Genius Dental Hydro Fusion</b>
Size (H x W x D)	Ø38.8×235 mm (Genius Dental Hydro Fusion) 128×118×198 mm (charging station with water hose)
Battery	22,2 V, Lithium-ion battery
Battery capacity	2200 mAh
Water tank capacity	0,5 l
<b>Mains adaptor</b>	
Manufacturer (name, address)	Shenzhen Wanyou Electronic Technology Co., Ltd. 3rd Floor, Zhenhua Electronics Factory, Shayi West Industrial Zone, Shayi Community, Shajing Street, Baoan District, SHENZHEN, Guangdong 518000 China
Model	M4-0510CU
Input voltage / AC frequency	100-240V ~ 50/60Hz, 0.3A Max
Output voltage / current / power	5V 1.0A, 5W
Average active efficiency	78.70%
Efficiency at low load (10 %)	35%
No-load power consumption	0.1W

## Customer service

If you have any questions about the device or spare parts/accessories, please contact your local retailer.

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