
Neck Beauty Device

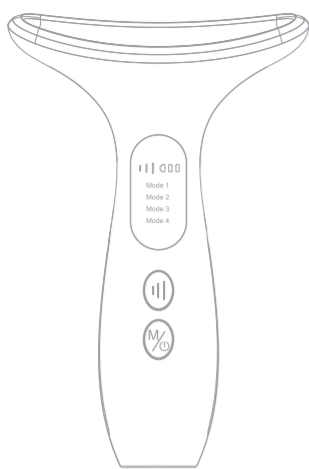Product Manual

Thank you very much for purchasing and using this product

Product Specifications	
Product name	Neck Beauty Device
Product number	HY38
Product Size	50*92*135mm
Rated power	5W
Rated voltage	3.7V/400mAh
Power adapter	DC5V =1A(USB)
Charging time	3 hours

Button definition

- ① Power on/off button: Long Press to Switch on/off Short Press to Switch Mode
- ② Gear key
- ③ Color light
- ④ DC charging port
- ⑤ Gear display
- ⑥ Battery display

Product Maintenance

- Please try to keep this product in a dry and room temperature environment or in a storage box.
- Do not overcharge to avoid damaging the battery life.
- Please wipe the product clean after charging and put it in the storage box if not used for a long time.
- The outside of the product can be wiped clean with alcohol or diluted detergent.
- During use, the device will become slightly hot, but it is normal for it not to be hot.
- It is forbidden to rinse the fuselage directly with water.
- If it is not used for a long time, please charge the battery at least once every 3 months.

Usage warning

- It is prohibited to use it with medical electrical appliances implanted in the body, such as pacemakers, as it may cause equipment failure or trigger danger.
- Do not use this product if your palms have skin problems such as wounds, eczema, blisters, and peeling.
- Do not use this product in places with high humidity such as bathrooms.
- Please do not use this product for other purposes than its intended use.
- Please do not disassemble, repair or modify this product without permission. May cause equipment malfunction or cause injury.
- Please do not give this product to children under six years old to avoid personal injury.

Packing list

Neck beauty device power cable manual

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